

TITEMS Monitor High Precision


– Type 4 Indicator

Titems Monitor HP (High Precision) are designed for use in steam sterilizers, providing an integrated response to defined combinations of time, temperature and presence of moist saturated steam.

These attributes are monitored to a point that is well in excess of biological spore kill where upon the indicator spot will change from blue to black.

Product Highlights:

- Type 4 performance
- Sterilant specific (only changes in steam)
- Easy interpretation
- Can monitor the mechanical aspect of the boiler, steam lines and sterilizer
- Can monitor the quality of the preparation of individual packs, the loading of the chamber and the use of the steam in other areas of the hospital
- Requires correct time & temperature with the presence of steam
- Accuracy through purity and quality of formulation

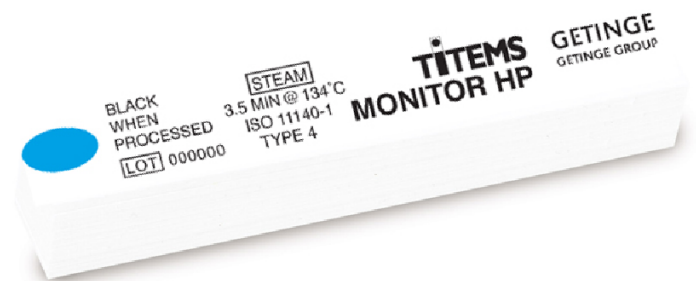

Storage:

- Store at room temperature between 10-38°C and normal humidity (10-70%)
- The expiry date is printed on the product packaging, 3 years from date of manufacture

Re-order Code	Description	Packaging
AUS410	TITEMS Monitor HP - Short Version	2000/Packet
AUS415	TITEMS Monitor HP - Short Version	125/Packet
AUS420	TITEMS Monitor HP - Long Version	1000/Paclet

Getinge is a global provider of innovative solutions for operating rooms, intensive care units, sterilization departments and for life science companies and institutions. Based on our firsthand experience and close partnerships with clinical experts, healthcare professionals and medtech specialists, we are improving the everyday life for people - today and tomorrow.

Australian Office - Level 2, 4 Talavera Road, Macquarie Park, Australia NSW 2113 • **Phone:** 1800 438 4643 • **Email:** sales.au@getinge.com

New Zealand Office - 34 Vestey Drive, Mt Wellington Auckland 1060 • **Phone:** 0800 1 438 4643 • **Email:** sales.nz@getinge.com