



Maquet Modueco Ceiling Supply Units

Facilitating effective patient care in
various compact healthcare environments

GETINGE 



Maquet Modueco

Facilitating effective patient care in various compact healthcare environments


Maquet Modueco ceiling supply unit offers high load capacity and flexibility, ensuring an optimal medical work environment even in space-constrained settings like:

- + Operating Rooms
- + ICUs and NICUs
- + Emergency rooms
- + CSSDs

With Maquet Modueco, benefit from:

- + Ergonomic solutions tailored for medical staff
- + Secured workstations for patients, medical staff and equipment
- + Upgradable medical workstation for future needs

Maquet Modueco Overview

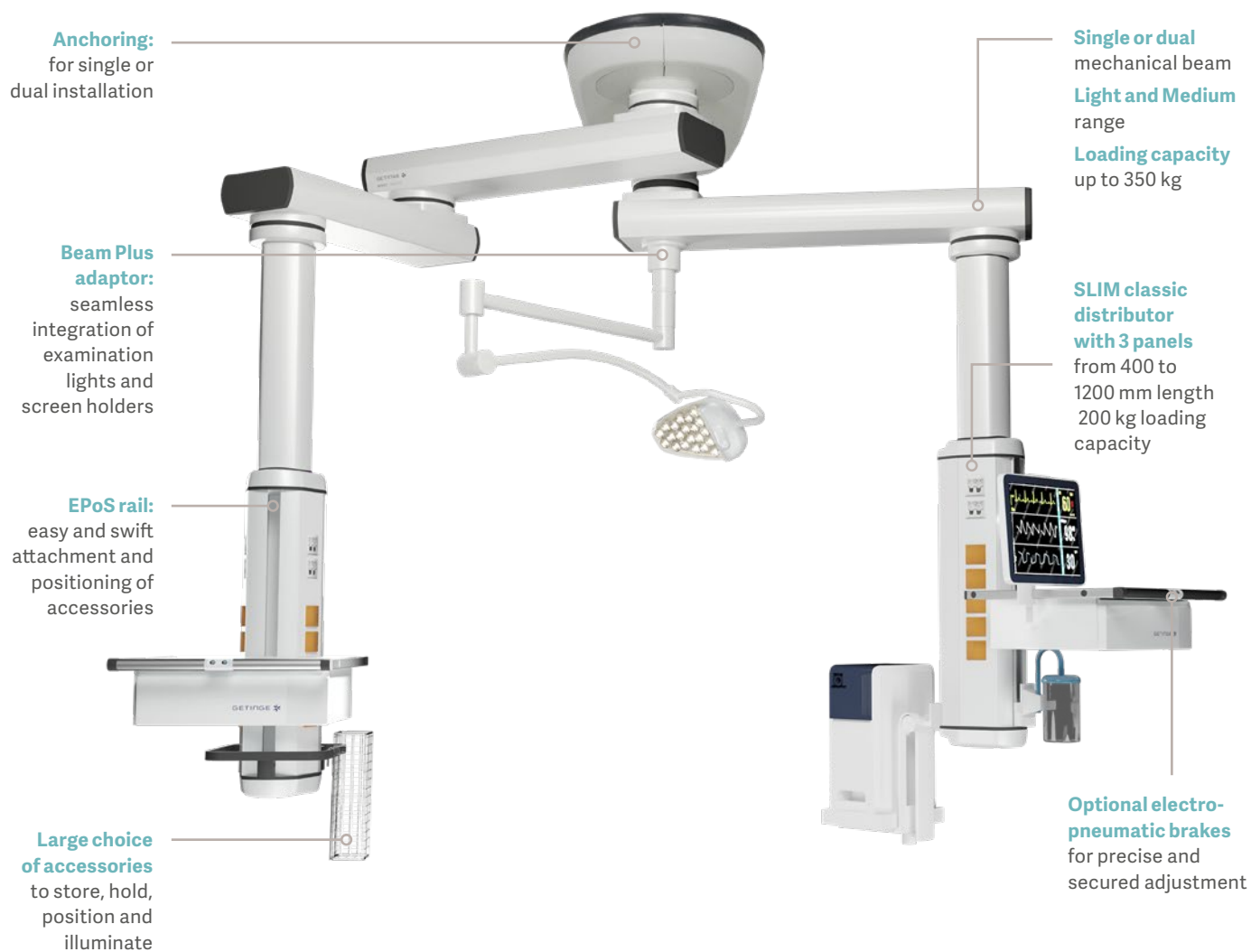

Maquet Modueco Tailored for various medical environments

In Intensive Care Units (ICU)

Maquet Modueco provides a **healing environment and tailored solution to enhance ergonomics**. The system allows for a clear organization between wet and dry sides, helping clinicians to achieve better space optimization while ensuring a **smooth integration of all critical care equipment** (ventilator, monitoring, IV pumps...).





In Operating Rooms

Whether in anesthesia, surgical, or endoscopic settings, Maquet Modueco provides effective support across all surgical workstations. **With a high load capacity** up to 350 kg and **flexible cable management solutions**, it enables to create effective minimally invasive surgery workstations.



In Neonatal Intensive Care Units (NICUs)

The design of the column allows for **more efficient use of space**. Thanks to its angled panels, **it offers an optimal access to sockets and outlets** while ensuring that treatment and monitoring technologies are within reach for medical staff to treat neonatal patients.



In Emergency Rooms

Time is critical. Medical professionals must act swiftly to keep patients alive and stable, without being hindered by equipment. The Slim classic distributor **maximizes space efficiency** for essential medical devices, ensuring patients remain stable and supported. It creates room for medical staff to treat or assess patients quickly and efficiently. In addition, the swivel beam and pneumatic brakes enable **quick and secure positioning** of various equipment.



In Central Sterile Services Department (CSSDs)

Maquet Modueco provides power supply and carrying capacity for **standardized and hygienic workstations**.

Thanks to its cable management solution and angled-panels, it ensures a **safe environment** while creating an **ergonomic working environment** for CSSD staff.



Easy and ergonomic
access to patient's head

Maquet Modueco Ergonomic solutions

Ergonomy plays a major role
in improving nurse satisfaction
and patient care.

Maquet Modueco has been
designed in order to facilitate
medical staff working conditions.



Maquet Modueco Secured workstations



Providing safer medical
working conditions is crucial
to improve patient care.



Electro-pneumatic brakes reduce the risk of unintentional
movement in order to protect patients, medical staff and
equipments. The brake release mechanism ensures easy
operation.



The combination of a **flexible, controllable braking system** and an **ergonomic operating handle** ensures **precise positioning, smooth and effortless movement**. **Simple to move**: brakes with compressed air are released when buttons are pressed.



Angled panels provide an **optimal access to electrical socket and gas outlets**.



Ambient lights available at the top of the beam, on the sides or at the bottom of the distributor **allow medical staff to keep on working without disturbing the patient in intensive care units**. In the operating room, It enhances visibility and creates a more comfortable environment for surgeons, especially in low-light laparoscopic procedures.



Hygiene with smart cable management.



Robust and reliable connection between the hose and the gas outlet. Flammable resistant cables.



Improving safety for patients and clinicians thanks to the integrated O₂ concentration monitor. It provides a **real time O₂ concentration** and is equipped with flammable resistant plastic cover.

Maquet Modueco

Upgradable medical workstations for future needs



+ A reservation can be included for future upgrades, helping to save installation time and costs.



+ The Beam Plus adaptor enables the integration of an examination light, a surgical light, or a flat-screen holder onto the beam right from the start of the project. Its independent rotation ensures easy and safe maneuverability. Ceiling supply unit and light share the same anchoring, no need for a second one, saving costs and space on the ceiling.



+ The distributor is equipped with the EPoS* system which includes three vertical rails, enabling:

- the flexibility to create ideal, customized workstations with the possibility to install accessories at the front and at the back sides,
- quick and easy attachment and positioning of accessories on every side and along the entire height.

* Without power



Ensuring your workstation can evolve along with your needs


Maquet Modueco Configure your workstation according to your needs




Column
configuration



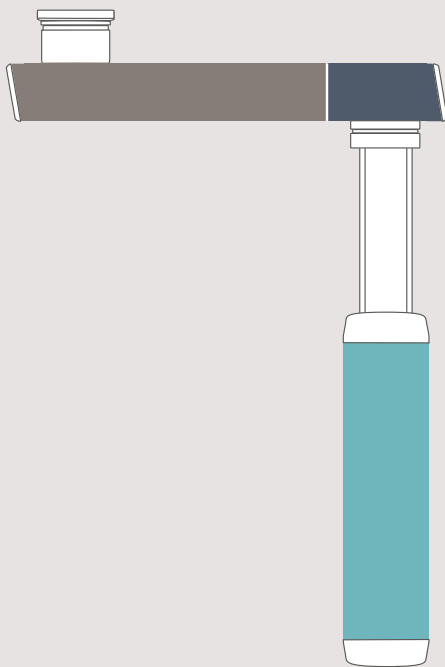

Single beam
configuration

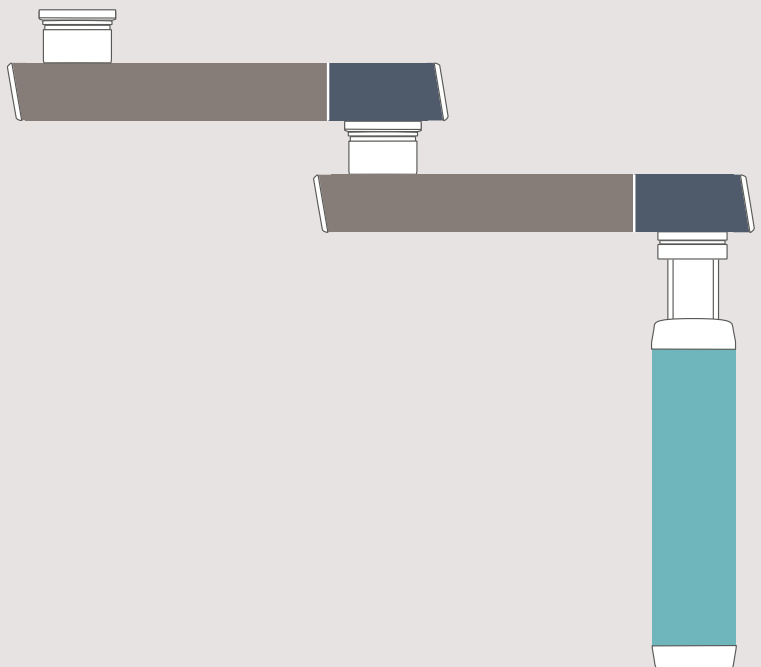

Double beam
configuration



Distributor length: from 400 mm up to 1200 mm
Beam length: from 600 mm up to 2400 mm
Column length: from 600 mm up to 1600 mm

Light/Medium 
Medium 

Maquet Modueco

A wide choice of accessories



- + Choice between standard, height-adjustable and strong IV poles



- + Control shelf handle



- + Cable winding wheel



- + Scanner holder with Easy Clamp



- + Ventilator docking holder



+ Ambient light with on/off lighting control



+ Ambu-bag / resuscitator holder with EasyClamp



+ Single or double drawers



+ Flexible wipes canister holder with EasyClamp

+ and even more accessories...



This information is intended for an international audience outside the US.

This information is aimed exclusively at healthcare professionals or other professional audiences and is for informational purpose only, is not exhaustive and therefore should not be relied upon as a replacement of the Instructions for Use, service manual or medical advice. Getinge shall bear no responsibility or liability for any action or omission of any part based upon this material, and reliance is solely at the user's risk.

Any therapy, solution or product mentioned might not be available or allowed in your country. Information may not be copied or used, in whole or in part, without written permission by Getinge.

This document is intended to provide information to an international audience outside of the US.

Manufacturer · Maquet (Suzhou) Co., Ltd. · No.158 Fangzhou Road, Suzhou, P.R. China 215024 · +86 (0) 512 62839880

© 2025 Getinge | Getinge and **GETINGE** * are trademarks or registered trademarks of Getinge AB, its subsidiaries or affiliates | DMS-0010225 | All rights reserved.