

Maquet Moduevo PLG-II SKY Ceiling Supply Unit Improving mobile medical device connectivity







Improving medical device connectivity Maquet Moduevo PLG-II SKY

Maquet Moduevo PLG-II SKY supplies OR mobile medical devices (incl. anesthesia machine, heart-lung machine or endoscope devices on a trolley) with power, gas and transmitting signals.

Comfortable and ergonomic working conditions

- Easily maneuverable thanks to its external handle for proper positioning during procedures and safe placement in parking position, enabling easier conditions for cleaning
- Peace of mind: intuitive equipment thanks to a large cable capacity with sockets and outlets on both sides as well as in the bottom panel.
- An optional built-in ambient light can be added in the bottom panel to provide illumination during minimally invasive procedures or to support anesthetist with the visibility needed to focus on the patient.
- Once the surgery is over, cleaning is made easier as the pendant allows head clearance thanks to the parking position and at the desired height.

• Height-adjustable down to 400 mm, PLG-II SKY is an ergonomically friendly ceiling supply unit that can be easily adapted to different devices dimensions. Depending on medical staff heights, it also allows an easy access to sockets and outlets.

Suitable for various operating rooms

- Thanks to a reach radius of 2,975 mm max., PLG-II SKY offers the flexibility to position the anesthesia machine either at the left foot side in case of neurosurgery or at the right head side for other types of procedures.
- Customizable solution due to the modularity of the height and length of the suspension tube, it can therefore be adapted to different room layouts.
- In Hybrid ORs, it offers more freedom in terms of set-up as it can be installed out of the ceiling rails.

Flexible

- Depending on your future needs, you can always reserve space for new data sockets thanks to a high cable capacity.
- Compatibility with the Beam Plus option, allowing the costeffective addition of a monitor, an examination or even a surgical light without an extra anchor.
- Maquet Moduevo PLG-II SKY's spring arm requires no electricity, reducing energy consumption.

+ Comfortable and ergonomic working conditions


- + Suitable for various operating rooms
- + Flexible



Explore Maquet Moduevo PLG-II SKY

Beam Plus option

streamlines the connection of a monitor or light without an additional anchor

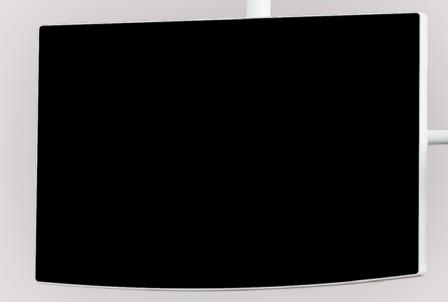

٠




without electricity

Intuitive external handle improves maneuverability





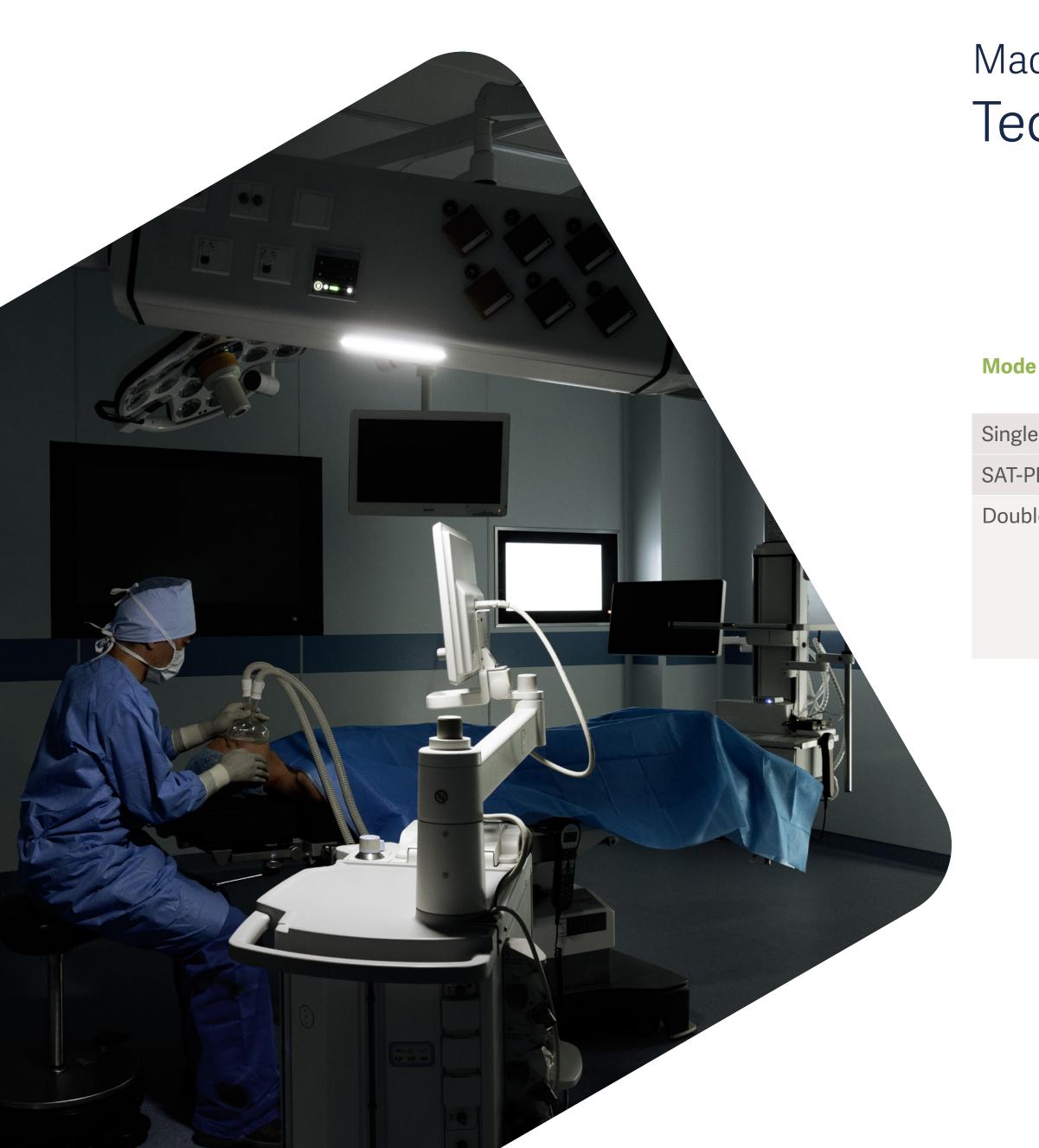

Maquet Moduevo PLG-II SKY Technical data

le	Configuration PLG-II SKY Distributor	Reach radius	Vertical movement
le beam mode	on anchor	1,475 mm/58"	
PLG mode	on SATELITE		
ble beam mode	with LL6-0	2,075 mm/82"	400 mm/16"
	with LL9-0	2,375 mm/94"	
	with LL12-0	2,675 mm/105"	
	with LL15-0	2,975 mm/117"	







With a firm belief that every person and community should have access to the best possible care, Getinge provides hospitals and life science institutions with products and solutions aiming to improve clinical results and optimize workflows. The offering includes products and solutions for intensive care, cardiovascular procedures, operating rooms, sterile reprocessing and life science. Getinge employs over 10,000 people worldwide and the products are sold in more than 135 countries

This document is intended to provide information to an international audience outside of the US.

Manufacturer: Maquet (Suzhou) Co., Ltd · No.158 Fangzhou Road · Suzhou Industrial Parl · 215024 Suzhou · China · +86 (0) 512 62839880

www.getinge.com

by DMS-00012226 · 11/2022