



OR Tables Accessories

2026 Catalog



Comprehensive expertise

for everyday clinical needs

A comprehensive range of high-quality accessories for the best possible patient positioning.

Getinge creates products in close collaboration with doctors and hospitals. This helps to optimize workflows, facilitate the work of the OR team, and improve the success of treatments for patients.

Getinge offers a comprehensive range of accessories for the best possible patient positioning. There are a variety of suitable accessories for different needs and types of surgery, be it infusion stands, or metal-

free and radiolucent multi-purpose plates. All accessories ensure safe and comfortable patient positioning and ergonomic working conditions for the surgical team. Nearly all accessories come equipped with special padding for optimum pressure distribution. Individual OR table configuration helps assist the needs of the patient and surgical team, thus ensuring a successful operation.



Content

00	Operating tables	06
Tested and proven around the globe: Maquet Operating Table Systems		
01	General accessories	16
Standard accessories that can be used in multiple surgical disciplines – from leg strap to trolley for accessories.		
02	Head rests	34
Head positioning accessories for every requirement, from the standard head rest to the radiolucent horseshoe head rest.		
03	Arm supports	45
Safe, gentle and flexible arm positioning using adjustable accessories and pressure-relieving pads.		
04	Back and lithotomy positioning	55
Cross-discipline solutions, from radiolucent multi-purpose plates and modular back plates through to seat plate extensions.		
05	Leg positioning	69
The wide range of leg positioning accessories offers optimum solutions for every discipline.		

About this catalog

This catalog gives an overview of the entire OR table accessory range from Getinge. The accessories are sorted according to practical use (positions) and surgical disciplines. Alongside the description of the accessory, there is also information about product variations plus a list of the corresponding OR tables. These OR tables are described in the following pages. The accessory matrix instantly shows which accessories may be mounted to a chosen OR table and whether any additional accessories are needed.

06	Orthopedics and traumatology	85
Innovative solutions for orthopedics and traumatology. Shape, material and pads are coordinated to ensure optimum X-ray results.		
07	Neurosurgery	110
Special accessories for safe and precise patient positioning designed to meet the special requirements of neurosurgery.		
08	Positioning aids	124
Gentle and preventive positioning of patients with the help of a variety of positioning aids. Simple to use, easy to clean and available in numerous sizes.		
09	Medical furniture	138
The Resist range of medical furniture consists of reliable and well proven accessories. Instrument tables, infusions stands, step stools and many more items are essential for the work of the OR team.		
10	Accessories matrix and index	148
The desired accessory can be located quickly and you can see, at a glance, which OR tables it is suitable for.		

Maquet Operating Table Systems

Tested and proven around the globe

Maquet Magnus

For extreme positions and optimum access

With inclination angles of up to 80° and tilt angles of up to 45° (may also be combined), Maquet Magnus enables extreme positions. These ensure optimum X-raying conditions and clear access to the operating field.

Thanks to a wide range of modules and accessories for all kinds of surgical intervention, Maquet Magnus is a highly versatile OR table that offers the greatest of flexibility in everyday clinical use.

The special features afforded by Magnus:

- Modular concept and components
- Flexible universal table top and carbon table tops for Hybrid ORs
- Backlit hand control with memory function
- Very low minimum height and very high maximum height
- 80° Trendelenburg low head / low feet
- 45° lateral tilt



Maquet Otesus

Explore new possibilities

Maquet Otesus features a number of key enhancements developed in collaboration with OR teams, while retaining many tried-and-trusted features of its proven predecessor – minimizing the need for training. The modular operating table system includes 4 columns, a smart universal hand control, a choice of user friendly transporters and a variety of table tops. Surgical staff can tailor the table to each surgery and patient by simply swapping a few modules.

Comfortable, hygienic cushions

Soft Integral Protective Comfortable (IPC) cushions offer excellent pressure relief performance that is comparable to Soft Foam Core (SFC) cushions. An inner core sandwich design connects 20 mm wear protection and 60 mm visco-elastic foam to a carrier plate, which is then covered with a spray skin. The hygienic, Velcro®-free pads are waterproof. Certain variants of the pads can be mechanically decontaminated and disinfected.

Patient transport between rooms

The Maquet Otesus Transporters allow surgical staff to maneuver patients with minimal physical effort. Even a motorized version – the transporter Ergodrive – is available. It simplifies the reconnection and disconnection of the table tops to the Maquet Otesus Column to streamline the workflow in the operating room.

All key operations at your fingertips

The Smart Control places all key table operations right at your fingertips. Bi-directional communication with the column provides real-time information on the current table configuration. The Smart Control continuously receives and displays the table set-up information, which saves time, enhances safety, and simplifies documentation.



Maquet Mobile Operating Tables

Versatile, tried and tested – real all-rounder

Maquet Betaclassic

High-quality and cost-efficient mobile operating table

A manually operated hydraulic power unit and ergonomically designed controls are standard on the robust Maquet Betaclassic. A few swift pumps of the pedal-operated selector enable precise adjustment of the table top position. The proven design ensures excellent table stability and unrestricted patient access for the surgical team.

Maquet Betaclassic provides dedicated solutions for military and disaster medicine, as it does not depend on any external power source. Thanks to the large castors, the OR table can be transported regardless of the terrain.



Maquet Lyra

A versatile table for all-around use

It's increasingly important to invest in reliable, multi-functional infrastructure. That's why Getinge developed the Maquet Lyra Mobile OR Table. Just one table can meet the needs of most surgical disciplines, and it is compatible with carbon-fiber accessories for excellent intra-operative diagnostics.

The Maquet Lyra Mobile OR Table is designed to work for almost all surgical disciplines, increasing the options for each OR. Plus, it's compatible with all existing Maquet Mobile OR Table accessories that you may already have, maximizing the return on your investment.

The Maquet Lyra Mobile OR Table has a robust design that ensures stability and safety for a broad range of patients. It has an overall load (patient plus accessory) of 360 kg (794 lbs) in central position, and 180 kg (397 lbs) without any restrictions – ideal for most patient populations.



Maquet Meera family

High performance, outstanding value

Hospitals of all sizes struggle with increasing cost pressures. With the Maquet Meera Mobile OR Table family, your hospital can cost-effectively deliver premium patient care.

Maquet Meera:

A universal mobile table for high-quality care

The mobile operating table Maquet Meera offers incredible flexibility to support a variety of procedures. The 310 mm motorized longitudinal shift delivers a large radiolucent area for gynecological or urological surgery, without an added seat plate extension. Maquet Meera is a stable OR table that meets the complex demands of bariatric surgery.



Maquet Meera ST:

A multifunctional table for various types of procedure

Maquet Meera ST offers the same benefits as the Maquet Meera, without the longitudinal shift. It ensures that your OR will benefit from a solid, functional table, without features that your surgeons do not require.



Maquet Meera CL:

Flexibility for many disciplines

Maquet Meera CL is built with the same stable components as the rest of the Meera range, but on a simplified base with a mechanical locking system. Specially designed castors enable smooth movement in all directions: intermediate operation, straight-line travel, and free-wheeling. Lateral movement capability is ideally suited for the needs of outpatient surgery centers.





—> Maquet Mobile Operating Tables

Maquet Corin

Meet your new teammate: intuitive, safe, smart and connected

Maquet Corin is an intelligent guided mobile OR table that can help to streamline processes in the daily routine of your surgical departments. Stay fully focused on patients knowing that your OR table is giving intuitive visual feedback that helps manage risks. Easily set up and position patients for optimized site access and ergonomics.

In short, Maquet Corin is your new teammate which helps you to:

- Facilitate teamwork and communication
- Ensure protection of people and equipment
- Prevent wrong-way errors with intuitive controls like feedback lights and acoustic signals
- Save time
- Stay flexible



GETINGE *
MAQUET | CORIN

OR TABLES ACCESSORIES

OR TABLES ACCESSORIES

Maquet Yuno II

Improves outcomes with advanced surgical workflows in orthopedics, traumatology and neurology

Maquet Yuno II ensures the positioning and access needed to achieve optimum results in open or minimally invasive interventions. The table can be configured for any discipline, increasing your options and maximizing the functionality of your OR. Easy handling ensures that Maquet Yuno II will simplify positioning for even the most complex techniques to help you focus on what's most important: your patients.

Minimally invasive surgery (MIS) can improve patient outcomes and save costs. Maquet Yuno II supports a variety of positions, bringing the benefits of MIS to

patients of various sizes. Radiolucent carbon-fiber accessories improve intraoperative diagnostics that are critical to achieving optimal outcomes.

When trauma patients arrive in your emergency department, there's no time to lose. Night or day, your OR must be ready. With the flexible Maquet Yuno II Mobile Operating Table, it's easy to quickly configure the table to meet the needs of the procedure and surgeon, even with limited staff on call.





01

Accessories for operating tables

General accessories

Each position is different, depending on the upcoming operation, the patient and the habits of the OR personnel.

Alongside the many special accessories, Maquet also offers standard OR table accessories that can be used regardless of the surgical discipline, from anesthesia screens and infusion stands through to body supports.

This chapter also lists general accessories, for example, leg straps for securing the patient to the OR table, side rail accessories, as well as transport and storage solutions.



1003.61A0
Radial setting clamp

- With self-clamping mechanism for the side rail
- Quick radial adjustment
- Simultaneous fixation of accessories, radial adjustment and clamping to the side rail

Technical information

Width x height x length:
70 x 140 x 120 mm

For attaching rods with
a diameter of 16 – 18 mm

EU version

Can be mounted to a side rail

Mechanical cleaning is possible



1003.23C0
Radial setting clamp

- For mounting accessories to the side rail
- Quick radial adjustment
- Simultaneous fixation of accessories and radial adjustment

Technical information

Width x height x length:
70 x 140 x 110 mm

For attaching rods with
a diameter of 16 – 18 mm

Can be mounted to a side rail

Mechanical cleaning is possible



1003.22C0
Standard clamp

- For mounting accessories to the side rail

Technical information

Width x height x length:
70 x 140 x 110 mm

For attaching rods with
a diameter of 16 – 18 mm

Can be mounted to a side rail

Mechanical cleaning is possible



1002.54A0
Anesthesia frame

- Individually adjustable
- Can be curved into any desired position

Technical information

Can be mounted to a side rail

2 standard clamps may be
ordered separately

Recommended combination with
• standard clamp 1003.22C0



1002.57A0
Anesthesia screen

- For holding OR drapes
- Insulated vertical rod protects the patient from grounding
- Height adjustable using standard clamp
- For placement and positioning of the patient's arm in conjunction with arm restraint cuff 1001.4600

Technical information

Height adjustment: 470 – 830 mm

Total height: 860 mm

Width of horizontal bar: 620 mm

Mounting to side rail
with standard clamp

Standard clamp may be
ordered separately

Recommended combination with
• anesthesia screen extension
1002.59A0



1002.59A0
Anesthesia screen extension

- Extension for anesthesia screen 1002.57A0
- Extension can be swivelled

Technical information

Length: 300 mm each

Including: set of 2

Only in conjunction with
• anesthesia screen 1002.57A0



1002.55A0
Tube holder

- For secure fixation of the breathing tube
- Flexible shaft enables infinite positioning
- May be used for a variety of tube diameters

Technical information

Can be mounted to a side rail

Standard clamp may be
ordered separately



1004.15A0
Lateral protection

- For securing the patient during transportation on the OR table
- Inner surface is padded with SFC
- Can be folded down

Technical information

Width x height x length:
45 x 485 x 725 mm

EU version

Can be mounted to a side rail

Integrated standard clamp



1005.88A0
Safety side rails

- For securing the patient during transportation on the OR table
- Can be folded down by 90°

Technical information

Width x height x length:
40 x 370 x 630 mm

Max. workload 20 kg in even
distribution corresponding to
a 200 kg patient weight

Integrated standard clamp

EU version

Material of device:
stainless steel, aluminum

Packing unit: 1 pair (2 pieces)

No mechanical cleaning

Can be mounted to a side rail



1005.89A0
Cover for safety side rails

- Cover for safety side rails
- Multiple use

Technical information

Width x height x length:
30 x 310 x 630 mm

Material of device: Mfoam

Packing unit: 1 piece

Can be mounted to the
safety side rail



1009.01C0
Infusion stand

- Ball joint enabling vertical alignment
- 4 metal hooks for holding
infusion containers

Technical information

Height adjustment: 870 – 1,550 mm

Max. load per hook: 1.25 kg

Can be mounted to a side rail

Integrated standard clamp



4009.3000
Infusion stand

- 2 metal hooks for holding
infusion containers

Technical information

Height adjustment: 875 – 1,300 mm

Max. load per hook: 2 kg

Standard clamp included

Can be set vertically in 15° angles
using the radial adjustment

Can be mounted to a side rail



1003.32C0
Mayo instrument stand

- Practical shelf for holding instruments
during surgery
- Height adjustable
- Flexible position of the shelf
- Angulation possible

Technical information

Shelf dimensions width x length:
260 x 450 mm

Height adjustment: 370 – 560 mm

Max. load: 3 kg

Integrated standard clamp with
a ball and socket joint for the rod

Can be mounted to a side rail



1002.98A0
Instrument table

- For holding instruments during surgery

Technical information

Interior dimensions of shelf surface
width x length: 215 x 380 mm

Max. load: 5 kg

Mounting to a pin mounting point of
an OR table or a seat plate extension



1003.44A0
Soiled instrument bowl

- For collecting OR waste
- Easy to move thanks to swivel mechanism

Technical information

Material: CrNi steel

Max. load: 5 kg

Mechanical cleaning is possible

Mounting to a pin mounting point



**1003.46A0
Trolley**

- Transportation and storage of a TUR rinsing set and a bucket

Technical information

Including: 1 CrNi steel bucket

Max. load: 20 kg

4 double swivel castors

Only in conjunction with

- TUR rinsing set 1003.45D0



**1003.45D0
TUR rinsing set**

- Hygienic intake and drainage of fluids
- Braked guide frame ensures safe use

Technical information

Comprises: 1 guide frame, 1 plastic basin, 1 drainage tube, 1 rinsing bowl, 1 sieve cap, 1 sieve, 2 tensioning bars

Mounting to OR table using guide rails

Recommended combination with

- trolley 1003.46A0



**1008.07E0
Clean catcher non sterile**

- Fully disposable uro-band and drainage bag eliminates cleaning dirty uro-band
- No more sagging bands

Technical information

Can be mounted to a side rail (EU and US)

Packing unit: 1 pack of 10 sets



**1008.07F0
Clean catcher sterile**

- Fully disposable uro-band and drainage bag eliminates cleaning dirty uro-band
- No more sagging bands

Technical information

Can be mounted to a side rail (EU and US)

Packing unit: 1 pack of 10 sets



**1002.40A0
Pin-joint arm for body supports**

- Very flexible positioning using 3 joints
- Quick locking using central handle screw
- High stability and safety

Technical information

For mounting of body supports 1002.11X0 and 1002.97A0

Can be mounted to a side rail

Clamp may be ordered separately



**1002.19C0
Fixture for body supports**

- Adjustment options for height and length
- Can be swivelled horizontally

Technical information

Width x height: 90 x 340 mm

Length: 280 – 380 mm

Standard clamp included

For mounting body supports 1002.11X0 and 1002.97A0

Can be mounted to a side rail



**1002.11A0
Back / buttocks support**

- Flexible positioning of the pad
- Pad can be rotated on the adapter

Technical information

Pad width x length: 120 x 170 mm

Material of pad: PUR

Combination with

- fixture for body supports 1002.19C0
- pin-joint arm 1002.40A0



**1002.11B0
Pubis / sternum support**

- Flexible positioning of the pad
- Pad can be rotated on the adapter

Technical information

Pad width x length: 85 x 85 mm

Material of pad: PUR

Combination with

- fixture for body supports 1002.19C0
- pin-joint arm 1002.40A0



**1002.11C0
Lateral support**

- Flexible positioning of the pad
- Pad can be rotated on the adapter
- Good pressure distribution thanks to slightly concave shape

Technical information

Pad width x length: 100 x 210 mm

Material of pad: PUR

Combination with

- fixture for body supports 1002.19C0
- pin-joint arm 1002.40A0



**1002.11D0
Supporting roll**

- Flexible positioning of the pad
- Pad can be rotated on the adapter

Technical information

Padded roll length: 200 mm

Ø 80 mm

Material of pad: PUR

Combination with

- fixture for body supports 1002.19C0
- pin-joint arm 1002.40A0



**1002.11E0
Buttocks support**

- For supporting the patient in lateral position
- Flexible positioning of the pad
- Pad can be rotated on the adapter

Technical information

Pad width x length: 240 x 150 mm

Material of pad: SFC

Combination with

- fixture for body supports 1002.19C0
- pin-joint arm 1002.40A0



**1002.97A0
Body support**

- For supporting the patient in lateral position
- Oval shape
- Pad can be rotated on the adapter

Technical information

Pad width x length: 120 x 140 mm (at widest point)

Thickness of pad: 35 mm

Material of pad: SFC

Recommended combination with

- fixture for body supports 1002.19C0
- pin-joint arm 1002.40A0



**1002.21C0
Lateral support**

- For lateral securing of the patient
- Good pressure distribution thanks to slightly concave shape

Technical information

Pad width x length: 100 x 215 mm

Material of pad: PUR

Overall height: 360 mm

Can be mounted to a side rail

Integrated standard clamp



**1002.36A0
Lateral support**

- For supporting the thorax in Beach Chair position
- Good pressure distribution thanks to slightly concave shape

Technical information

Pad width x length: 100 x 215 mm

Including: set of 2 (right and left)

Material of pad: PUR

Can be mounted to a side rail

Integrated standard clamp



**1002.17A0
Double lateral support**

- Two individually adjustable pads for stable support, especially in lateral position

Technical information

Pad width x length: 85 x 85 mm

Including: standard clamp

Material of pad: PUR

Can be mounted to a side rail



**1002.09C0
Shoulder support**

- Support of the shoulder area, especially in Trendelenburg position
- Individually adjustable to patient's anatomy

Technical information

Length of vertical bar: 190 mm

Size of pad width x length: 100 x 210 mm

Comprises: 2 shoulder supports and 2 standard clamps

Material of pad: PUR

Can be mounted to a side rail

Available versions:

- 1002.09C0
Length of horizontal bar: 190 mm



**1001.75A0
Table width extension, long**

- Good pressure distribution thanks to SFC padding
- With side rail for mounting of additional accessories

Technical information

Width x length: 100 x 300 mm

Including: set of 2

Integrated fixation clamp

EU side rail

Can be mounted to a side rail



1001.76A0
Table width extension, short

- Good pressure distribution thanks to SFC padding
- With side rail for mounting of additional accessories

Technical information

- Width x length: 100 x 200 mm
- Including: set of 2
- Integrated fixation clamp
- EU side rail
- Can be mounted to a side rail



1002.94A0
X-ray cassette holder

- For horizontal positioning and securing of an X-ray cassette underneath the OR table top
- For X-ray cassettes with a short side longer than 180 mm

Technical information

- Diameter of support plate: 200 mm
- Cassette holder height x length: 402 x 497 mm
- Can be mounted to a side rail
- Integrated standard clamp



1002.93A0
X-ray cassette slide

- For accomodating an X-ray cassette
- Width can be adjusted to dimensions of the cassette
- Push rod can be folded in

Technical information

- Variable width: 360, 390, 400 mm
- Length: 506 mm
- Suitable for use with**
 - 1118.01, 1118.06, 1118.10
 - 1131, 1132, 1133
 - 1160.10, 1160.30
 - 1433
 - 7000, 7200, 7100, 7700
 - 1150 table tops except 1150.16/20



1008.11A0
Radiation protection

- Flexible X-ray protective curtain
- Lead equivalent value: 0.5 mm
- 2 segments for the lower body mounted on a pivotal arm
- Integrated side rail for mounting accessories
- 3 segments for the upper body
- Upper body protection may be adapted individually

Technical information

- Weight: 16.9 kg
- Suitable for EU and US side rails
- No mechanical cleaning
- Can be mounted to a side rail
- Integrated standard clamp
- Kenex catalog no. 312/DS-0412/6



1009.69A0
Sensor Drive

- Wireless control unit for OR tables with Autodrive function
- Precise speed regulation using the throttle grip
- Easy movement of the OR table
- Forward and backward motion is possible

Technical information

- Maximum speed: 3 km/h
- Can be mounted to a side rail
- Recommended combination with**
 - charging units 1009.70A0 or 1009.71A0/B0
- Suitable for use with**
 - tables with Autodrive function
 - 1133.22B5, 7200.01B2, 1433.02B0
 - 7700, 7700.01B0, B2



1003.2500
Push bars

- Comfortable moving of operating tables
- Handle with plastic cover

Technical information

- Width x height x length: 50 x 200 x 400 mm
- Including: set of 2
- Can be mounted to a side rail
- Integrated standard clamp



1180.49A0
Side rail adapter pair

- To mount the catheter holder and other side rail accessories

Technical information

- EU side rail
- Maximum load per adapter: 10 kg
- Suitable for use with**
 - 1180.16X0 CF table tops



1180.47A0
Catheter holder

- Catheter holder for secure storage of catheters during surgical interventions
- The side guard prevents catheters from falling down
- Can easily be mounted to a side rail adapter or an OR table top

Technical information

- Width x height x length: 600 x 190 x 1,215 mm
- Maximum load: 6 kg
- Combination with**
 - the pair of side rails 1180.49A0
- Suitable for use with**
 - 1180.16X0 CF table tops
 - 1160.16AX CF table tops



1005.40A0
Locking device for side rails

- For stabilisation of the seat plate, e.g. for neurosurgical interventions
- Side rail accessories can be mounted as usual

Technical information

- Length: 400 mm
- EU side rail
- Can be mounted to a side rail
- 2 standard clamps integrated



1004.91A0
Side rail extension

- For mounting side rail accessories

Technical information

- Length: 490 mm
- EU side rail
- Can be mounted to a side rail
- 1 standard clamp integrated



1001.54A0
Side rail extension

- Side rail extension with two integrated standard clamps

Technical information

- Length: 690 mm
- EU side rail
- Can be mounted to a side rail
- 2 standard clamps integrated
- Recommended combination with**
 - OR table top 1150.20
- Can be combined with**
 - 1008.11A0



1004.91B0
Side rail extension

- For mounting side rail accessories

Technical information

- Length: 200 mm
- Material: CrNi steel
- EU side rail
- 1 standard clamp integrated
- May only be mounted to the end of the side rail



1005.22A0
Extension rail

- Extremely long side rail for mounting accessories to the carbon-fiber table top

Technical information

- Length: 1,700 mm
- EU side rail
- 2 standard clamps integrated
- Only for use in combination with**
 - accessory adapters 1150.70A0 or 1150.72A0
 - 1150.16A0



1001.56A0
Body strap

- Body strap with metal closure
- Continuously adjustable
- 2 half straps, one long, one short

Technical information

- Width x length: 120 x 1,500 mm
- Material of strap: polyester fabric
- Strap is washable (95°C)
- Variant 1001.56A1 for Magnus
- Can be mounted to a side rail with a metal strap holder



1001.56B0
Body strap

- Body strap with Velcro® closure
- Continuously adjustable
- 2 half straps, one long, one short

Technical information

- Width x length: 120 x 1,500 mm
- Strap is washable (60°)
- Variant 1001.56B1 for Magnus
- Can be mounted to a side rail with a metal strap holder



1001.59A0
Body strap

- Infinite adjustment
- Secure closure with metal clamp

Technical information

- Width x length: 120 x 1,500 mm
- Material of strap: polyester fabric
- Strap is washable (60°)
- Can be mounted to a side rail with a metal strap holder



1001.59B0
Body strap

- Infinite adjustment
- Secure closure with hook and loop fastener

Technical information

- Width x length: 120 x 1,500 mm
- Material of strap: polyester fabric
- Strap is washable (60°)
- Can be mounted to a side rail with a metal strap holder



1008.12A0
Patient safety strap, medium

- Body strap with metal closure
- Continuously adjustable from 450 – 1,140 mm
- Gel padded center section

Technical information

- Width x height x length: 127 x 12 x 2,430 mm
- Maximum patient weight: 227 kg
- Material of strap: nylon
- No mechanical cleaning
- Can be mounted to a side rail
- 2 clamps integrated
- Suitable for use with all OR tables



1001.58A0
Body strap

- Body strap with Velcro® fastener
- Continuously adjustable

Technical information

- Width x length: 120 x 2,200 mm
- Material of strap: polyester fabric
- Strap is washable (60°)
- Suitable for use with all OR tables



1002.18A0
Support bar

- Height-adjustable
- Stable stand
- Adaptation to the floor is possible via integrated ball joint in the base stand

Technical information

- Adjustable height: 740 – 1,205 mm
- For mounting in round tube interface (Ø 18 mm), e.g. radial setting clamp
- Recommended combination with**
 - pediatric operating table top 1005.78C0



1004.07A0
Column protection device

- Protection against damage at the column cover
- Upper part removable
- Can be mounted to the basis segment of a mobile OR table with three point stand

Technical information

- Width x height x length: 70 x 132 x 440 mm
- Material of device: stainless steel
- Packing unit: 1 set of 4 items
- Mechanical cleaning is possible

Can be combined with

- Meera 7100, 7200
- Corin 7700

Suitable for use with

- Corin 7700.01XX
- Meera 7100.01XX, 7200.01XX



1132.40A0
Accessory trolley

- Space saving storage of accessories for mobile operating tables
- Mounting points for head, back and leg plates
- 1 side rail for mounting side rail accessories
- Convenient transport and storage of accessories

Technical information

- Width x height x depth: 820 x 1,550 x 805 mm
- Up to 3 wire baskets may be mounted
- 4 double swivel castors, 2 with brakes
- Mechanical cleaning is possible



1003.48A0
Accessory trolley

- Space-saving storage of accessories for OR tables
- 2 side rails provided
- Mounting options for leg holders and anesthesia screen

Technical information

- Width x height x depth: 900 x 1,450 x 576 mm
- Including: 1 wire basket 0802.0783
- Up to 4 wire baskets can be mounted
- 4 double swivel castors, 2 with brakes
- Mechanical cleaning is possible



0802.0783
Wire basket

- For accommodating and storing small accessories

Technical information

Width x height x length:
275 x 135 x 835 mm

Mechanical cleaning is possible

Combination with

- accessory trolleys 1003.48A0, 1132.40A0, 1007.06A0



1007.06A0
Accessory trolley

- Space saving storage of accessories for the Magnus OR Table

Technical information

Width x height x length:
700 x 1,390 x 900 mm

1 metal basket mounted

2 side rails for mounting accessories

4 double swivel castors,
2 with brakes

Recommended combination with

- wire basket 0802.0783

02 Accessories for operating tables

Head rests



The head rest plays an important role in the anatomically correct positioning of the patient. For this purpose, Maquet offers a wide variety of head positioning accessories: from head rest adapters and simple head rests through to radio-lucent horseshoe head rests.

All standard head rests and horseshoe head rests have special pads. SFC materials and gel pads ensure optimum distribution of pressure. This enables the OR team to concentrate on the operation with the patient comfortably and safely positioned. Easy Click interfaces enable quick and easy mounting of the head positioning accessories. The requirements of anesthesia and surgical departments are equally met: The anesthesia department requires easy operation with lateral side rails for accessories, whereas surgeons require precise adjustment options.



1130.54B0
Connection bracket


- For mounting head positioning accessories, e.g. for neurosurgery

Technical information

Mounting to a small pin mounting point

Square mount for additional head positioning accessories, e.g. fixture 1002.65A0

Round tube mount for DORO® Adjustable Base Unit 1005.50A0

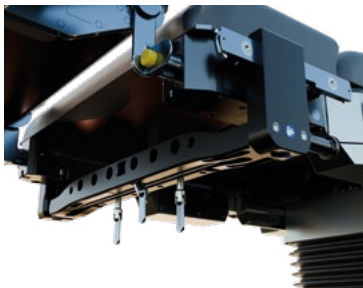

1004.05A0
Connection bracket

- For EU side rail
- For head positioning accessories
- Fast and easy installation from one side via "Quick Lock" by a single person
- Width-adjustable:
 - Versatile for various head options
- Double-lock prevents accidental release
- Allows installation even when the patient is already intubated

Technical information

Can be mounted to a side rail

Weight: 4.5 kg

Maximum patient weight: 250 kg

Compatible with:
DORO® Adjustable Base Unit - skull clamp holders from Integra Mayfield® and Mizuho Sugita® Head frame system as well as for accessories with square rod mounts



1180.36A0
Connection bracket

- Enables attachment of a variety of head positioning accessories
- Easy mounting with Easy Click technology

Technical information

Interface: Magnus basic unit

Round tube mount for DORO® Adjustable Base Unit 1005.50A0

Square mount for head positioning accessories, e.g. fixture 1002.65A0

Suitable for use with

- motorized joint module 1180.11A0
- 1180.32A0



1002.07A0
Head rest adapter

- For mounting standard head rests and horseshoe head rests
- Height and inclination are adjustable
- Securing mechanism to prevent unintentional pull-out

Technical information

Mounting to a square interface

Recommended combination with

- connection bracket, e.g. 1130.54B0



1007.16A0
Head rest adapter

- For mounting standard head rests and horseshoe head rests
- Height and inclination can be adjusted
- Can be rotated
- Securing mechanism to prevent unintentional pull-out

Technical information

Length of height-adjustable bar: 340 mm

Length of side rail: 200 mm

Mounting to square mount

Combination with

- connection bracket 1130.54B0



1002.10A0
Adjustable head rest adapter

- For mounting standard head rests and horseshoe head rests
- High stability and safety
- Flexible and easy adjustment by using 3 joints
- Quick locking with only one central handle screw

Technical information

Maximum patient weight: 250 kg

Width x length:
200 x 438 mm extended

Safety lever to prevent unintentional pull-out from interface

Mounting to OR tables, table tops and accessories with a round pipe mounting point

1002.66A0
Horseshoe head rest with one-handed adjustment

- For placement and positioning of the patient's head
- Can be used in prone and supine position
- Quick one-handed adjustment
- Infinite adjustment
- SFC pads

Technical information

Including: fastening strap and pad

Mounting to a square mounting point

Width x height x length:
320 x 120 x 460 mm

Safety lever to prevent unintentional pull-out from interface

Locking lever can be slid in and out

Combination with

- head rest adapter or trapezoidal shaped head plate (e.g. 1160.35AC)



1002.65A0
Fixture

- Connection element for head rests, e.g. horseshoe or head plate to the OR table
- Side rail for mounting of head rests
- Adjustment by using 3 individually adjustable joints

Technical information

Length: 430 mm (extended)

Width of side rail: 250 mm

Mounting to a square mounting point

Safety lever to prevent unintentional pull-out from interface

Combination with

- connection brackets, e.g. 1130.54B0

1002.81D0
Head rest

- Good access to the patient's head thanks to the extra-slim design
- Infinite inclination adjustment
- Controlled adjustment via braking mechanism

Technical information

Width x length: 210 x 260 mm

Material of pad: PUR

Adjustment angle: ± 50°

Height adjustment: 80 mm

Mounting to a square mounting point



1002.73A0
Flat head rest

- Radiolucent with the exception of the mounting point region
- Support plate can be shifted longitudinally and is fixed with 2 handle screws

Technical information

Width x length: 230 x 290 mm

Material of pad: SFC

Mounting to a fixture with side rail mounting point, e.g. 1002.65A0

1130.81A0
Head rest adapter

- For accommodating head rests and connection brackets to the lower back plate

Technical information

1 pair (right, left)

Mounting to a pin mounting point

Suitable for use with

- 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



1002.72A0
Horseshoe head rest – two-section

- Width of the two sections can be adjusted
- Can be used in prone and supine positions
- 2 seamless SFC pads

Technical information

Width x height x length:
230 x 120 x 260 mm
(total width of two sections)

Including: 1 fastening strap incl. forehead pad

Material of pad: SFC

Mounting to fixture (e.g. 1002.65A0) or head rest adapter with a side rail



1132.81A0
Head rest adapter

- For accommodating head rests and connection brackets on the lower back plate

Technical information

1 pair (right, left)

Mounting to a pin mounting point

Suitable for use with

- 1118
- 1132



1130.64G0
Head rest

- For positioning of the patient's head
- Quick adjustment using a trigger bar
- Equipped with X-ray cassette holder

Technical information

Adjustment angle: + 47° / - 61°

Width x height x length:
508 x 160 x 250 mm
(without side rails)

Maximum patient weight: 380 kg

Material of pad: SFC

Mounting to a small
pin mounting point

EU side rail

Recommended combination with

- extension plate 1131.31BC
- head rest adapter 1130.81A0 or 1132.81A0

Available versions:

- 1130.64C0 with SFC padding, without side rail
- 1130.64G0 with SFC padding, with side rail



1130.67B0
Head rest

- Easy adjustment using pneumatic support
- Adjustment can be set via a release lever
- Equipped with X-ray cassette holder

Technical information

Adjustment angle: + 30° / - 45°

Width x height x length:
540 x 160 x 295 mm

Maximum patient weight: 380 kg

Material of pad: SFC

Mounting to a small
pin mounting point

EU side rail

Recommended combination with

- extension plate 1131.31BC
- head rest adapter 1130.81A0 or 1132.81A0



1130.53B0
Head rest

- Controlled adjustment via braking mechanism
- Additional adjustment options thanks to double joint – head plate can be inclined and raised
- Enabling comfortable lateral positioning
- Equipped with X-ray cassette holder

Technical information

Adjustment angle: + 33° / - 90°

Width x height x length:
508 x 160 x 250 mm

Maximum patient weight: 135 kg

Mounting to a small
pin mounting point

Recommended combination with

- extension plate 1131.31BC
- head rest adapter 1130.81A0 or 1132.81A0

Available versions:

- 1130.53B0 with SFC padding
- 1130.53D0 with PUR padding



1160.64AC
Head rest

- Head rest with inclination
- Colour-coded adjustment levers
- IPC pad with integrated base plate – Velcro®-free fixation – available
- Side rail mounted on each outside edge
- Fitted with cassette holder for X-ray cassette slide (1002.93A0)
- Easy mounting with Easy Click technology

Technical information

Adjustment angle: + 84° / - 60°

Width x height x length:
590 x 139 x 300 mm

Maximum patient weight: 454 kg

Mounting to a pin mounting point

1160.64AC suitable for
mechanical cleaning

EU side rail

Suitable for use with

- 1160.10/30
- mobile OR tables 1133.22, 1433, 7000, 7100, 7200, 7700

Available versions:

- 1160.64AC with IPC padding
- 1160.64BC with SFC padding



- 1160.53AC**
Head rest
- Head rest can be inclined and raised
 - Colour-coded adjustment levers
 - IPC pad with integrated base plate – Velcro®-free fixation – available
 - Side rails mounted on each outside edge
 - Controlled adjustment via braking mechanism
 - Fitted with cassette holder for X-ray cassette slide (1002.93A0)
 - Easy mounting with Easy Click technology

Technical information

Large adjustment angle: +84° / -60°

Width x height x length:
590 x 139 x 300 mm

Maximum patient weight: 380 kg

Mounting to a pin mounting point

1160.53AC suitable for mechanical cleaning

EU side rail

Suitable for use with

- 1160.10/30
- mobile OR tables, 1133.22, 1433, 7000, 7100, 7200, 7700

Available versions:

- 1160.53AC with IPC padding
- 1160.53BC with SFC padding



- 1180.53A0**
Head rest
- Additional adjustment options thanks to double joint – head plate can be inclined and raised
 - Anatomical positioning of the patient's head in lateral and supine positions
 - Easy mounting with Easy Click technology

Technical information

Large adjustment angle: +117° / -43°

Width x length: 550 x 350 mm

Maximum patient weight: 380 kg

Interface: Magnus basic unit

Material of pad: SFC

EU side rail

Recommended combination with

- extension plate 1180.32A0

Suitable for use with:

- 1180.10A0



- 1180.50A0**
Head rest
- Quick adjustment of inclination via release lever below the plate
 - Controlled adjustment via braking mechanism
 - Easy mounting with Easy Click technology

Technical information

Large adjustment angle: +36° / -60°

Width x length:
580 x 360 mm

Maximum patient weight: 380 kg

Interface: Magnus basic unit

Material of pad: SFC

EU side rail surrounding the plate

Recommended combination with

- back plate 1180.31A0

Suitable for use with:

- 1180.10A0



- 1180.46A0**
Extension plate
- For extending the table top
 - Metal-free radiolucent

Technical information

Width x length: 460 x 390 mm

Maximum patient weight: 250 kg

Material of pad: SFC

Mounting to screw-fitted mounting point

Recommended combination with

- carbon-fiber plate 1180.45A0



- 1002.15A0**
Head support for shoulder operations
- Side elements can be folded down
 - Can be used for a variety of head circumferences
 - Stable fixation of the head

Technical information

Width x height x length:
200 x 156 x 332 mm

Maximum patient weight: 250 kg

Including: 2 fixture straps

Soft pads are fixed with Velcro® and can be replaced

Mounting to the side rail mounting point

Recommended combination with

- shoulder plates



1150.71C0
Head rest

- Trapezoidal-shaped head rest
- Metal-free radiolucent
- Lightweight

Technical information

Width x height x length:
340 x 65 x 200 mm

Maximum patient weight: 245 kg

Material of pad: SFC

Mounting to a dovetail
mounting point

Recommended combination with

- carbon-fiber table top 1150.16B0

Available versions:

- 1150.71B0 with SFC padding,
combination with 1150.16B0;
width x length: 460 x 200 mm
- 1150.71C0,
combination with back plates
1150.37B0, 1150.45BC;
width x length: 340 x 200 mm

03 Accessories for operating tables

Arm supports



Correct support for the patient’s arm plays an important part in a successful procedure. A poorly positioned arm may lead to nerve damage for the patient. Even in short interventions, it may result in pressure points on the skin.

Maquet accessories with SFC padding minimise this risk. In addition, the exact adjustment mechanisms ensure that the supports can be adapted to suit both patient and OR team so that the patient's arm can be securely and comfortably positioned.

These arm-support accessories meet the requirements of anesthesia and surgical departments alike. The anesthesia team benefits from quick and easy adjustment of the arm, while the work of the OR team is facilitated by using the arm protector 1002.25A0. This arm protector secures the patient’s arm and, at the same time, offers the surgeons something to lean on and thus relieve the strain on their spine.



**1001.44E0
Arm support**

- Individually adjustable
- Can be rotated and adjusted in length and height
- Hollow shape of pad supporting secure positioning
- Secure fixation of the arm with 2 fixture straps

Technical information
Length of arm rest: 300 mm
Including: 2 fixture straps
Spare part: fastening strap 3100.7614
Can be mounted to a side rail
Integrated standard clamp
Available versions:
• 1001.44E0 with SFC padding



**1001.44F0
Arm support**

- Individually adjustable
- Can be rotated and adjusted in length and height
- Hollow shape of pad supporting secure positioning
- Secure fixation of the arm with 2 fixture straps

Technical information
Length of arm rest: 450 mm
Including: 2 fixture straps
Spare part: fastening strap 3100.7614
Can be mounted to a side rail
Integrated standard clamp
Available versions:
• 1001.44D0 with PUR padding
• 1001.44F0 with SFC padding



**1008.17A0
KTEK Armboard Swivel, large**

- Individually adjustable
- Can be rotated and adjusted in length and height
- Foam pad with integrated padded fixation part for the arm

Technical information
Width x height x length: 317 x 16 x 450 mm
Weight: 2.8 kg
Maximum patient weight: 227 kg
For EU and US side rails
Material of device: stainless steel, aluminum
Material of pad: Polyurethane foam
Comprises: 1 piece incl. pad with fixation part
No mechanical cleaning
Can be mounted to a side rail
Clamps to be ordered seperately:
• 1003.22C0 standard clamp
• 1003.23C0 radial setting clamp
Available versions:
• 1008.17A0
• 1008.17B0



1008.17B0
KTEK Armboard Swivel, small

- Individually adjustable
- Can be rotated and adjusted in length and height
- Foam pad with integrated padded fixation part for the arm

Technical information
Width x height x length: 165 x 16 x 450 mm
Weight: 2.7 kg
Maximum patient weight: 227 kg
For EU and US side rails
Material of device: stainless steel, aluminum
Material of pad: Polyurethane foam
Comprises: 1 piece incl. pad with fixation part
No mechanical cleaning
Can be mounted to a side rail
Clamps to be ordered seperately: • 1003.22C0 standard clamp • 1003.23C0 radial setting clamp
Available versions: • 1008.17A0

SU01.04990001
Swivel arm board

- Pressure distribution thanks to SFC padding
- Can be swivelled horizontally gradually



Technical information
Width x height x length: 145 x 575 x 210/190 mm
Weight: 2.8 kg
Maximum patient weight: 135 kg
Swivel range: 0 – 180°, in 22.5° steps
Including: standard clamp
For EU side rails
Can be mounted to a side rail
Available versions: • SU01.04990001 with padding height 80 mm • SU01.04990002 with padding height 60 mm



1001.6000
Arm support

- One-hand operation
- Adjustable both horizontally and vertically
- Self-locking adjustment lever
- Lightweight accessory

Technical information
Width x height x length: 180 x 200 x 600 mm
Weight: 3 kg
Maximum patient weight: 135 kg
Adjustment angle: +30° / -40°
Including: 1 standard clamp and 2 fixture straps
With SFC pad
Spare part: fastening strap 3101.7614
Can be mounted to a side rail

1008.15A0
KTEK Trigger Armboard

- One-hand operation
- Freely adjustable
- Integrated snap clamp can be mounted on side rails in any position
- Radiolucent



Technical information
Width x height x length: 140 x 75 (pad) x 660 mm
Weight: 1.8 kg
Maximum patient weight: 227 kg
For EU side rails
Material of device: stainless steel, aluminum, carbon steel (?)
Material of pad: soft foam
Comprises: 1 armboard, 1 pad, 2 soft armboard straps
Spare parts: soft arm board straps (Kyra prod. no. 8310)
No mechanical cleaning
Can be mounted to a side rail
Integrated snap clamp



1002.30A0
Trimano Fortis

- Support arm can be used in various surgical disciplines
- Simple to operate with the Click and Move system
- Flexible positioning also during the operation
- Easy mounting of accessories
- Immediately ready for use without power supply or additional compressed air connection

Technical information

Width x height x length:
110x 390x 900 mm

Weight: 6.4 kg

Maximum patient weight: 250 kg

Sterile cover sets available

Can be mounted to a side rail

Combination with

- adapter Trimano Fortis 1002.31A0
- arm support 1002.32A0



1002.32A0
Arm support

- Accessory for Trimano Fortis
- Easy mounting with Easy Click technology
- Secure positioning of the arm

Technical information

Width x height x length:
125 x 90 x 310 mm

Weight: 0.7 kg

Maximum patient weight: 250 kg

Including: 1 fastening strap

Material of pad: SFC

Variant:

arm rest 1002.43C0 for Trimano 3D



1002.49A0
Arm rest

- Very flexible adjustment using 3 joints
- Quick locking using central handle screw
- Colour coding identifies the closure of the handle screw
- Rounded pad for protecting either the bend of the elbow or the axilla

Technical information

Dimensions of arm board pad
width x length: 130 x 520 mm

Extended length of arm rest: 780 mm

Weight: 4 kg

Maximum patient weight: 250 kg

Including: 2 fastening straps

Adjustable arm: 780 mm

Can be mounted to a side rail

Clamp may be ordered separately:

- 1003.22C0
- 1003.23C0
- 1003.61A0



1005.30A0
Fisso arm support lateral

- Flexible positioning by using 3 articulations
- Central locking mechanism for quick fixation
- Pad shaped ergonomically
- Pad mounted directly to carrier plate

Technical information

Height adjustment joint arm: 460 mm

Length of the pad: 450 mm

Maximum patient weight: 10 kg

Fixation clamp included



1001.4600
Arm restraint cuff

- For fixation of the patient's lower arm to an anesthesia screen
- Enables raised position of the arm
- Secure fixation by hook and loop fastener

Technical information

Width x length: 120 x 350 mm

Weight: 0.5 kg

Maximum patient weight: 135 kg

Strap width x length: 120 x 350 mm

Strap is washable (60 °C)

Only in conjunction with

- anesthesia screen 1002.57A0



1001.63A0
Multi-purpose plate for arm surgery

- Metal-free and radiolucent in the middle section
- Horizontal position of the hand operating table can also be maintained with Trendelenburg position
- Support bar can be mounted on both sides of the table end
- Additional accessories may be mounted on the side rail
- Radiolucency from finger tip to finger tip when using two hand operating tables in combination with a carbon-fiber table top and adapter 1001.64X0

Technical information

Width x height x length:
380 x 25 x 845 mm

Weight (plate only): 5 kg

Maximum patient weight: 454 kg

Height adjustment of support bar:
700 – 1,160 mm

Comprises: multi-purpose plate,
1 support bar 3115.1069, 1 T-adapter

EU side rail

Mounting to a side rail
with radial setting clamp

Clamp may be ordered separately:

- 1003.22C0
- 1003.23C0
- 1003.61A0

Recommended combination with

- adapter for carbon-fiber table top 1001.64X0



1005.95A0
Carbon-fiber arm surgery table

- Carbon-fiber carrier plate
- Soft foam pad
- Center piece is radiolucent
- CF support bar can be mounted on both sides of the table end

Technical information

- Height adjustment: 720 – 1,150 mm
- Width x length: 525 x 820 mm
- Weight: 5.6 kg
- Maximum patient weight: 270 kg
- EU side rail
- Material of device: carbon fiber, aluminum, stainless steel
- Comprises: 1 multi purpose plate incl. pad, 1 support bar
- No mechanical cleaning
- Can be mounted to a side rail
- 2 fixation clamps included



1001.64X0
Adapter for CF table top

- For mounting the hand operating table 1001.63A0 to a carbon-fiber multi-purpose plate or OR table top
- Radiolucent

Technical information

- Including: 1 side rail
- EU side rail

Available versions:

- 1001.64A0 for CF plate 1150.45BC, 1001.64B0 for CF plate 1180.45A0 and OR table tops 1180.16A1, A2, A5, 1160.16 A1, A0
- 1001.64C0 for CF plate 1132.45AC and OR table tops 1150.16B0, 1180.16A6



1001.42A0
Arm support

- Radiolucent through a full 360°
- Easy to slip under the SFC padding
- Radiolucency from spine to hand is possible, e.g. by combining with carbon-fiber multi-purpose plate 1150.45BC

Technical information

- Posturing plate width x length: 500 x 1,000 mm
- SFC pad width x length: 250 x 670 mm
- Comprises: 1 arm board, 1 SFC pad, 2 fastening straps
- Spare part: fastening strap 3101.8434



1001.40B0
Arm protector, curved

- Easy to slip under the SFC padding
- The patient's arm is protected safely
- The curved shape gives extra space for the arm

Technical information

- Width x height x length: 275 x 156 x 400 mm
- Weight: 1 kg
- Material: Plexiglas, transparent
- No mechanical cleaning

Variant:

- 1001.40A0

Compatible with

- all OR table tops with SFC padding with exception of 1150.16B0, 1180.10, 1150.13B0, 1180.13A0



1001.40A0
Arm protector

- Easy to slip under the SFC padding
- The patient's arm is protected safely

Technical information

- Width x height x length: 225 x 210 x 400 mm
- Material: Plexiglas, transparent
- No mechanical cleaning

Compatible with

- all OR table tops with SFC padding with exception of 1150.16B0, 1180.10, 1150.13B0, 1180.13A0



1002.25A0
Arm protector

- Arm protector for positioning and securing the patient's arm
- Sturdy stainless steel guard protecting the arm
- Inlay: SFC padding
- Relieves strain on the back muscles by allowing the surgeon to rest against it

Technical information


Width x height x length:
250 x 190 x 450 mm

Including: 2 straps to secure the arm

Material of pad: SFC

Can be mounted to a side rail

Integrated standard clamp

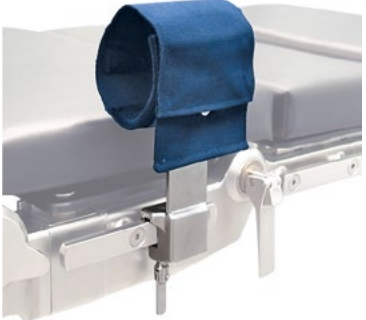

1002.24C0
Arm restraint cuff

- For securing the arm
- Fixation of arm strap with Velcro® fastener
- Strap can be removed and is washable

Technical information

Width x height x length:
70 x 50 x 410 mm

Weight (plate only): 0.7 kg

Spare part: strap 3101.0234

Can be mounted to a side rail


Integrated standard clamp

Strap is washable (60 °C)

04

Cross-discipline accessories

Back and lithotomy positioning

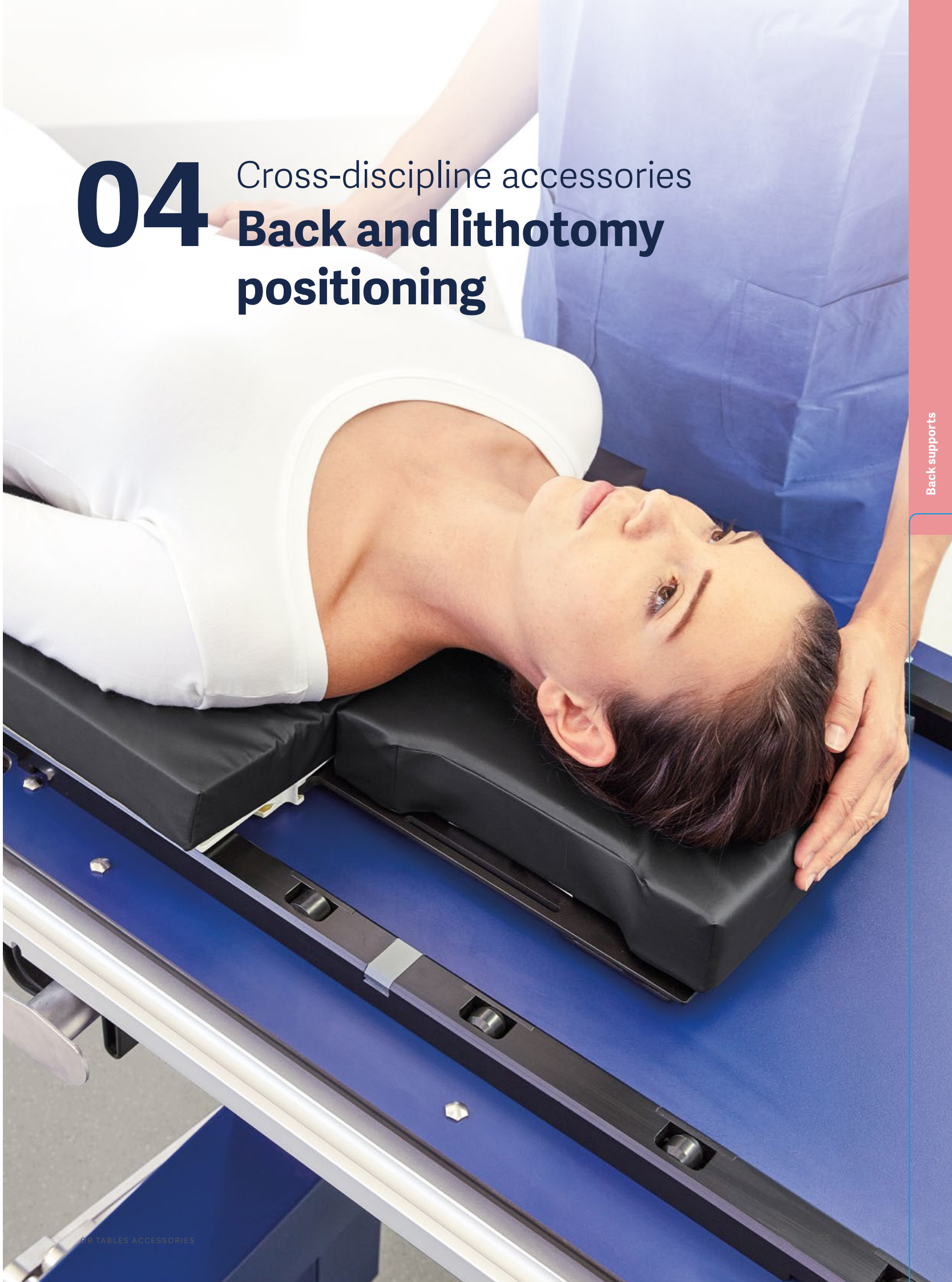

The range of back plates supplied by Maquet is extremely versatile. They include from radiolucent multi-purpose plates, modular back plates for shoulder operations, extension plates for very large patients, and trapezoidal shaped plates that improve ergonomic access to patients. Seat plate extensions, for example, for laparoscopy or urological interventions, complete the range.

All accessories for supine positions have something in common: high quality and durability. As surgical procedures have changed over the years, appropriate OR tables and accessories have been developed in line with the new requirements. They are designed to promote optimum workflow by way of their high modularity, adaptability and short changeover times.



- 1160.45XC**
Carbon-fiber insertable plate
- Radiolucent through a full 360°
 - Very long X-ray window depending on the version and attached head rests

Technical information
Mounting to interface for insertion of 1160 table tops and mobile OR table Corin
Material: carbon fiber
With SFC pad
Recommended combination with • accessory trolley 1003.86A0
Combination with • adapter for CF plates 1001.64B0 • accessory adapter 1180.72A0
Suitable for use with • 1160.10/30 • 7.700.01
Available versions: • Insertable CF plate 1160.45AC without head rest interface, length: 1,850 mm • Insertable CF plate 1160.45BC with head rest interfaces for screw-fitted and push-fit head rests, length: 1,520 mm



- 1433.45BC**
Carbon-fiber insertable plate
- Radiolucent through a full 360°
 - Very long X-ray window depending on the version and attached head rests

Technical information
Interface for insertion on 1433 of Yuno
With SFC pad
Recommended combination with • accessory trolley 1003.86A0
Combination with • adapter for CF plates 1001.64B0 • accessory adapter 1180.72A0
Suitable for use with • Yuno 1433.01/02
Available versions: • 1433.45AC insertable plate (long version), without head rest interface, length: 1,830 mm • 1433.45BC insertable plate (short version), with head rest interfaces for screw-fitted and push-fit head rests, length: 1,535 mm



1003.86A0
**Accessory trolley for
carbon-fiber insertable plates**

- Accessory trolley for storage and transport of carbon-fiber clip-on plates and related pads

Technical information

Width x length x height:
700 x 850 x 1,350 mm

4 double swivel castors,
2 castors with brakes

Suitable for

- carbon-fiber plates for insertion
- 1433.45AC/BC for Yuno II
- 1160.45AC/BC for 1160.10/30 table tops



1150.45BC
Carbon-fiber multi-purpose plate

- Radiolucent through a full 360°
- CT compatible

Technical information

Width x length: 600/395 x 900 mm

Material of pad: SFC

Dovetail interface for mounting specific head-positioning accessories, e.g. horseshoe head rest 1002.04B0

Mounting to a pin mounting point

Recommended combination with

- accessory adapter 1150.72B0

Suitable for use with

- 1150.19/25/30



1180.45A0
Carbon-fiber multi-purpose plate

- Radiolucent through a full 360°
- Pressure distribution thanks to SFC padding
- Easy mounting with Easy Click technology

Technical information

Interface: Magnus basic unit

Material of pad: SFC

Radiolucent area width x length:
460 x 650 mm

Screw-fitted head rest mounting point for additional accessories

Suitable for use with

- 1180.10A0



1433.33AC
Carbon-fiber back plate

- Radiolucent through a full 360°
- Attachment of accessories for positioning of the head possible
- Pressure distribution thanks to SFC padding
- Easy mounting with Easy Click technology

Technical information

Width x length: 530 x 495 mm

Material of pad: SFC

Screw-fitted and push-fit head rest mounting points for additional accessories, e.g. accessory adapter 6005.25A0 or skull clamp holder 1005.48B0

Mounting to a pin mounting point

Suitable for use with

- 1131, 1132, 1133, 1433, 7200, 7100, 1160.10/30, 7700



1180.37A0
Carbon-fiber back plate

- Radiolucent through a full 360°
- Short side rail for mounting side rail accessories
- Easy mounting with Easy Click technology

Technical information

Interface: Magnus joint module

Material of pad: SFC

EU side rail

Screw-fitted and push-fit head rest mounting points for additional accessories, e.g. accessory adapter 6005.25A0

Suitable for use with

- motorized joint module 1180.11A0



1132.45AC
Back plate

- Radiolucent through a full 360°
- Pressure distribution thanks to SFC padding
- Accessory adapter with side rails can be mounted

Technical information

Width x length: 520 x 1,155 mm

Material of pad: SFC

Mounting to a pin mounting point

Recommended combination with

- accessory adapter 1150.72A0
- accessory stand 1003.87A0

Suitable for use with

- 1118.06/10, 1131, 1132, 1133, 1433, 7000, 7100, 7200, 1160.10/30, 7700



1180.44AC
**Carbon-fiber multi-purpose plate
for lithotomy position**

- Radiolucent through a full 360°
- Extra wide sacral cut-out
- Attachment of leg plates or leg holders possible
- Easy mounting with Easy Click technology
- Pin mounting point for leg plates 1133.53, 1133.58, 1118.50B0

Technical information

Overall dimensions width x length:
576 x 852 mm

Radiolucent area width x length:
400 x 530 mm

Interface: Magnus basic unit

Material of pad: SFC

EU side rail (short)

Suitable for use with

- 1180.10A0



1003.87A0
Accessory stand for carbon-fiber back plate

- Space-saving storage of back plate 1132.45AC and related SFC pad

Technical information

Width x height x length:
540 x 900 x 700 mm

4 double swivel castors,
2 with brakes



1150.72A0
Accessory adapter

- For mounting side rail accessories to a carbon-fiber OR table top and carbon-fiber back plates
- Radiolucent in the middle segment

Technical information

EU side rail

Available versions:

- 1150.72A0 suitable for use with table top 1150.16B0 (length of side rail: 400 mm)
- CF plate 1132.45AC
- CF table top 1180.16A6
- 1150.72B0 suitable for use with CF plate 1150.45BC (length of side rail: 300 mm)



1150.70A0
Accessory adapter

- For mounting side rail accessories to OR table top 1150.16B0
- Radiolucent

Technical information

Length of side rail: 340 mm

EU side rail

Suitable for use with

- 1150.16B0
- 1180.16A6
- 1132.45AC



1180.72X0
Accessory adapter

- For mounting accessories to carbon-fiber plates
- Multi-purpose plate 1180.45A0
- Table tops 1180.16X0
- Radiolucent CF bar in the middle segment

Technical information

Length of side rail: 300 mm

EU side rail

Available versions:

- 1180.72A0 suitable for 1180.16A1/A2/A3/A5, CF plates 1180.45AC, insertable CF plates 1160.45XC, 1433.45XC and 1160.16A0/A1/A2 CF table tops
- 1180.72B0 suitable for transfer table top 1180.12A0



1180.31A0
Back plate

- For modular configuration of a Magnus Table Top
- Excellent pressure distribution thanks to hybrid pad
- Easy mounting with Easy Click technology

Technical information

Width x length:
580 x 365 mm (with side rail)

Interface: Magnus joint module

EU side rail

Suitable for use with

- motorized joint module 1180.11A0



1160.32AC
Extension plate

- For configuration of the operating table top
- Mounting to a pin mounting point
- Easy mounting with Easy Click technology

Technical information

Width x length x height:
590 x 250 x 158 mm (with side rail)

EU side rail

Suitable for use with

- 1160.10/30,
- mobile OR tables 1118.06/10, 1131, 1132,1133, 1433, 7200, 7100, 7000, 7700

Available versions

- 1160.32AC with IPC padding
- 1160.32BC with SFC padding
- 1160.32DC with IPC pad, washable



1160.33XC
Extension plate, short

- For configuration of the operating table top
- Especially for procedures on head and neck
- Mounting of retrofit kit for skull clamp holders possible
- Easy mounting with Easy Click technology

Technical information

Width x height x length:
590 x 158 x 160 mm

Weight: 5.6 kg

Maximum patient weight: 545 kg

EU side rail

Material of pad: SFC, IPC

Mounting to a pin mounting point

Suitable for use with

- 1160.10/30
- Mobile OR tables 1118.06/10, 1131, 1132,1133, 1433, 7200, 7100, 7000, 7700

Available versions

- 1160.33BC with SFC padding, EU side rail
- 1160.33DC with IPC padding, EU side rail , washable

Only in conjunction with

- 1160.32XC extension plate
- 1160.55XC seat plate extension
- 1131.55BC seat plate extension



1160.36A0
Retrofit kit for standard holder for DORO® skull clamps

- For combination with a standard holder for DORO® skull clamps and DORO® Adjustable Base Unit 1005.50A0

Technical information

Weight: 0.9 kg

Maximum patient weight: 250 kg

Comprises: 1 pair

Only in conjunction with

- 1160.33XC extension plate (may be ordered separately)



1160.36B0
Retrofit kit for standard holder for Sugita Skull Clamps

- Retrofit kit for combination with a standard holder for Sugita Skull Clamps

Technical information

Weight: 0.8kg

Maximum patient weight: 250 kg

Comprises: 1 pair

Only in conjunction with

- 1160.33XC extension plate (may be ordered separately)



1150.31B0
Back plate

- Back plate with claw and small pin interface

Technical information

Width x length: 540 x 300 mm

Adjustment angle + 90°/-60°

EU side rail

Mounting to claw mounting point

Suitable for use with

- 1150.20/30

Available versions

- 1150.31B0 with SFC padding
- 1150.31D0 with PUR padding



1160.35AC
Back plate

- Trapezoidal-shaped plate for close access to head and shoulder area
- Enabling procedures on the head, e.g. ENT procedures
- Additional square mounting point for head positioning accessories
- Easy mounting with Easy Click technology
- Square mount for head positioning accessories

Technical information

Width x length: 546 x 278 mm

Mounting to a pin mounting point

Suitable for use with

- 1160.10/30
- mobile OR tables 1118.06/10, 1133.22, 1433, 7000, 7100, 7200, 7700

Available versions

- 1160.35AC with IPC padding, washable
- 1160.35BC with SFC padding



1180.38A0
Universal extension plate

- For modular configuration of a Magnus OR Table Top
- Excellent pressure distribution thanks to hybrid pad
- Easy mounting with Easy Click technology

Technical information

Width x length: 580 x 550 mm (with side rail)

Interface: Magnus basic unit

Material of pad: SFC

EU side rail

Suitable for use with

- basic unit 1180.10A0



1180.32A0
Extension plate

- For modular configuration of a Magnus OR Table Top
- Excellent pressure distribution thanks to hybrid pad
- Easy mounting with Easy Click technology

Technical information

Width x length: 580 x 230 mm (with side rail)

Interface: Magnus basic unit

Material of pad: SFC

EU side rail

Suitable for use with

- basic unit 1180.10A0



1131.31BC
Extension plate

- For extending the table top
- Guide rails for X-ray cassettes
- Easy mounting with Easy Click technology

Technical information

Width x length: 590 x 250 mm

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1118, 1131, 1132.21/22, 1133
- 1433
- 7000, 7100, 7200, 7700
- 1160.10/30



1180.11A0
Motorized joint module

- High maneuverability provided by two separately controlled motors
- Excellent pressure distribution thanks to hybrid pad
- Easy mounting with Easy Click technology

Technical information

Width x length: 580 x 300 mm

Interface: Magnus basic unit

Material of pad: SFC

Including: 2 joint modules, 1 padded plate 1180.11A1

EU side rail

Suitable for use with

- 1180.10



1180.11B0
Motorized joint module

- High maneuverability provided by two separately controlled motors
- Improved pressure distribution thanks to hybrid pad
- Easy mounting with Easy Click technology
- Padded plate with sacral cut-out

Technical information

Width x length: 580 x 345 mm

Interface: Magnus basic unit

Material of pad: SFC

Including: 2 joint modules, 1 padded plate 1180.11B1

EU side rail

Suitable for use with

- 1180.10



**1180.11A1
Padded plate**

- Enables the mounting of back plates to the motorized joint module
- Improved pressure distribution thanks to hybrid pad
- Simple to mount

Technical information

Width x length: 520 x 320 mm

Material of pad: SFC

For combination with

- motorized joint module 1180.11X0

Suitable for use with

- 1180.10



**1180.11B1
Padded plate**

- Enables the mounting of leg plates to the motorized joint module
- Improved pressure distribution thanks to hybrid pad
- Simple to mount
- With sacral cut-out

Technical information

Width x length: 520 x 345 mm

Material of pad: SFC

For combination with

- motorized joint module 1180.11X0

Suitable for use with

- 1180.10



**1180.11C1
Urology seat plate**

- For use in conjunction with joint module 1180.11X0
- Padded plate with sacral cut-out enabling good access to the surgical field
- Plate is fitted with guide rails for insertion of a TUR rinsing set
- Good pressure distribution thanks to hybrid pad

Technical information

Width x length: 520 x 345 mm

Material of pad: SFC

Mounting to motorized joint module 1180.11X0

Recommended combination with

- TUR rinsing set 1003.45D0



**1180.55A0
Universal module**

- Extension of the table top
- With sacral cut-out
- Pressure distribution thanks to SFC padding
- Easy mounting with Easy Click technology

Technical information

Width x length: 582 x 160 mm

Interface: Magnus basic unit

Material of pad: SFC

EU side rail

Suitable for use with

- 1180.10



**1160.55AC
Seat plate extension**

- For extending the table top e.g. for X-ray and TUR applications
- Seat plate with gyn. cut-out
- Fitted with guide rails for insertion of X-ray cassettes
- Easy mounting with Easy Click technology
- For use with TUR procedures, please order rail set 1160.56A0

Technical information

Width x height x length: 590 x 158 x 250 mm

Maximum patient weight: 454 kg

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1160.10/30
- 1118.06, 1118.10, 1133.22, 1433.02, 7000, 7100, 7200, 7700

Available versions:

- 1160.55AC with IPC padding
- 1160.55BC with SFC padding
- 1160.55DC with IPC padding, washable

May be ordered separately:

- 1160.56A0 guide rails for TUR rinsing set, pair



**1150.55HC
Seat plate extension**

- For extending the table top
- Especially for X-ray and TUR applications
- With sacral cut-out
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 573 x 260 mm

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1150.19/30

Variant:

- 1150.55DC washable



1131.55BC
Seat plate extension

- For extending the table top for radiolucency of the pelvic area
- Guide rails for X-ray cassettes or TUR rinsing set
- With gyn. cut-out
- Easy mounting with Easy Click technology

Technical information

Width x length: 590 x 250 mm

Including: 2 guide rails for combination with TUR rinsing set 1003.45D0

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1118, 1131, 1132.21/22, 1133, 1433, 7000, 7100, 7200, 7700, 1160.10/30



1130.55B0
Seat plate extension

- For extending the table top for radiolucency of the pelvic area , e.g. for urological interventions
- Guide rail for X-ray cassettes or TUR rinsing set

Technical information

Width x length: 573 x 260 mm

EU side rail

Mounting to a pin mounting point

Recommended combination with


- TUR rinsing set 1003.45D0

Suitable for use with

- 1118.01 and 1118.03

Available versions:

- 1130.55A0 with PUR pad
- 1130.55B0 with SFC pad

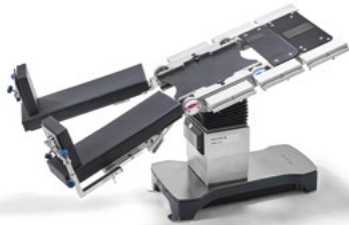

1001.74B0
Bariatric OR table set

- For extension of the table top area
- For placement and positioning of adipose patients
- Set of 3 pair table width extensions and 2 extra wide and extra thick SFC pads

Technical information

Maximum patient weight: 545 kg

Width x length of table width extenders: 133 x 300 mm

Width x height x length of SFC pads:
a) 760 x 100 x 870 mm
b) 760 x 100 x 250 mm

EU side rail

Material of pad: SFC

Comprises 3 pair of width extenders and 2 pads

Mounting to side rails or table top

Integrated fixation clamps

Only in conjunction with

- 7700, 1160.10, 1160.32BC, 1160.55BC

Available versions:

- 1001.74B0 SFC pad



1005.78C0
Pediatric operating table top

- For placement and positioning of infants and small children
- Good access thanks to slim construction
- Side rail accessories can be mounted as usual

Technical information

Width x height x length: 410 x 120 x 910 mm

Including: 1 table top and 1 support bar

Material of pad: SFC

2 EU side rails

Mounting to a pin mounting point

Suitable for use with

- OR tables 1118, 1131, 1132, 1133, 7000, 7100, 7200, 7700, 1160.10/30

05 Cross-discipline accessories

Leg positioning



Safe and comfortable leg positioning supports optimum treatment outcomes for the patient. At the same time, lightweight leg plates or transfer boards considerably ease the frequent changeover of OR tables, as is usual in, for example, gynaecology.

For every type of leg position, Maquet offers suitable accessories. From simple leg plates and flexible leg plates in four segments, through to leg holders for gynecological and urological interventions. Special carbon-fiber materials meet the high demands of imaging procedures.



- 1003.80A0
Leg holder
- Adjustable with pneumatic support
 - Foot and calf support is automatically adapted when the position of the leg is adjusted
 - One-hand operation

Technical information
Width x length x height: 1,027 x 240 x 185 mm
Net weight: 6 kg
Maximum patient weight: 250kg
1 pair supplied
1 pair of blade clamps (1003.24A0) is included in the delivery



- 1003.81A0
Knee crutches, standard
- To position and support the leg in lithotomy position
 - Side wings for stable support of the calf and knee area

Technical information
Width x length x height: 220 x 385 x 566 mm
Net weight: 5.1 kg
1 pair supplied
Suitable for • patient weight up to 180 kg
Only in combination with • leg holder (1003.80A0)



- 1003.82A0
Knee crutches, adipositas
- To position and support the leg in lithotomy position
 - Side wings for stable support of the calf and knee area
 - Particularly suitable for obese patients

Technical information
Width x length x height: 260 x 385 x 566 mm
Net weight: 6.9 kg
1 pair supplied
Suitable for • patients up to 250 kg
Only in combination with • leg holder (1003.80A0)



- 1003.24A0
Clamp with blade interface
- For mounting accessories to the side rail
 - Especially for mounting accessories with blade interface
 - Enabling close rotation point of the leg holder to the hip joint

Technical information
Width x length x height: 60 x 85 x 130 mm
Maximum load: 49.2 kg



1008.08A0
KTEK 250 Stirrup – Junior

- To position and support the leg in lithotomy position
- For small patients
- One-hand operation
- Boot without side wing

Technical information

Width x height x length:
28 x 40 x 100 mm

Weight: 4.6 kg per leg holder

Maximum patient weight: 113 kg

Lithotomy range +85 / - 55°,
abduction range +25 / -9°

Material of pad: soft foam

Material of strap: silicone strap

Packing unit: 1 pair (2 pieces)

No mechanical cleaning

No clamps integrated

Recommended combination with

- KTEK Stirrup Cart 1008.02B0

Only in combination with

- secure release clamps 1008.05A0
for EU side rail
(may be ordered seperately)



1008.09A0
KTEK Comfort 800 Stirrup

- To position and support the leg in lithotomy position
- Adjustable with pneumatic support
- Foot and calf support is automatically adapted when the position of the leg is adjusted
- One-hand operation
- Boot with side wing for supporting calf and knee area

Technical information

Width x height x length:
28 x 40 x 100 mm

Weight: 6.6 per leg holder

Maximum patient weight: 363 kg

Lithotomy range +85 / - 55°,
abduction range +25 / -9°

Material of pad: soft foam

Material of strap: silicone strap

Packing unit: 1 pair (2 pieces)

No mechanical cleaning

No clamps integrated

Can be mounted to a side rail with a
blade clamp

Recommended combination with

- KTEK Stirrup Cart 1008.02B0

Only in combination with

- secure release clamps 1008.05A0
for EU side rail
(may be ordered seperately)



1008.05A0
Secure release clamp

- Clamp with blade interface
- For mounting the KTEK Stirrups
to the side rail

Technical information

Width x height x length:
54 x 54 x 76 mm

Packing unit: comprises 1 clamp

Suitable for

- combination with KTEK Stirrups



1008.02B0
KTEK Stirrup Cart

- For storage and transport of 2 KTEK Stirrups
- With small storage bin for accessories

Technical information

Width x height x length:
549 x 549 x 1035 mm

Weight: 5 kg

1 piece supplied

No mechanical cleaning

Combination with

- KTEK Stirrups and
secure release clamps

Only in conjunction with

- KTEK Stirrups



1002.2300
Leg holder

- For holding the leg with heel loops
- Can be adapted for different leg lengths
- Can be rotated

Technical information

Height: 865 mm

Comprises: 1 leg holder

Loops are washable

Spare part: loop 3101.0334

Can be mounted to a side rail

Standard clamp may be
ordered separately



1001.65A0
Leg holder

- To position and support the leg
in lithotomy position
- Goepel version
- Good positioning comfort ensured
by sufficient positioning space and
appropriate geometry
- Adjustment of calf plate via ball joint

Technical information

Comprises: 1 leg holder,
1 support strap

Inlay: SFC padding

Spare part: strap 3101.2754,
pad 9093.4593

Can be mounted to a side rail
with a clamp

Radial setting clamp (1003.23C0)
may be ordered separately



**1001.73A0
Leg holder**

- Extra wide leg holder
- To position and support the leg in lithotomy position
- Particularly suitable for obese patients

Technical information

Comprises: 1 leg holder, 2 fastening straps

Inlay: SFC padding

Can be mounted to a side rail with a setting clamp

Radial setting clamp (1003.23C0) may be ordered separately



**1008.10A0
Knee crutch lift assist**

- Gas spring supported adjustment
- Inlay: soft foam pad

Technical information

Width x height x length: 356 x 152 x 720 mm

Weight: 4.6 kg

Maximum patient weight: 227 kg

Material of device: stainless steel, aluminium

Material of pad: soft foam

Packing unit: 1 pair (2 pieces)

No mechanical cleaning

Can be mounted to a side rail with a blade clamp

No clamps integrated

Only in combination with

- secure release clamps 1008.05A0 for EU side rail (may be ordered seperately)



**1130.65A0
Transfer board**

- For placement and positioning of legs
- Radiolucent
- Lightweight accessory
- Quick and easy mounting and dismounting

Technical information

Width x length (with pins): 505 x 700 mm

Material of pad: SFC

Mounting to a pin mounting point

Suitable for use with

- 1118.01, 1118.06, 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



**1132.65A0
Transfer board**

- For pre- and post-operative positioning of the legs in lithotomy position
- Lightweight accessory for back-friendly working
- Quick and easy to mount and remove

Technical information

Width x length (without pins): 493 x 930 mm

Mounting to a pin mounting point

Suitable for use with

- 1118.01, 1118.06, 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



**1180.57A0
Transfer board**

- For pre- and post-operative positioning of the legs in lithotomy position
- Lightweight accessory for back-friendly working
- Quick and easy to mount and remove thanks to Easy Click technology

Technical information

Width (with pad) x length: 577 x 925 mm

Interface: Magnus joint module

Recommended combination with

- transfer board pad 1180.56A0

Suitable for use with

- motorized joint module 1180.11B0
- universal module 1180.55A0/A1



**1180.56A0
Transfer board pad**

- SFC pad for improved positioning comfort of the legs

Technical information

Recommended combination with

- transfer board 1180.57A0



**1180.51A0
Adapter clamp with transfer board**

- Lightweight transfer board enables quick and easy attachment and detachment
- For pre- and post-operative positioning of legs in lithotomy position

Technical information

Transfer board width x length: 510 x 900 mm

Interface: Magnus joint module 1180.11B0

Comprises: 2 adapter clamps, 1 transfer board

Material of pad: SFC

EU side rail

Suitable for use with

- motorized joint module 1180.11B0



1130.70AC
Maquetmatic leg plate

- Anatomical positioning of the legs using a four-segment leg plate
- Enables good positioning comfort in lithotomy position
- Leg plates can be folded independently and abducted outwards
- Thigh and calf plates can be removed
- Knee and calf crutches can be mounted in place of thigh and calf plates

Technical information

Width x length: 785 x 940 mm

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Recommended combination with

- knee crutches 1005.01B0
- filler piece 1130.71B0

Suitable for use with

- 1118.01, 1118.06, 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 7700
- 116010/30



1130.71B0
Filler piece

- Pressure distribution thanks to SFC padding

Technical information

Only in conjunction with

- Maquetmatic 1130.70AC



1005.01B0
Knee crutches

- To support and position the legs of a patient

Technical information

Comprises: 1 pair of knee crutches and 2 support straps

Spare part: strap 3101.2754, pad 9093.4593

Combination with

- Maquetmatic leg plates



1180.48B0
Lightweight leg plate

- Easy access thanks to abductable leg plates
- Lightweight leg plate supports back-friendly working
- Pressure distribution thanks to SFC padding

Technical information

Width
- with pad: 540 mm
- with side rail: 585 mm
Length: 922 mm

Interface: Magnus joint module

Material of pad: SFC

EU side rail

Suitable for use with

- motorized joint module 1180.11B0



1118.60B0
Pair of leg plates

- Gentle adjustment with pneumatic support
- With joint for abduction
- Recommended for use with operating tables without motorized leg adjustment

Technical information

Width x length: 570 x 750 mm

Adjustment angle:
abduction 90° / downwards 90°

Maximum patient weight: 250 kg

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1118.06, 1118.10
- 7000.01, 7100, 7200
- 1133.22
- 1433.02
- 1160.10/30



1150.73BC
Pair of leg plates

- Multi-functional four-piece leg plate
- Can be used for knee-elbow position
- Thigh and calf plate can be adjusted separately
- Calf plate can be removed
- Particularly stable

Technical information

Width x length: 548 x 923 mm

Maximum patient weight: 380 kg

Weight: 11 kg / plate

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1150.19/30
- Available versions:**
- 1150.73BC with SFC padding
- 1150.73DC with PUR padding



1118.50B0
Pair of leg plates

- Gentle adjustment with pneumatic support
- Abduction of leg plates via the leg plate mounting point
- Recommended for use with operating tables without motorized leg adjustment

Technical information

Width (with side rail) x length:
570 x 750 mm

Weight: 11 kg per plate

Maximum patient weight: 135 kg

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Recommended combination with

- 1118.01B0/K0, 1118.06B0/K0

Suitable for use with

- 1118.01, 1118.06, 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 1160.10/30



1007.42AC
Pair of leg plates

- Good pre- and post-operative positioning comfort for the patient
- Increased working comfort for medical staff during pre- and post-operative period
- Removable filler pieces enable easy adaptation to 1-hole and 3-hole CF plate
- Dedicated leg plate for Yuno II

Technical information

Width x height x length:
253 x 166 x 1,320 mm

Weight: 14 kg (pair)

Maximum patient weight: 250 kg

EU and US side rail

Material of pad: SFC

Mounting to a pin mounting point

Recommended combination with

- Yuno II
- Carbon-fiber extension plate 1433.66AC/BC

Suitable for use with

- 1433.02



1133.67BC
Carbon-fiber pair of leg plates

- Lightweight 2-piece leg plate
- Radiolucent through a full 360°
- Can be abducted
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 235 x 910 mm

Maximum patient weight: 225 kg

Material of pad: SFC

Material of plate: carbon fiber

Mounting to a pin mounting point

Suitable for use with

- 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



1160.50AC
Pair of leg plates

- Double joint abduction enables good access between the legs of the patient
- Weight-optimized leg plate allows back- friendly mounting to OR table

Technical information

Width x length: 580 x 850 mm

Weight: 19.5 kg (pair)

Maximum patient weight: 545 kg

Material of pad: IPC, SFC

EU side rail

Cleaning options:
1160.50DC mechanical cleaning

Mounting to a pin mounting point

Suitable for use with

- 7000, 7100, 7200
- 7700
- 1433.02
- 1160.10/30

Available versions:

- 1160.50AC with IPC pad
- 1160.50BC with SFC pad
- 1160.50DC with IPC pad, washable



1133.58BC
Leg plate

- One-piece leg plate
- Pressure distribution thanks to SFC padding
- Quick and easy to mount

Technical information

Width x height x length:
585 x 175 x 651 mm

Material of pad: SFC

Maximum patient weight: 454 kg

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1118.06, 1118.10
- 1131, 1132, 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



1133.53BC
Pair of leg plates

- Double joint abduction enables good access between the legs of the patient
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 250 x 680 mm

Weight: 9.5 kg / plate

Maximum patient weight: 454 kg

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1118.06, 1118.10
- 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



1150.64B0
Pair of leg plates

- For positioning of legs in conjunction with procedures with traction
- Lightweight leg plates
- Quick unlocking for detachment

Technical information

Width x length: 250 x 680 mm

Weight: 9 kg (pair)

Maximum patient weight:
225kg, 1150.64E0 = 250kg

Material of pad: SFC

Mounting to square mount
of OR table or traction device

Packing unit: 1 pair

Combination with

- 1118, 1132, 1133, 7000, 7100, 7200 with extension device 1419.X0
- 1150.19/30 with extension device 1150.59AC
- 1180.10 with extension device 1180.19A0
- 1150.64E0 only in combination with adapter plate for traction 1433.41AC or support plate for extension 1160.59AC

Suitable for use with

- 1150.20

Available versions:

- 1150.64B0 with SFC padding 80 mm
- 1150.64E0 with SFC padding 40 mm



1180.54A0
Pair of leg plates

- Multi-functional four-piece leg plate
- Can be used for knee-elbow position
- Thigh and calf plate can be adjusted separately
- Calf plate can be removed
- Particularly stable
- Simple mounting using Easy Click technology

Technical information

Width with pad: 520 mm,
width with side rail: 580 mm,
length: 922 mm

Interface: Magnus joint module

Weight: 12.5 kg / plate

Maximum patient weight 380 kg

EU side rail

Suitable for use with

- motorized joint module 1180.11B0

Available versions:

- 1180.54A0 with PUR pad
- 1180.54B0 with SFC pad



1180.41A0
Pad for rectal positioning

- Body support for knee-elbow position
- For anatomical positioning during rectal interventions
- Supports the pelvis and reduces abdominal pressure

Technical information

Width x length x height:
627 x 380 x 395 mm

Material of pad: SFC

Combination with

- joint module 1180.11X0
- leg plates 1180.54X0



1133.73BC
Pair of leg plates

- Multi-functional four-piece leg plate
- Can be used for knee-elbow position
- Thigh and calf plate can be adjusted separately
- Calf plate can be removed

Technical information

Width x length: 583 x 850 mm

Weight: 12 kg per plate

Maximum patient weight: 454 kg

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Suitable for use with

- 1133
- 1433
- 7000, 7100, 7200
- 7700
- 1160.10/30



1133.73AC
Pair of leg plates

- Multi-functional four-piece leg plate
- Can be used for knee-elbow position
- Thigh and calf plate can be adjusted separately
- Calf plate can be removed

Technical information

Width x length: 853 x 882 mm

Weight: 12 kg per plate

Maximum patient weight: 454 kg

Material of pad: IPC

Suitable for use with

- 1133.22
- 1433.02
- 7000, 7100, 7200
- 7700
- 1160.10/30



1001.86C0
Foot plate

- Foot support in a standing position
- Axially correct setting using rotating joint
- In standing position additional support of the leg with a leg strap is needed, e.g. leg strap 1001.53A0
- Set of 2

Technical information

Width x height x length:
335 x 380 x 315 mm

Foot plate width x length:
300 x 300 mm

Weight: 7.3 kg

Maximum patient weight: 250 kg

Material of pad: SFC

Can be mounted to a side rail

Integrated standard clamp



1001.66A0
Leg support

- For placement and positioning of the patient's leg during skin disinfection
- Foot holder can be moved
- Height adjustable using the clamp

Technical information

Bar width x length: 410 x 710 mm

Height of pad: 80 mm

Material of pad: PUR

Can be mounted to a side rail with a clamp

Standard clamp may be ordered separately



1008.06B0
KTEK Nissen Straps, standard

- Securing of the patient's leg in the thigh area
- Pressure relief thanks to soft padding

Technical information

Width x height x length:
125 x 12 x 2430 mm


Product weight: 2 pieces = 9kg

Weight: 4.5 kg/ strap

Maximum patient weight: 227 kg

Packing unit: 1 pair including 2 side rail clamps

No mechanical cleaning

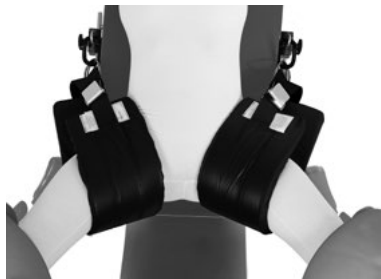

1008.06D0
KTEK Nissen Straps, bariatric

- Securing of the patient's leg in the thigh area
- Pressure relief thanks to soft padding

Technical information
Width x height x length: 126 x 13 x 2430 mm
Product weight: 2 pieces = 9kg
Weight: 12 kg per plate
Maximum patient weight: 363 kg
Packing unit: 1 pairs and 4 clamps
No mechanical cleaning



1001.57A0
Leg strap

- For fixing the patient's legs on a multiple-section leg plate


Technical information
Width x length: 120 x 1,300 mm
Comprises: 1 leg strap
Material of strap: polyester fabric
Strap is washable (60 °C)
Compatible with
• 2-piece and 4-piece leg plates



1001.4700
Leg restraint cuff

- For securing the lower leg on a split leg plate

Technical information
Comprises: 1 leg restraint cuff
Strap is washable (60 °C)
Compatible with
• 2-piece and 4-piece leg plates

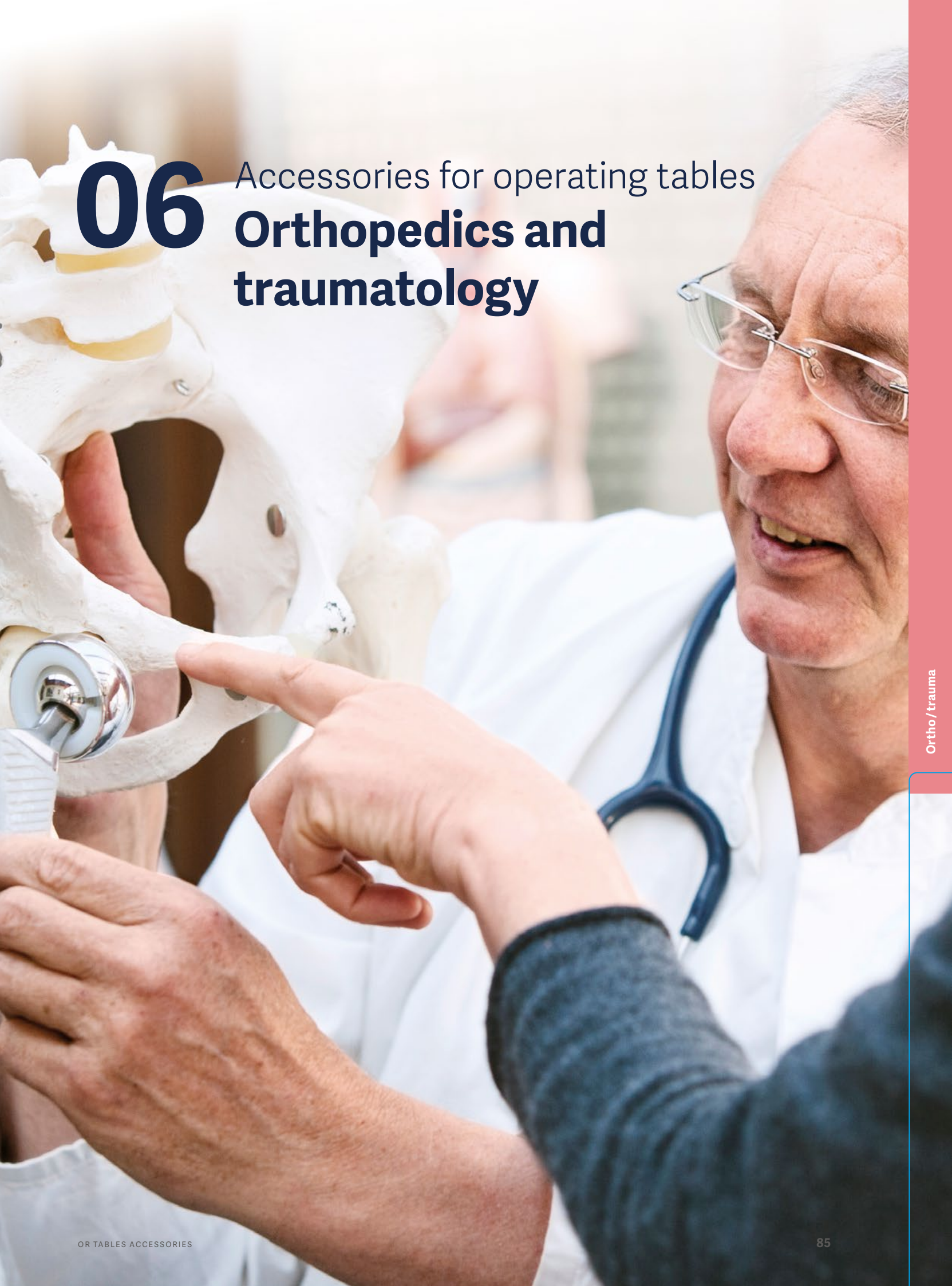

06 Accessories for operating tables

Orthopedics and traumatology

Modern orthopedics and traumatology require particularly demanding positions and accessories that are coordinated for interventions.

Maquet’s accessories for orthopedics and traumatology cover the full range of standard applications. In conjunction with experts and users, innovative solutions are developed which support new surgical procedures. More advanced imaging procedures, such as 3D X-ray technology, again increase the demands placed on OR tables and accessories in orthopedics. Maquet accessories for orthopedics and traumatology meet the requirements for optimum X-ray results: from high-quality carbon-fiber materials through to the shape of the modules and to the radiolucency of the pad material.



1003.60D0
Clip-on device for minimally invasive surgery

- For hip arthroscopy and minimally invasive hip endoprosthesis
- Radiolucent CF bar
- High stability thanks to wide tracking castors
- Back-friendly handling thanks to pneumatic support
- Adjustment of CF bar horizontally and vertically possible
- Can be adapted to the length of the patient’s leg and foot position by infinite sliding of the skid
- Mounting to a table top or extension device

Technical information

Length of carbon-fiber bar: 1,465 mm

Comprises: 1 carbon-fiber bar, 1 sled with screw tension device and rotating and tilting clamp, 1 hypomochlion

Recommended combination with

- traction boot 1003.75A0 and/or foot plate 1001.97A0

Suitable for use with

- 1150.20, table tops for extension 1160.59, 1433.41
- traction devices 1419.01X0 1150.59, 1180.19



1180.19A0
Extension device

- Quick conversion of a universal Magnus Table Top into a table top for traction
- Simple mounting of bars
- Good C-arm access
- Simple transportation and mounting using the trolley
- Space-saving storage of all accessories in the trolley

Technical information

Interface: Magnus basic unit

Including: 1 trolley, 2 telescopic bars (1 short, 1 long), 2 bar extensions, 1 screw tension device, 1 foot plate mount, 1 countertraction post, 1 side rail extension, 1 rotating and tilting clamp, 1 traction stirrup clamp with rotation, 1 support bar, 1 radial setting clamp

Material of pad: SFC

EU side rail

Recommended combination with

- traction boot 1003.75A0
- foot plate 1001.97A0
- pair of leg plates 1150.64B0

Suitable for use with

- 1180.10



1150.59AC
Extension device

- Quick conversion of a universal OR table to a traction table
- Easy mounting to OR table
- Good C-arm access
- Combination with leg plates (1150.64B0) possible
- Simple transportation and mounting using the trolley
- Space-saving storage of all accessories in the trolley

Technical information

Comprises: 1 pair of adapters, 1 support frame, telescopic bars (1 short, 1 long), 1 screw tension device, 2 bar extensions, 1 foot plate mount, 1 countertraction post, 1 side rail extension, 1 rotating and tilting clamp, 1 traction stirrup clamp with rotation, 1 radial setting clamp, 1 support bar, 1 trolley

Material of Pad: SFC

EU side rail

Recommended combination with

- traction boot 1003.75A0

Suitable for use with

- 1150.19/30



1160.59BC
Support plate for extension

- Table top for orthopedic procedures
- Quick and easy conversion of a universal OR table top into a table top for procedures with traction

Technical information

Mounting to interface for insertion

Material of pad: SFC

Recommended combination with

- pair of leg plates 1150.64E0
- traction boot 1003.75A0
- foot plate 1001.97A0

Suitable for use with

- 1160.10/30
- 7700

Available versions:

- 1160.59AC without accessories
- 1160.59BC including accessories: 1 trolley , 1 table top for traction, 1 set of traction bars, telescopic bars (1 short, 1 long), 1 counter-traction post, 1 padded roll, 1 left adapter bar, 1 right adapter bar, 1 bar elongation, 1 screw tension device, 1 foot plate mount, 1 side rail extension, 1 rotating and tilting clamp, 1 traction stirrup clamp



1419.01X0
Traction device

- Quick conversion of a normal OR table to an extension table
- Attachment of leg plates during preparation phase possible
- Simple mounting of traction bars
- Good C-arm access
- Space-saving storage of all accessories in the trolley
- Back-friendly mounting using the trolley

Technical information

Comprises: 1 trolley with basket, 1 pair of adapters, telescopic bars (1 short, 1 long), 1 screw tension device, 1 foot plate mount, 2 support bars, 2 radial setting clamps, 1 countertraction post, 1 padded roll for countertraction post, 1 side rail extension

Material of pad: SFC

Recommended combination with

- traction boot 1003.75A0
- foot plate 1001.97A0
- leg plates 1150.64X0

Available versions:

- 1419.01B0 suitable for 7000, 7100, 7200 on reverse mounting point only
- 1419.01JC suitable for 7000, 7100, 7200 on normal and reverse mounting point
- 1419.01C0 suitable for 1118.01, 1118.06, 1118.10
- 1419.01HC suitable for 1133.22 on reverse mounting point



1433.41AC
Adapter plate for extension

- Simple mounting of traction bars
- Good C-arm access
- Combination with existing extension accessories possible
- Easy to use

Technical information

Material of pad: SFC

2 removable seat plate segments

Mounting to interface for insertion

Suitable for use with


- 1433.01 and 1433.02

Available versions:

- 1433.41AC without accessories for traction
- 1433.41BC with accessories for traction: telescopic bars (1 short, 1 long), 1 screw tension device, 1 foot plate mount, 2 support bars, 2 radial setting clamps, 1 countertraction post, 1 side rail extension, 2 bar extensions, 1 rotating and tilting clamp, 1 traction stirrup clamp with rotation

Recommended combination with

- 1160.57A0 accessory trolley (may be ordered separately)

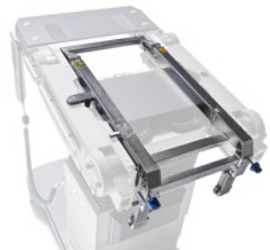**1003.72AC**
Traction adapter

- Easy conversion from Maquet Corin or Maquet Otesus to a specialized orthopedic setup
- Stainless steel frame with transponder technology
- Insertable interface for extension plate 1433.66 without traction accessories

Technical information

L x W x H: 680 x 535 x 250 mm

Suitable for

- 1160.10, 1160.30, 7700 in combination with 1433.66AC/BC

Material of device

- Stainless steel

**1160.57A0**
Accessory trolley

- Space-saving storage of accessories for traction

Technical information

Including: 1 wire basket

Width x length: 850 x 650 mm

4 double swivel castors,
2 with brakes

Recommended combination with

- 1160.59AC and 1433.41AC

**1007.40B0**
Traction bar with ball joint

- For hip arthroscopy, minimally invasive hip procedures and Direct Anterior Approach (DAA)
- For femur fractures
- Adjustment of bar: horizontally and vertically
- Can be adapted to the length of the patient's leg by sliding the slider
- Carbon-fiber bar is radiolucent

Technical information

L x W x H : 2284 x 331 x 573 mm

Suitable for

- 1160.10, 1160.30, 7700 in combination with 1003.72AC

Comprises

- Slider

**1433.66AC**
Carbon-fiber extension plate with 3 holes

- Radiolucent through 360°
- Good access to pelvis thanks to variable positioning of the countertraction post
- 2 removable seat plate segment items possible to be mounted: countertraction post
- Countertraction post included

Technical information

Width x height x length:
515 x 80 x 1,150 mm

Weight: 6.8 kg

Maximum patient weight: 250 kg

Material: carbon fiber, with SFC pad

No mechanical cleaning

Mounting to interface for insertion

Suitable for use with

- 1433.01/02

**1007.41B0**
Traction bar

- For hip arthroscopy and femur fractures
- Adjustment of bar: horizontally
- Can be adapted to the length of the patient's leg by sliding the slider
- Carbon-fiber bar is radiolucent

Technical information

L x W x H : 2258 x 287 x 351 mm

Suitable for

- 1160.10, 1160.30, 7700 in combination with 1003.72AC

Comprises

- Slider

**1433.66BC**
Carbon-fiber extension plate with 1 hole

- Radiolucent through 360°
- Pressure distribution thanks to SFC padding
- 2 removable seat plate segment items possible to be mounted: countertraction post
- Countertraction post included

Technical information

Width x length x height:
515 x 1,150 x 80 mm

Weight: 6.9 kg

Maximum patient weight: 250 kg

Material: carbon fiber, with SFC pad

No mechanical cleaning

Mounting to interface for insertion

Suitable for use with

- 1433.01/02

**1007.46B0**
Trolley for traction bars

- For storage and transportation
- Enables easy mounting of traction bars to Maquet Otesus & Maquet Corin

Technical information

Comprises: 1 trolley

Only in combination with

- traction bars 1007.41B0 and 1007.40B0

**1007.40A0**
Traction bar with ball joint

- For hip arthroscopy and minimally invasive hip procedures
- For femur fractures
- Adjustment of bar: horizontally and vertically
- Can be adapted to the length of the patient's leg by sliding the slider
- Carbon-fiber bar is radiolucent

Technical information

Length: 2,471 mm

Material of bar: carbon-fiber

Recommended combination with

- screw tension device 1007.43A0
- accessory trolley for traction bars 1007.47A0

Suitable for

- Yuno II 1433.02

**1007.41A0**
Traction bar

- For hip arthroscopy and femur fractures
- Adjustment of bar: horizontally
- Can be adapted to the length of the patient's leg by sliding the slider
- Carbon-fiber bar is radiolucent

Technical information

Length: 2,454 mm

Material: carbon fiber

Recommended combination with

- screw tension device 1007.43A0
- vertical adjustment 1007.51X0
- screw tension device 1003.7300
- accessory trolley for traction bars 1007.47A0

Suitable for

- Yuno II 1433.02

**1007.46A0**
Trolley for traction bar

- For storage and transportation
- Enables easy mounting of traction bars to Yuno II

Technical information

Comprises: 1 trolley

Only in combination with

- traction bars 1007.40A0 and 1007.41A0

**1007.44A0**
Sliding clamp with side rail

- The side rail of the clamp can be used to mount accessories to a traction bar

Technical information

EU side rail

Suitable for

- traction bars 1007.40A0 and 1007.41A0

**1007.47A0**
Accessory trolley

- Space-saving storage and transportation of accessories for Yuno II

Technical information

Width x height x depth:
750 x 1,450 x 820 mm

4 double swivel castors,
2 with brakes

Suitable for

- Yuno II accessories

**1433.42A0**
Femur hook

- To hold the thigh bone for minimally invasive hip procedures
- Anatomical design

Technical information

Comprises: femur hook (left and right), 1 support arm, 1 clamp with square mount

Material: stainless steel

All components can be sterilized

Mounting to a side rail

Suitable for use with

- 1433.01 and 1433.02

**1003.55A0**
Traction bar, carbon-fiber version

- 360° radiolucent with exception of the mounting point area

Technical information

Length: 1,374 mm

Including: 1 pair

Material: carbon fiber

Recommended combination with

- sliding clamp 1003.58B0 for attachment of screw tension device

Suitable for use with

- 1150.20, extension devices 1419.01X0, 1150.59AC, 1180.19A0
- adapter plate for traction 1433.41AC
- support plate for extension 1160.59AC

**1003.58B0**
Sliding clamp

- To mount a screw tension device to a carbon-fiber bar
- With side rail mounting point

Technical information

To be used in combination with CF traction bar 1003.55A0 only

EU side rail

Combination with

- screw tension device 1003.3700
- screw tension device 1007.43A0 plus vertical adjustment

**1004.85B0**
Countertraction post for femur

- For positioning the patient's thigh during traction in lateral position
- Radiolucent curved perinial support
- Perinial arc is height-adjustable using the hand crank

Technical information

Comprises: 1 countertraction post for femur and 1 extension plate

Material of pads: PUR

Mounting to lateral square mount

Suitable for use with

- table top 1150.20B0/D0
- extension devices 1419.01X0, 1150.59AC, 1160.59AC, 1180.19A0, 1433.41AC
- adapter clamp 1007.50AC



1007.50AC
Adapter clamp

- For combination with additional accessories
- Lateral positioning of the patient during femur treatment is possible

Technical information

- Suitable for Yuno II
- Mounting to a pin mounting point

Only in conjunction with

- extension plates 1433.66BC (1 hole) and 1433.66AC (3 holes)
- femur countertraction post 1004.85X0



1003.57B0
Seat plate extension, carbon fiber

- Easy access to pelvis thanks to variable positioning of the countertraction post
- Radiolucent through 360°
- Lateral carbon-fiber supporting roll for supporting the femur
- Pressure distribution thanks to SFC padding

Technical information

- Material: carbon fiber, PUR, SFC

Recommended combination with

- clip-on device for MIS interventions 1003.60D0
- extension device 1180.19X0

Suitable for use with

- 1150.20
- extension devices 1150.59AC, 1180.19A0
- adapter plate for traction 1433.41AC
- support plate for extension 1160.59AC



1004.8700
Leg plate clamp

- To accommodate leg plates on a traction bar
- Support of anatomical leg positioning

Technical information

- Set of 2
- Mounting to standard traction bars of traction table tops and traction devices

Recommended combination with

- leg plate 1150.64X0



1433.67A0
Tibia adapter

- For mounting the countertraction post for tibia and fibula 1003.50X0

Technical information

- Including: SFC padding suitable for extension plate with 1 hole

Only in conjunction with

- CF extension plates 1433.66AC/BC



1003.50C0
Countertraction post for tibia and fibula

- For positioning and securing the leg in the event of calf fractures
- Easy repositioning of the bone using the height and inclination adjustment of the padded roll
- Radiolucent in the padded roll region

Technical information

- Interface: extension table top
- Including: 1 padded roll
- Material of pad: PUR

Suitable for use with

- 1150.20
- extension devices 1419.01X0, 1150.59AC, 1180.19A0
- adapter plate for traction 1433.41AC
- support plate for extension 1160.59AC

Available versions

- 1003.50A0 with horizontal guide (50 x 50 mm) for telescopic bar with 40 mm clear inner dimension
- 1003.50C0 with horizontal guide (40 x 40 mm) for telescopic bar with 35 mm clear inner dimension



1004.93B0
Condyle fixation

- For exact fixation of the knee
- Individually adjustable thanks to asymmetrical lateral discs

Technical information

- Comprises: 2 pieces
- Material of Pad: PUR

Only in conjunction with

- countertraction post for tibia and fibula 1003.50X0



1003.5100
Guiding bar

- For correction of the tensile direction on the calf

Technical information

Suitable for use with

- 1150.20
- extension devices 1419.01X0, 1150.59AC, 1180.19A0
- adapter plate for traction 1433.41AC
- support plate for extension 1160.59AC



1004.42B0
Joint-supporting arm

- For mounting a telescopic bar and a screw tension device, supporting the execution of traction
- Precise setting thanks to horizontal bar that can be angled

Technical information

Only in conjunction with

- countertraction post for tibia and fibula 1003.50C0

Available versions

- 1004.42B0 for combination with telescopic bar 35 mm
- 1004.42A0 for combination with telescopic bar 40 mm



1003.34A0
Rotating and tilting clamp

- For mounting a traction boot or a foot plate
- Precise setting of the boot using rotation and tilt function

Technical information

Mounting to the side rail mounting point

Combination with

- screw tension device 1003.3700
- foot plate mount 1003.49B0
- screw tension device 1007.43A0



1003.54A0
Bar extension

- To shorten or lengthen the telescopic bar during traction positions
- Precise adaptation of the traction bar to the size of the patient, e.g. for children

Technical information

Recommended combination with

- screw tension device 1003.3700

Suitable for use with

- 1150.20
- extension devices 1419.01X0, 1150.59AC, 1180.19A0
- adapter plate for traction 1433.41AC
- support plate for extension 1160.59AC



1003.35A0
Traction stirrup clamp with rotation

- For mounting and positioning of tension bars for traction
- Precise setting of the bar using the rotation function

Technical information

Mounting to the side rail mounting point

Combination with

- screw tension device 1003.3700
- foot plate mount 1003.49B0
- screw tension device 1007.43A0



1003.3700
Screw tension device

- For the execution of traction
- Manual fine adjustment of traction using a hand crank
- Adjustable via ball-and-socket joint and vertical telescopic bar
- Side rail for mounting of further accessories

Technical information

Combination with

- table tops for traction or traction devices
- rotating and tilting clamp 1003.34A0 possible
- accessories to accommodate the foot
- traction stirrup clamp with rotation 1003.35A0 possible



1003.49B0
Foot plate mount

- For holding the healthy leg during traction positions with the help of further accessories

Technical information

Recommended combination with

- table tops for traction or traction devices
- foot plate or traction boot



1007.43A0
Screw tension device

- Easy execution of fine traction using a hand crank
- Adjustment of the leg by a ball and socket joint
- Precise setting of the leg using integrated rotation and tilt function
- Scale for good orientation while using the rotating function
- Ergonomic rotation possible thanks to star handle design
- Free rotation possible

Technical information

Recommended combination with

- traction bars 1007.40A0 and 1007.41A0 in combination with vertical adjustment 1007.51A0/B0
- accessories to accommodate the foot e.g. traction boot 1003.75A0 or foot plate 1001.97A0



1007.52A0
Bar extension

- To lengthen traction bar 1007.40A0

Technical information

Width x height: 248 x 278 mm

Combination with

- traction bar 1007.40A0 only



1007.51B0
Vertical adjustment, short

- For accommodation and vertical adjustment of a screw tension device
- For femoral fractures
- Easy to mount

Technical information

Length: 500 mm

Only in conjunction with

- screw tension device 1007.43A0

Recommended combination with

- traction bar 1007.41A0



1007.51A0
Vertical adjustment, long

- For accommodation and vertical adjustment of a screw tension device
- For hip arthroscopy and femoral fractures
- Easy to mount

Technical information

Length: 800 mm

Only in conjunction with

- screw tension device 1007.43A0

Recommended combination with

- traction bar 1007.41A0
- telescopic bar 1003.59V0



1003.59V0
Telescopic bar

- Combination with a traction device
- To be used instead of a telescopic bar in order to reach higher traction forces e.g. for hip arthroscopy

Technical information

Combination with

- 1007.51A0 necessary
- traction devices 1419.01X0, 1150.59AC, 1180.19A0
- adapter plate for traction 1433.41AC
- support plate for extension 1160.59AC

Suitable for use with

- 1150.20

Variant:

- Set of traction bars 1003.59C0



1004.86B0
Universal support

- For holding the healthy leg during extension positions
- Secure fixation using 2 fixture straps

Technical information

Including: 2 fixture straps

Material of pad: PUR

Spare part: strap 3100.7614

Recommended combination with

- screw tension device 1003.3700



1004.8800
Side rail clamp

- For mounting accessories to the side rail with the help of a setting clamp

Technical information

EU side rail

Mounting to standard metal traction bars



1007.45A0
Support plate

- To support the patient's leg
- Can be mounted to the sliding clamp 1007.44A0 using a clamp (e.g. 1003.23C0)

Technical information

Plate width x length x height: 252 x 290 x 465

Height adjustment: 465 mm


Material of pad: SFC

Recommended combination with

- sliding clamp 1007.44A0 and a standard or radial setting clamp

Suitable for use with

- 1433.02

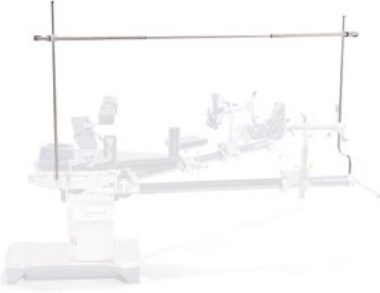

1007.48A0
Frame for sterile drape

- For holding sterile drape coverage
- Easy mounting to OR table
- Length of the horizontal bar is adjustable

Technical information

Comprises: 2 vertical bars and 1 horizontal bar

Adjustment range of horizontal bar: 1,500 – 2,800 mm

1 setting clamp may be ordered separately

Combination with

- traction bars 1007.40/41A0
- traction devices and orthopedic table tops possible



1007.49A0
Extension frame

- For holding a patient's leg in lateral and prone position
- Flexible positioning of the leg possible thanks to the adjustment range
- Easy to mount

Technical information

Comprises: 1 horizontal and 1 vertical bar

Only in conjunction with

- screw tension devices 1003.3700 and 1007.43A0

Recommended combination with

- foot plate 1001.97A0
- traction stirrup clamp with rotation 1003.35A0

**1003.3600**
Pelvis plaster cast device

- For supporting the pelvis during the application of a plaster cast
- Can be used for both children and adults
- Easy to apply plaster thanks to height-adjustable countertraction post

Technical information

Comprises: 2 removable pelvis rests, 1 adult, 1 child

Mounting to square mount of the OR table

Suitable for use with

- traction device 1180.19 and 1419
- OR table top 1150.20
- OR table 1425 and 1433.02
- support plate for extension 1160.59BC
- 1433.41AC

**1001.97A0**
Foot plate

- Secure fixation of the foot with hook and loop fastener via Velcro® straps
- Can be adapted to the width of the foot
- Pad to protect the heel

Technical information

Set of 2

Spare part: strap 9300.0943, rubber lining 9010.2444, heel pad 9010.2454

Recommended combination with

- rotating and tilting clamp 1003.34A0
- screw tension device 1003.3700
- foot plate mount 1003.49B0
- screw tension device 1007.43A0

**1003.75A0**
Traction boot

- Traction boot for positioning and securing of foot and calf during orthopedic procedures
- Shell is extremely robust
- 2 handles on each side for easy adjustment
- Size 37 – 45

Technical information

Comprises: 1 traction boot, 1 pack of 10 x 1003.76A0

Mounting to the side rail mounting point

Recommended combination with

- screw tension device 1003.3700, 1007.43A0
- foot plate mount 1003.49B0

**1003.76A0**
Disposable inlay and mounting brackets, big size

- Inlay for traction boot
- Skin-friendly material
- For boots size 37 – 45

Technical information

Packing unit: 1 pack of 10 inlays, each with 2 brackets

Suitable for use with

- traction boot 1003.75A0 only

Variant:

- 1003.76C0 for traction boot 1003.75C0

**1003.77A0**
Reusable inlay

- Reusable inlay for traction boot
- Skin-friendly material
- Size: 37 – 45

Technical information

Set of 2

Packing unit: 1 set of 2 pieces

Combination with

- disposable mounting brackets 1003.78A0

Suitable for use with

- traction boot 1003.75A0 only

**1003.78A0**
Disposable mounting brackets

- Stabilizes the patient's foot
- Quick and easy to mount

Technical information

Packing unit: 1 pack containing 15 pairs of brackets

Comprises: 1 pack containing 15 pairs of brackets

Recommended combination with

- reusable inlay 1003.77A0

Suitable for use with

- traction boot 1003.75A0 only

**1003.75C0**
Traction boot, small

- Traction boot for positioning and securing of foot and calf during orthopedic procedures
- Shell is extremely robust
- 2 handles on each side for easy adjustment
- For small size boots 32 – 37

Technical information

Comprises: 1 traction boot, 1 pack of 10 x 1003.76C0

Mounting to the side rail mounting point

Recommended combination with

- screw tension device 1003.3700, 1007.43A0
- foot plate mount 1003.49B0
- disposable inlay and brackets 1003.76C0

**1001.8700**
Foot plate, adult

- For holding the foot during extension positions

Technical information

Supplied individually

Only in conjunction with

- traction boot 1001.8800

Recommended combination with

- rotating and tilting clamp 1003.34A0
- screw tension device 1003.3700
- foot plate mount 1003.49B0
- screw tension device 1007.43A0



1001.8800
Traction boot, adult

- Secure fixation of the foot using 3 buckles on the leather cuff
- Pressure distribution thanks to interior padding

Technical information

Set of 2

Only in conjunction with

- foot plate 1001.8700



1001.9000
Foot plate, child

- For holding the foot during extension positions

Technical information

Supplied individually

Only in conjunction with

- traction boot for children 1001.9100

Recommended combination with

- rotating and tilting clamp 1003.34A0
- screw tension device 1003.3700
- foot plate mount 1003.49B0
- screw tension device 1007.43A0



1001.9100
Traction boot, child

- Secure fixation of the foot using 3 buckles on the leather cuff
- Pressure distribution thanks to interior padding

Technical information

Set of 2

Only in conjunction with

- foot plate 1001.9000



1001.92C0
Armpit support

- Radiolucent countertraction post
- Padded axilla support

Technical information

Material of pad: PUR

Can be mounted to a side rail

Integrated standard clamp

Only in conjunction with

- elbow support 1004.23C0



1004.23C0
Elbow support

- For repositioning of elbow fractures
- Range of movement adapted to the anatomy of the patient

Technical information

Can be mounted to a side rail

Integrated standard clamp

Combination with

- armpit support 1001.92C0
- Weinberger arm restraint cuff 1001.48A0



1004.9800
Humerus positioning device

- For repositioning of humerus fractures
- Supporting the implementation of intra-medullary nails for humerus fractures
- Axilla correct positioning of the forearm using inclined side rail

Technical information

Comprises: 1 standard clamp, 1 vertical bar, 1 adjustable horizontal bar with a side rail

Can be mounted to a side rail

Recommended combination with

- humerus countertraction post 1004.8000
- Weinberger arm restraint cuff 1001.48A0



1001.48A0
Weinberger arm restraint cuff

- Secure fixation of the hand

Technical information

Material: hand plate and finger fixation made of plastic

Can be sterilized

Recommended combination with

- humerus positioning device 1004.9800



1004.8000
Humerus counter traction post

- For fitting of intra-medullary nails for humerus fractures
- Radiolucent pad

Technical information

Can be mounted to a side rail

Clamp may be ordered separately:

- 1003.22C0 standard clamp
- 1003.23C0 radial setting clamp
- 1003.61A0 radial setting clamp

Recommended combination with

- humerus positioning device 1004.9800
- Weinberger arm restraint cuff 1001.48A0



1001.45C0
Upper arm posturing plate

- For holding the upper arm in prone position
- Radiolucent, good C-arm access
- May be fitted parallel or at a right-angle

Technical information

- Width x length: 210 x 260 mm
- Height-adjustable: 100 mm
- Material of pad: PUR
- Standard clamp supplied
- Can be mounted to a side rail



1002.30A0
Trimano Fortis

- Support arm can be used in various surgical disciplines
- Simple to operate with the Click and Move system
- Flexible positioning also during the operation
- Easy mounting of accessories
- Immediately ready for use without power supply or additional air connection

Technical information

- Sterile cover sets available
- Can be mounted to a side rail
- Combination with**
 - adapter Trimano Fortis 1002.31A0
 - arm support 1002.32A0



1002.31A0
Adapter Trimano Fortis

- Adapter for mounting accessories to Trimano Fortis 1002.30A0
- Can be sterilized

Technical information

- Combination with**
 - arm support 1002.43A0/B0
 - Trimano elbow countertraction post 1002.45A0
 - Trimano finger fixation 1002.48A0
 - disposable sterile arm rest 1002.52A0
- Variant:**
 - adapter 1002.42A0 for Trimano 3D (1002.41A0)



1002.48A0
Trimano finger fixation

- Secure fixation of the hand during surgery
- Mounting using Easy Click technology

Technical information

- Material: hand plate and finger fixation made of plastic
- Can be sterilized
- Sterile set 1002.44D0 may be ordered separately
- Suitable for use with**
 - adapter Trimano Fortis 1002.42A0 and 1002.31A0



1002.43B0
Arm rest for beach chair position

- Arm rest for Trimano Support Arm
- For supporting and securing the arm in beach chair position

Technical information

- Length: 270 mm
- Material: steel
- Can be sterilized
- Only in conjunction with**
 - adapter 1002.42A0 for Trimano 1002.41A0
 - adapter Trimano Fortis 1002.31A0 for Trimano 1002.30A0
 - sterile set 1002.44A0



1002.43A0
Arm rest for lateral position

- Accessory for Trimano Support Arm
- For lateral position

Technical information

- Length: 270 mm
- Material: steel
- Can be sterilized
- Only in conjunction with**
 - adapter 1002.42A0 for Trimano 1002.41A0
 - adapter Trimano Fortis 1002.31A0 for Trimano 1002.30A0
 - sterile set 1002.44A0



1002.45A0
Trimano elbow countertraction post

- Accessory for Trimano support arm
- For fitting of intra-medullary nails for humerus fractures
- Radiolucent
- Simple mounting using Easy Click technology

Technical information

- Including: padded roll with 60 mm and 90 mm diameter
- Length of countertraction post: 226 mm
- Material: PUR and carbon-fiber
- Sterile sheath 1002.44C0 may be ordered separately
- Combination with**
 - adapter Trimano Fortis 1002.42A0 and 1002.31A0



1002.52A0
Disposable arm rest

- Accessory for Trimano 3D and Trimano Fortis
- Light and secure fixation of the arm during surgery
- For single use

Technical information

- Comprises: disposable arm rest, bandage, sterile cover
- Packing unit: Set of 6 (comprises: disposable arm rest, bandage and sterile cover)

**1180.34A0**
Shoulder module

- Back plate for shoulder surgery
- Easy access to the shoulder region thanks to removable shoulder modules
- Head mount can be swivelled and inclined
- Inclination adjustment via motorized mounting point
- Pressure distribution thanks to SFC padding
- Side rails enable attachment of additional accessories

Technical information

Width x length: 582 x 760 / 870 mm

Interface: Magnus joint module

Including: 1 head rest adapter

EU side rail

Recommended combination with

- head support 1002.15A0
- trolley 1180.34T0

Suitable for use with

- motorized joint module 1180.11A0

**1180.34T0**
Trolley

- For storage and transportation of shoulder module 1180.34A0
- Easy mounting of the back plate to Magnus joint module
- Storage options for leg plates, extension plates and head rests

Technical informationWidth x height x length:
850 x 910 x 790 mm4 double swivel castors,
2 with brakes**1007.18B0**
Back plate for shoulder operations

- Good access to the shoulder region thanks to removable shoulder modules
- Head mount can be swivelled and inclined
- Inclination adjustment by means of a hand crank
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 575 x 550 mm

Including: 1 head rest adapter

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Recommended combination with

- head support 1002.15A0
- trolley 1007.20A0

Suitable for use with

- 1118, 1131, 1132, 1133, 1433, 7000, 7100, 7200, 7700

**1007.20A0**
Trolley

- For transportation and easy mounting of back plates for shoulder surgery
- Compact storage of back plates
- Attachment of a wire basket for storage of accessories possible

Technical informationWidth x height x length:
835 x 725 x 735 mmFor accommodating shoulder plate
1433.34AC order retrofit kit
3102.3089**Suitable for use with**

- shoulder plates
1007.18B0, 1007.19BC, 1007.21B0

May be ordered separately:

- wire basket 0802.0783

**1132.34A0**
Back plate for shoulder operations

- Removable shoulder modules enable good access to the shoulder region
- Head mount can be swivelled and inclined
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 530 x 290 mm

Including: 1 head rest adapter

Material of pad: SFC

Mounting to a pin mounting point

Recommended combination with

- head support 1002.15A0

Suitable for use with

- 1118.06/10, 1131, 1132, 1133, 1160.10/30, 1433, 7000, 7100, 7200, 7700

**1433.34AC**
**Carbon-fiber back plate
for shoulder operations**

- Removable shoulder modules enabling good access to shoulder region
- Side rails on shoulder modules can be removed
- Lateral supports can be flexibly adapted to patient's width
- Inclination adjustment via motorized mounting point
- Simple mounting of arm supports by way of additional side rails

Technical informationDimensions of the pad
width x length: 503 x 540 mmIncluding: 2 lateral supports,
2 side rails

Material of pad: SFC

Transportation and storage possible
on trolley 1007.20A0 with retrofit kit
3102.3089

EU side rail

Mounting to a pin mounting point

Recommended combination with

- head rest adapter 1002.10A0
- head support 1002.15A0

Suitable for use with

- 1131, 1132.21/22, 1133, 1160.10/30, 1433, 7100, 7200, 7700

**1007.19BC**
Back plate for shoulder operations

- Good access to the shoulder region thanks to removable shoulder modules
- Head mount can be swivelled and inclined
- Inclination adjustment via motorized mounting point
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 600 x 530 mm

Including: 1 head rest adapter

Material of pad: SFC

EU side rail

Mounting to a pin mounting point

Recommended combination with

- head support 1002.15A0
- trolley 1007.20A0

Suitable for use with

- 1150.19/25/30

**1007.21B0****Back plate for shoulder operations**

- Easy access to the shoulder region thanks to removable shoulder modules
- Head mount can be swivelled and inclined
- Inclination adjustment via motorized mounting point
- Pressure distribution thanks to SFC padding

Technical information

Width x length: 600 x 550 mm

Including: 1 head rest adapter

EU side rail

Mounting to a pin mounting point

Recommended combination with

- head support for shoulder operations 1002.15A0
- trolley 1007.20A0

Suitable for use with

- 1131
- 1132.21/22
- 1133
- 1433
- 7100, 7200, 7700
- 1160.10/30

**1009.02C0****Meniscus positioning device**

- Flexible positioning of the leg using the radial setting clamp
- Radiolucent padded roll

Technical information

Width: 367 mm

Overall height: 600 mm

Material of pad: PUR

Can be mounted to a side rail

Radial setting clamp 1003.23C0 may be ordered separately

**1004.84A0****Knee positioning device**

- Assists the surgeon during procedures on the knee
- Motorized adjustment operated via foot switch
- Secure fixation of the thigh
- Can be adapted to various thigh diameters
- Pressure relief thanks to extra soft gel pad

Technical information

Comprises: 1 positioning device, 1 footswitch, 5 straps

Voltage: 220 V

Spare parts: strap 3101.6183, gel pad 4006.1700

Can be mounted to a side rail

**1004.84B0****Knee positioning device**

- Secures the thigh during procedures on the knee
- Manual adjustment using the radial setting clamp
- Secure fixation of the thigh
- Can be adapted to various thigh diameters
- Pressure relief thanks to extra soft gel pad

Technical information

Comprises: 1 positioning device, 1 radial setting clamp, 1 gel pad, 1 strap

Spare parts: strap 3101.6183, gel pad 4006.1700

Can be mounted to a side rail

**1008.14C0****Arthroscopic leg holder**

- Secures the thigh during procedures on the knee
- Adjustable width for various diameters of the leg
- Secure fixation of the thigh through secure-lok strap
- Pressure relief thanks to disposable pad

Technical information

Width x height x length: 64 x 190 x 305 mm

Weight: 1.8 kg without clamp

Maximum patient weight: 227 kg

Material of device: stainless steel, Polycarbonate

Packing unit: 1 piece without inlay

Comprises: 1 piece with secure lock strap, 1 round post clamp, without disposable inlay

No mechanical cleaning

Can be mounted to a side rail

Integrated round post clamp (EU)

Recommended combination with

- 1008.14B0

**1008.14B0****Disposable inlay for arthroscopic leg holder**

- Disposable pad for arthroscopic leg holder
- Antimicrobial coating
- Non-sterile

Technical information

EU, US side rails

Material of pad: disposable

Packing unit: 1 unit/box = 12 pieces

Single use


Can be mounted to arthroscopic leg holder

Recommended combination with

- 1008.14C0 arthroscopic leg holder

07 Accessories for operating tables

Neurosurgery

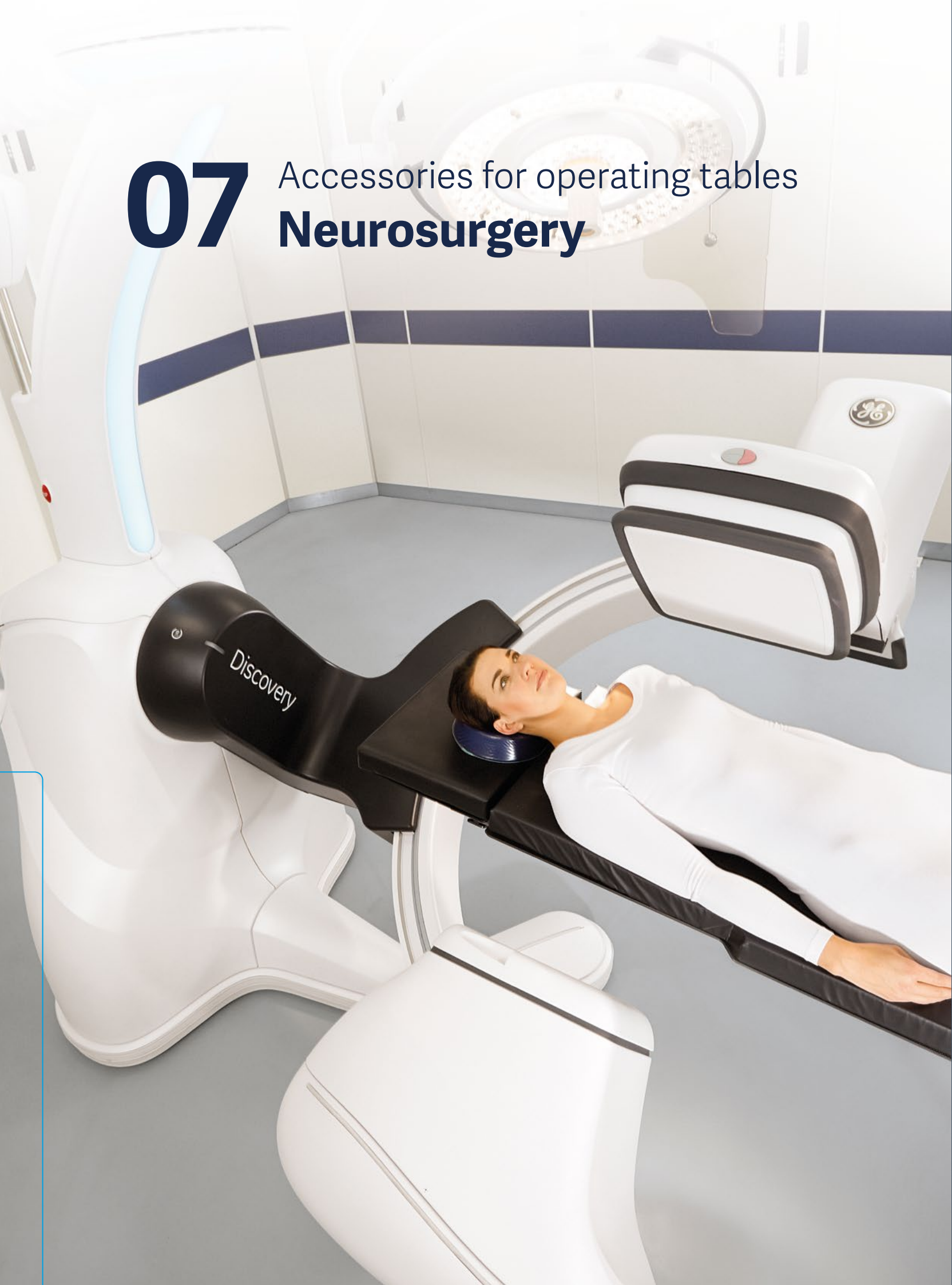

Stable and anatomically precise patient positioning is important for all surgical disciplines, none more so than neurosurgery because any unintentional movement of the OR table could influence the position of the patient and have fatal consequences.

Special accessories such as the cross-bar attachment 1005.2300 and securing functions like the Neuro-Lock, a feature of the Magnus OR Table System, provides a secure and precise positioning for neurosurgical interventions.

All Maquet OR Tables are designed to give the surgeon optimum access to the patient. The use of modern imaging procedures, for example, 3D X-rays, places increasing demands on the OR table accessories for neurosurgery. With specially designed and radiolucent accessories as well as high-quality pads, our products facilitate optimum X-ray results.



1005.50A0
DORO® Adjustable Base Unit

- The adjustable base unit is part of a system to position and secure the patient's head for head and neck surgery
- To fix a skull clamp
- For positioning a patient in prone and supine positions
- For use in non-imaging neuro-surgical applications
- Adjustable base unit including wrench for width adjustment

Technical information

Width x height: 547 x 213 mm

Weight: 3.2 kg

Material of device: aluminum

Mounting to round tube mounting point e.g. connection brackets 1130.54B0 and 1180.36A0 and cross-bar attachment 1005.2300

Combination with

- DORO® Swivel Adaptor 1005.51A0
- DORO® QR3 Skull Clamp 1005.52C0

Suitable for use with

- 1160.36A0 retrofit kit



1005.52C0
DORO® QR3 Skull Clamp

- Unique shape grants free space to place the intubation tube, even in prone position.
- No play between archholder and skull clamp base when locked.
- Minimum play between skull clamp base and extension arm ensures high stability and reduced risk to lost navigation.
- Easy to read pressure scale in metric and imperial units.
- Quick rail interface to securely attach retractor systems and navigation tracking devices
- With receptacles for skull pins for rigid cranial fixation

Technical information

Height : 284 mm

Weight: 1.8 kg

Material: light alloy and CrNi steel

Pins may be ordered separately: DORO® Reusable pins for adults 1005.53B0, DORO® Reusable pins for children 1005.53A0, DORO® Sterile Disposable Skull Pins, adult 1005.67A0

Recommended combination with

- reusable or disposable pins

Only in conjunction with

- DORO® Adjustable Base Unit 1005.50A0
- DORO® Swivel Adaptor 1005.51A0



1005.53B0
DORO® Reusable Skull Pins, adult

- Pins for fixing the patient's head in a skull clamp
- For multiple use
- Can be sterilized

Technical information

Material: stainless steel

Set of 3

Combination with

- skull clamps 1005.52C0 / 1005.49B0 / 1005.42A0 / 1005.48B0 / 1005.53X0 / 1005.67A0

Available versions:

- DORO® Reusable pins for children 1005.53A0



1005.67A0
DORO® Sterile Disposable Skull Pins, adult

- For single use
- For adults

Technical information

Set of 36

Material: stainless steel

Combination with

- skull clamps 1005.52C0 / 1005.49B0 / 1005.42A0 / 1005.48B0 / 1005.53X0 / 1005.67A0



1005.56A0
DORO®Swivel Horseshoe Head Rest, adult

- For anatomical head support
- With extension bar for traction on the spine
- Can be used in prone and supine positions
- Pressure relief thanks to extra-soft gel pad

Technical information

Material of pad: Gel

Only in conjunction with

- DORO® Adjustable Base Unit 1005.50A0
- DORO® Swivel Adaptor 1005.51A0



1005.51A0
DORO® Swivel Adaptor

- For connecting the adjustable base unit with a skull clamp
- Pivotal by 360°
- Flexible adjustment of the head rest system to the position of the patient

Technical information

Material of device: aluminum

Only in conjunction with

- DORO® QR3 Skull Clamp 1005.52C0
- DORO® Adjustable Base Unit 1005.50A0



- 1005.2300
Crossbar Attachment**
- Secure fixation of the skull clamp to patients in seated position
 - Can be adapted to the OR table width

Technical information
Can be mounted to a side rail
2 clamps may be ordered separately: 1003.22C0, 1003.23C0, 1003.61A0
Recommended combination with
• DORO® Adjustable Base Unit 1005.50A0



- 1005.49B0
DORO® Skull Clamp Radiolucent***
- With receptacles for skull pins for rigid cranial fixation
 - Metal-free radiolucent

**Phase-out planned in Q2 2026*

Technical information
Maximum patient weight: 250kg
Pins must be ordered separately, e.g. disposable pins 1005.67A0 or reusable pins 1005.53B0
Combination with
• DORO® Swivel Adaptor Radiolucent 1005.64A0
• DORO® Transitional Member 1005.63A0, 1005.63B0
• skull clamp holder 1005.48B0



- 1005.48B0
Skull clamp holder**
- For mounting a radiolucent skull clamp
 - Flexible adjustment via 3 joints
 - Metal-free radiolucent
 - Easy to mount using Easy Click technology

Technical information
Width x height x length: 399 x 200 x 360 mm
Maximum patient weight: 250kg
Adjustment range: +135°/-110° for joint 1, 360° all other joints
Mounting to push-fit head rest mounting point
Combination with
• carbon-fiber back plates 1189.37A0, 1433.33AC
• carbon-fiber insertable plates 1433.45BC, 1160.45BC
• skull clamp radiolucent 1005.49B0 / 1005.42A0 / 1005.48B0 / 1005.53X0 / 1005.67A0
• 1180.16A0/A1/A2/A3/A5 CF table tops with pin mounting point
• 1160.16A0, A2
Variant:
• 1005.48A0 for adaptation of skull clamp (Mayfield® Infinity XR2 or skull clamp A-2002 of company Integra)



- 1005.42A0
DORO® LUCENT skull clamp**
- Reduced x-ray artefacts
 - Reduced geometry – but increased head space
 - Unique release mechanism
 - Supportive markings
 - Designed for all patients (from adults to kids)
 - Automated cleaning
 - Pediatric set: 1005.42B0

Technical information
Combination with
• 1005.48B0 skull clamp holder
• 1005.64A0 DORO® swivel adaptor radiolucent
• 1005.72A0 DORO® interface skull clamp radiolucent
• 1005.64A0, 1005.48B0



- 1005.71A0
DORO® Table Adaptor Maquet, radiolucent**
- For connecting a DORO® head rest system for surgery on the head with the OR table
 - Metal-free radiolucent

Technical information
Mounting to dovetail mounting point
Combination with
• accessory adapter 6005.25A0
• carbon-fiber multi-purpose plate 1150.45BC
• further head positioning accessories



- 1005.63A0
DORO® Transitional Member Radiolucent, long**
- Long version
 - Combination with a DORO® Skull Clamp
 - Metal-free radiolucent

Technical information
Combination with
• DORO® Table Adaptor Maquet, radiolucent 1005.71A0
• DORO® Swivel Adaptor Radiolucent 1005.64A0
• DORO® Transitional Member Radiolucent, short 1005.63B0



- 1005.63B0
DORO® Transitional Member Radiolucent, short**
- Short version
 - Combination with a DORO® Skull Clamp
 - Metal-free radiolucent

Technical information
Combination with
• DORO® Table Adaptor Maquet, radiolucent 1005.71A0
• DORO® Swivel Adaptor Radiolucent 1005.64A0
• DORO® Transitional Member Radiolucent, long 1005.63A0



1005.64A0
DORO® Swivel Adaptor Radiolucent

- For mounting to a DORO® Skull Clamp
- Radiolucent

Technical information

Combination with

- DORO® Transitional Member Radiolucent 1005.63A0/B0
- Skull clamp radiolucent 1005.49B0 1005.42A0



1005.72A0
DORO® Interface Skull Clamp Radiolucent

- For mounting a radiolucent skull clamp to an aluminum swivel adapter
- Metal-free radiolucent

Technical information

Combination with

- DORO® Swivel Adaptor 1005.51A0
- Skull clamp radiolucent 1005.49B0 1005.42A0



1002.03A0
Horseshoe head rest

- For anatomically correct head support
- Can be used in prone and supine positions
- Metal-free radiolucent
- Pressure relief thanks to extra soft gel padding
- Individual adaptation to the patient's anatomy by using 3 rotating joints
- Mounting to push-fit head rest mounting point

Technical information

Width x height x length:
399 x 176 x 371 mm

Maximum patient weight: 250 kg

Combination with

- CF back plate 1180.37A0
- CF back plate 1433.33AC
- carbon-fiber insertable plates 1160.45BC and 1433.45BC

Suitable for use with

- carbon-fiber table tops 1180.16A1/A2/A3/A5
- 1160.16A0, A2



1007.12A0
Guiding roller for head-side traction

- Optimum patient positioning due to length and height adjustment

Technical information

Mounting to the side rail of the head rest adapter

Combination with

- horseshoe head rest 1002.72A0
- fixture 1002.65A0



1002.83A0
Head rest that can be angled

- Metal-free radiolucent
- Easy mounting with Easy Click technology

Technical information

Width x height x length:
455 x 155 x 324 mm

Product weight: 2.8 kg

Maximum patient weight: 250 kg

Height: max. 155 mm

Length: max. 428 mm

Inclination: +90° / - 45°

Material of pad: SFC

Mounting to push-fit head rest mounting point

Recommended combination with

- CF back plates 1180.37AC, 1433.33AC
- carbon-fiber insertable plates 1433.45BC, 1160.45BC

Suitable for

- carbon-fiber table tops 1180.16A0/A1/A2/A3/A5
- 1160.16A0, A2



1002.82A0
Rigid head rest

- Metal-free radiolucent
- Mounting with Easy Click technology

Technical information

Maximum patient weight: 250 kg

Product weight: 2.1 kg

Width x height x length:
420 x 155 x 303 mm

Material of pad: SFC

Mounting to push-fit head rest mounting point

Recommended combination with

- CF back plates 1180.37A0, 1433.33AC
- CF insertable plates 1433.45BC, 1160.45BC

Suitable for use with

- carbon-fiber table tops 1180.16A0/A1/A3/A5 with pin mounting point
- 1160.16A0, A2
- Compatible with CF table top 1180.16A2



1002.04B0
Horseshoe head rest

- For anatomically correct head support
- Metal-free radiolucent
- Pressure relief thanks to extra soft gel padding
- Individual adaptation to the patient's anatomy by using 2 rotating joints

Technical information

Maximum patient weight: 200 kg

Large adjustment angle ± 50°

Mounting to dovetail mounting point

Combination with

- 1180.37A0 and 1433.33BC with accessory adapter 6005.25A0
- carbon-fiber multi-purpose plates 1160.45BC and 1433.45BC with accessory adapter 6005.25A0
- 1180.45BC
- 1150.30BB0/D0 with carbon-fiber multi-purpose plate 1150.45BC

Suitable for use with

- 1150.16B0



1007.24A1
Universal frame for patient positioning

- For prone, supine and lateral positions
- 360° radiolucent
- Mirror to check the intubation position
- Easy access to surgical field thanks to slim construction
- Slim frame construction enables good X-ray access
- Simple transportation and mounting using the trolley

Technical information

Height adjustable: 600 – 900 mm

Including: 1 trolley, 1 pair of adapters, 1 head rest adapter, 1 mirror, 2 arm supports, 1 Plexiglas head rest, 1 pair of adjustable side rails, 1 side rail clamp, 1 filler piece

Square mount for head positioning accessories

Round tube mount for standard holder

Mounting to a pin mounting point

Suitable for use with

- 1118.01B0/K0, 1118.06B0/K0, 1118.10B0/ K0
- 1133.22B4/B5
- Meera: 7000.01B0, 7100,01B0/B2, 7200.01B0/B2
- 1150.30B0/D0
- 1180.10A0: adapter 1007.27A0 has to be ordered separately
- 7700.01A0, B0, A2, B2
- 1160.16/30
- Maquet Yuno II: 1433.02B0



6005.25A0
Accessory adapter

- For mounting accessories with a dovetail mounting point to a carbon-fiber plate
- Metal-free radiolucent
- Suitable for use in MRI

Technical information

Maximum patient weight: 180 kg

Including: 2 screws

Material: carbon fiber

Combination possible with

- carbon-fiber insertable plates 1433.45BC, 1160,45BC
- CF multi-purpose plate 1180.45A0
- CF back plate 1180.37A0
- CF back plate 1433.33AC



1007.27A0
Adapter for universal frame for Magnus

- For mounting the universal frame for patient positioning 1007.24A1 to the Magnus OR Table System

Technical information

1 pair

Only in conjunction with

- motorized joint module 1180.11A0



1007.26A0
Short support plate

- For supine and lateral positions
- Good access to the head area
- Metal-free radiolucent

Technical information

Length: 850 mm

Material of pad: SFC

Combination with

- universal frame for patient positioning 1007.24A1



1007.26B0
Long support plate

- For prone, supine and lateral positions
- Enables very large imaging window
- Metal-free radiolucent

Technical information

Length: 1,350 mm

Material of pad: SFC

Combination with

- universal frame for patient positioning 1007.24A1



1007.28B0
Adjustable positioning elements

- For positioning the patient in prone position on the universal frame
- Can be adjusted in three planes: vertically, horizontally, laterally
- Small body contact area enables free X-ray channel to the spine

Technical information

Comprises: 2 chest elements, 2 hip elements, 1 chest pad

Material of pad: SFC

Recommended combination with

- disposable covers 1007.29A0

Combination with

- universal frame for patient positioning 1007.24A1



1007.25A0
Side rail clamp

- For the adaptation of additional accessories to the side rail mounting point

Technical information

Can be mounted to a bar

EU side rail

Combination with

- universal frame for patient positioning 1007.24A1



1007.30A0
Adapter for radiolucent spine surgery frame

- Enables the adaptation of radiolucent spine surgery frame
- Metal-free radiolucent

Technical information

Combination with

- universal frame for patient positioning 1007.24A1



1008.01A0
KTEK Curve Spine Frame

- Especially suitable for spinal surgery
- Curve enables necessary space between vertebrae
- Curve can be adjusted by using the hand crank
- Width of curve elements can be adjusted
- Radiolucent

Technical information

Width x height x length: 444 x 254 x 720 mm

Weight: 11 kg

Maximum patient weight: 272 kg

Hand crank included



1008.02A0
KTEK Curve Cart

- Storage and transport of the KTEK Curve Spine Frame
- Storage bin for the hand crank in the upper area of the cart
- Small storage bin for accessories in the base segment

Technical information

Width x height x length: 549 x 1035 x 549 mm

Weight: 5 kg

May be ordered separately

Combination with

- 1008.01A0 KTEK Curve Spine Frame



1008.03A0
Curve comfort covers

- Disposable comfort covers for use with curve spine frame

Technical information

Packing unit: 1 unit = 6 pieces

Comprises: cosmetic foam pad covers, head block, lateral rod pad, arm cushions

Disposable

- Combination with 1008.01A0 KTEK Curve Spine Frame



1008.03B0
Curve standard cover

- Disposable standard covers for use with curve spine frame

Technical information

Packing unit: 1 unit = 6 pieces

Comprises: standard foam pad covers, head block, lateral rod pad, arm cushions

Disposable

- Combination with 1008.01A0 KTEK Curve Spine Frame



1008.03C0
Head block

- Disposable head block to position the head in prone position

Technical information

Material of pad: foam

Packing unit: 1 unit = 6 pieces

Comprises: 1 head block

Disposable

- Combination with 1008.01A0 KTEK Curve Spine Frame



- 1007.04D0**
Kneeling frame
- For positioning the patient in kneeling position
 - Flexible positioning using height-adjustable calf support
 - Manual adjustment of calf support using a hand crank
 - Mounting to OR table is made comfortable by using a trolley
 - Mounting to Magnus joint module

Technical information
Maximum patient weight: 135 kg/ in combination with 1180.10: 250kg
Interface: Magnus joint module 1180.11A0
Material of pad: SFC
EU side rail
Recommended combination with • trolley 1007.07D0 • sitting bracket with pelvis supports 1007.05A0
Available versions: • 1007.04BC suitable for 1131.12B0, 1133.22B4/B5, 7100.01B0/B2, 7200.01B0/B2, 1160.10A0/B0/D0, 1160.30A0/B0/D0, 7700.01A0/B0/A2/B2 • 1007.04D0 suitable for 1180.10A0



- 1007.07A0**
Trolley
- For support and transportation of kneeling frame
 - Facilitates back-friendly mounting of the kneeling frame

Technical information
Including: 1 wire basket
Available versions: • 1007.07A0 for kneeling frame 1007.04BC • 1007.07D0 for kneeling frame 1007.04D0



- 1007.54A0**
Trolley
- For support and transportation
 - Enables easy mounting of kneeling frame 1007.53A0

Technical information
Combination with • kneeling frame 1007.53A0



- 1007.53A0**
Kneeling frame
- For positioning a patient in knee-elbow position
 - Height-adjustable calf support
 - Manual adjustment of calf support using a hand crank
 - Mounting to OR table is made comfortable by using the trolley

Technical information
Material of pad: SFC
EU side rail
Mounting to a pin mounting point
Recommended combination with • trolley 1007.54A0 • sitting bracket with pelvis supports 1007.05A0
Suitable for use with • 1118.01B0/K0 • 1118.06B0/K0 • 1133.22B4/B5 • 7000.01B0 • 7100.01B0/B2 • 7200.01B0/B2 • 1160.10A0/B0/D0 • 1160.30A0/B0/D0 • 1433.02B0 • 7700.01A0, B0, A2, B2



- 1007.05A0**
Sitting bracket with pelvis support
- For buttocks support during surgery
 - Lateral supports secure the thighs
 - Width can be adjusted

Technical information
Comprises: 1 sitting bracket with padded roll, 2 lateral pelvis supports
Material of pads: PUR
2 setting clamps may be ordered separately: 1003.22C0, 1003.23C0
Combination with • kneeling frame 1007.04D0/1007.04BC and 1007.53A0



- 1004.75A0**
Wrist support
- For supporting the surgeon's hands and arms during surgery
 - Wrist support can be flipped 180°
 - Height-adjustable

Technical information
Width x height: 355 x 240 mm
Including: 1 standard clamp
Mounting to a side rail mounting point
Recommended combination with • head rest adapter 1002.07A0 • fixture 1002.65A0 • head rest adapter 1007.16A0/B0

08

Accessories for operating tables

Positioning aids



Positioning injuries are potential complications of surgical operations. Positioning aids can reduce this risk and help to prevent decubitus during the surgery and on the ward.

In addition to the medical aspect, positioning aids are also ideal for making patients more comfortable on the OR table. When required to lie in a supine position for a long period, patients often complain of backache. Gently raising the legs with the help of a pad, for example, the leg rest 1001.94A0, will take the pressure off the spine.

The positioning aids listed here are simple to use and easy to clean. They are high quality and very durable. For all positions, there are complementary gel cushions or pads in various sizes.

Head rests



4006.01A0
Oasis Elite Closed Head Ring Adolescent

- To support and position the patient's head
- For adolescents, closed version

Technical information

Height: 34 mm

Ø 140 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.02A0
Oasis Elite Closed Head Ring Adult

- To support and position the patient's head
- For adults, closed version

Technical information

Height: 50 mm

Ø 200 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.03A0
Oasis Elite Open Head Ring Adolescent

- To support and position the patient's head
- For adolescents, small version
- Open on one side

Technical information

Height: 34 mm

Ø 140 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.04A0
Oasis Elite Open Head Ring Adult

- To support and position the patient's head
- For adults
- Open on one side

Technical information

Height: 54 mm

Ø 200 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.19A0
Oasis Elite Prone Head Rest, large

- For adults
- For support, positioning and protection of the head in prone position
- Lateral cut-outs offer space for ventilation tubes

Technical information

Width x height x length:
240 x 140 x 280 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.20A0
Oasis Elite Prone Head Rest, medium

- For adolescents
- For support, positioning and protection of the head in prone position
- Lateral cut-outs offer space for respiration tubes

Technical information

Width x height x length:
273 x 135 x 220 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet




4006.21A0
Oasis Elite Supine Head Rest Adult

- For adults
- For support, positioning and protection of the head in supine position

Technical information

Width x height x length:
237 x 90 x 275 mm

Material: Silicone gel

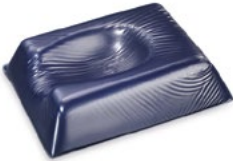

4006.22A0
Oasis Elite Supine Head Rest Adolescent

- For adolescents
- For support, positioning and protection of the head in supine position

Technical information

Width x height x length:
230 x 65 x 275 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.23A0
Azure Head Protector

- To support and position the patient

Technical information

Width x height x length:
250 x 50 x 260 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



1000.5700
Head ring

- For placement and positioning of the patient's head

Technical information

Height: 40 mm

Ø 210 mm

Material: PUR

Operating table pads



4006.05A0
Universal operating table pad

- To cover the flat surface of the operating table
- To ensure even distribution of patient weight

Technical information

Width x height x length:
520 x 10 x 520 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.06A0
Oasis ¾ Length Table Pad

- Long version
- To cover the flat surface of the operating table
- To ensure even distribution of patient weight

Technical information

Width x height x length
520 x 10 x 1,150 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.07A0
Oasis Full Length Table Pad

- Extra-long version
- To cover the flat surface of the operating table
- To ensure even distribution of patient weight

Technical information

Width x height x length:
520 x 10 x 1,800 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.08A0
Oasis Torso Section

- To cover a specific area of the operating table
- With sacral cut-out

Technical information

Width x height x length:
500 x 10 x 990 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.09A0
Oasis Perineal Table Pad

- To cover a specific area of the operating table
- With sacral cut-out

Technical information

Width x height x length:
520 x 10 x 520 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.24A0
Azure Sacral Protector, small

- To support the patient's sacral region
- Ensuring more even weight distribution

Technical information

Width x height x length:
380 x 25 x 400 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet

Arm supports



4006.12A0
Oasis Armboard Pad

- Small version
- For placement and support of the extended arm of a patient

Technical information

Width x height x length:
115 x 10 x 395 mm

Material: Silicone gel



4006.13A0
Oasis Armboard Pad

- Large version
- For placement and support of the extended arm of a patient

Technical information

Width x height x length:
200 x 10 x 600 mm

Material: Silicone gel

Recommended combination with

- arm posturing device 1001.44D0



4006.26A0
Azure Arm Protector, short

- Small version
- For placement and support of the patient's arm

Technical information

Width x height x length:
140 x 30 x 330 mm

Material: Silicone gel



4006.27A0
Azure Arm Protector, long

- Large version
- For placement and support of the patient's arm

Technical information

Width x height x length:
180 x 35 x 460 mm

Material: Silicone gel

Recommended combination with

- arm posturing device 1001.44D0

Leg rests and heel protector pads



4006.14A0
Oasis Crutch Stirrup Pads

- Additional padding to increase positioning comfort

Technical information

Width x height x length:
220 x 10 x 320 mm

Including: 2 pieces

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.15A0
Oasis Crutch Stirrup Pads

- Additional padding to increase positioning comfort

Technical information

Width x height x length:
245 x 10 x 395 mm

Set of 2

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.17A0
Oasis Thigh Pad

- To relieve pressure

Technical information

Width x height x length:
130 x 10 x 520 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.28A0
Azure Leg Protector

- For placement and support of the patient's leg
- Ensures more even distribution of patient weight

Technical information

Width x height x length:
180 x 40 x 515 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.16A0
Oasis Elite Heel Pads

- To support the leg at the ankle

Technical information

Width x height x length:
100 x 70 x 180 mm

Set of 2

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.29A0
Azure Heel Protectors

- For placement and support of the heel

Technical information

Width x height x length:
175 x 30 x 190 mm

Set of 2

Material: Silicone gel

Suitable for

- warming in a warming cabinet

Leg rests and heel protector pads



4006.10A0
Oasis Universal Patient Positioner

- Medium version
- To support the patient’s position

Technical information

Width x height x length:
100 x 60 x 290 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.11A0
Oasis Universal Patient Positioner

- Large version
- To support the patient’s position

Technical information

Width x height x length:
120 x 70 x 345 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



4006.25A0
Oasis Elite Flat Bottom Chest Roll

- To support the patient’s position

Technical information

Width x height x length:
140 x 110 x 505 mm

Material: Silicone gel

Suitable for

- warming in a warming cabinet



1003.7400
Pad

- Semicircular
- For placement and positioning of the patient

Technical information

Width x height x length:
190 x 115 x 478 mm

Material: PUR



4006.18A0
Oasis Elite Lateral Positioner

- To support the patient lying on their side on the OR table
- To prevent the occurrence of nerve damage and pressure points

Technical information

Width x length: 500 x 720 mm, height up to 160 mm

Material: Silicone gel



4006.3200
Lateral positioning pad

- For placement and support of the patient in lateral position

Technical information

Width x height x length:
520 x 230 x 413 mm

Material: Silicone gel / SFC

Fixation to side rails of the operating table possible with the help of Velcro® straps



1000.77A0
Tunnel cushion

- For placement and positioning of the patient’s legs, especially in lateral position

Technical information

Width x height x length:
400 x 220 x 650 mm

Material of pad: SFC



1000.6900
Plexus cushion

- For positioning the patient, especially supporting the neck and shoulder region in supine position

Technical information

Width x height x length:
500 x 100 x 440 mm

Material of pad: SFC



**1001.94A0
Legrest**

- For positioning the patient, in particular the patient's legs

Technical information

Width x height x length:
450 x 200 x 600 mm

Material of pad: SFC



**1000.5600
Pad**

- For positioning the patient, especially during surgical procedures on the spine

Technical information

Outer dimensions
width x height x length:
460 x 100 x 840 mm

Abdominal cut-out width x length:
200 x 370 mm

Material of pad: SFC



**4006.3100
Prone positioning pad**

- For positioning the patient, especially during surgical procedures on the spine

Technical information

Dimensions cut-out:

Width chest: 250 mm
Width pelvis: 190 mm
Length: 580 mm

Outer dimensions:

Width x length x height:
480 x 1,014 x 140 mm

Material of pad: Special foam



**1000.68C0
Pad**

- For positioning the patient, especially during surgical procedures on the spine

Technical information

Thorax pad width x height x length:
500 x 160 x 250 mm

Pelvis side width x height x length:
500 x 160 x 340 mm

Material of pads: SFC



**1000.24A0
Padded roll for countertraction post**

- To pad a countertraction post for positioning the patient

Technical information

Ø (external) 180 mm

Height: 190 mm

Material: SFC

Combination with

- countertraction post, e.g. 1004.90B0



**1000.78A0
Vacuum mattress**

- Microspace filling adapts smoothly to the body shape
- With sacral cut-out
- Latex-free
- Radiolucent
- MRI and CT compatible

Technical information

Width x length: 750 x 1,250 mm

Length between shoulder and sacral cut-out: 1,000 mm

Filling: Microspace System

Cover: PUR

Non-slip foil 1101.6864 available



09

Accessories for operating tables **Medical furniture for the operating room**

The Resist range of medical furniture consists of reliable and well proven accessories. Instrument tables, infusions stands, step stools and many more items are essential for the daily work of the OR team.

Important technical features are hygienic design, corrosion resistance, durability, tip-proof design, resistance to disinfectants and robustness. The dirt-repellent stainless-steel surfaces are designed to maximize hygiene.



4403.9302
Instrument stand

- Instrument table for placing medical instruments during surgery
- Removable non-pivoting shelf
- Manual height adjustment

Technical information

- Shelf dimensions width x length: 600 x 400 mm
- Height adjustment: 800 – 1,200 mm
- Max. load: 30 kg
- Material: CrNi steel
- 4 double swivel castors (Ø 75 mm), 2 castors with brakes



4464.98A0
Mayo instruments table

- Instrument table designed to hold sterile goods containers, instrument trays, instruments, swabs and other small items required for the operation
- Deep-drawn seamless shelf
- Convenient height adjustment via foot pedal

Technical information

- Hydraulic height adjustment: 990 – 1,355 mm
- Max. load: 30 kg
- 3 double swivel castors (Ø 75 mm)
- Available versions:**
 - 4464.98A0, shelf dimension: 600 x 400 mm
 - 4464.98B0, shelf dimension: 700 x 500 mm



4464.9X00
Instrument stand

- Instrument table to hold containers for sterile goods, instrument trays, instruments, swabs and other small items required for the operation
- Deep-drawn seamless shelf
- Table top can be swivelled from +90 to -90°
- Table top fixation with a lever below the shelf
- Convenient height adjustment via foot pedal
- Low foot segment

Technical information

- Hydraulic height adjustment: 960 – 1,355 mm
- Max. load: 30 kg
- 3 double swivel castors (Ø 50 mm)
- Available versions:**
 - 4464.9500, shelf dimensions: 600 x 400 mm
 - 4464.9600, shelf dimensions: 700 x 500 mm



4809.9003
Bowl stand

- For collecting used consumables, for example medications, liquids, towels and similar products
- 2 removable bowls
- Handle for easy maneuvering

Technical information

- Height: 825 mm
- Base Ø: 635 mm
- 2 bowls (6 l each)
- Max. load per bowl: 4 kg
- Material: CrNi steel
- 5 double swivel castors (Ø 50 mm)



4404.32X3
Instrument table

- Hygienic accessory table with 2 removable shelves
- Easy to maneuver in the OR
- Resistant to disinfectants
- Finely polished and dirt-repellent surface

Technical information

- Shelf dimensions width x length: 625 x 425 mm
- 2 removable shelves
- Max. load per shelf: 20 kg
- 4 double swivel castors, 2 of which with brakes
- Available versions:**
 - 4404.32A3, height 800 mm, distance between shelves: 500 mm
 - 4404.32B3, height 900 mm, distance between shelves: 600 mm
 - 4404.32C3, height 1,000 mm, distance between shelves: 700 mm



4405.32X3
Instrument table

- Hygienic accessory table with 2 removable shelves
- Easy to maneuver in the OR
- Resistant to disinfectants
- Finely polished and dirt-repellent surface

Technical information

- Shelf dimension width x length: 720 x 625 mm
- 2 removable shelves
- Max. load per shelf: 20 kg
- 4 double swivel castors, 2 of which with brakes
- Available versions:**
 - 4405.32A3, height 800 mm, distance between shelves: 500 mm
 - 4405.32B3, height 900 mm, distance between shelves: 600 mm
 - 4405.32C3, height 1,000 mm, distance between shelves: 700 mm



4406.32X3
Instrument table

- Hygienic accessory table with 2 removable shelves
- Easy to maneuver in the OR
- Resistant to disinfectants
- Finely polished and dirt-repellent surface

Technical information

Shelf dimensions width x length:
940 x 625 mm

2 removable shelves

Max. load per shelf: 20 kg

4 double swivel castors,
2 of which with brakes

Available versions:

- 4406.32A3, height 800 mm, distance between shelves 500 mm
- 4406.32B3, height 900 mm, distance between shelves 600 mm
- 4406.32C3, height 1,000 mm, distance between shelves 700 mm



4402.9002
Hand and arm positioning table

- For positioning the patient's hand and arm
- Table top pivotable $\pm 45^\circ$
- Table top can be tilted $\pm 29^\circ$
- Plate is radiolucent
- Can be fixed in a comfortable position for patient and surgeon

Technical information

Dimensions of table top
width x length: 450 x 450 mm

Height adjustment: 650 – 950 mm

Max. load: 25 kg

4 castors ($\varnothing$ 75 mm),
all with brakes

Only in combination with

- SFC pad 4000.05A0



4000.05A0
Pad

- SFC pad is used in combination with hand and arm posturing table 4402.9002
- Radiolucent

Technical information

Width x height x length:
450 x 30 x 450 mm

Velcro® strap below the pad for
attachment to hand and arm table



4902.3503
Disposal trolley

- Mobile unit for collecting waste or fluids in the operating theatre
- Removable bucket

Technical information

Height to upper edge of bucket:
465 mm

Kick-ring base: $\varnothing$ 466 mm

Capacity of bucket: 15 liter

Max. load: 20 kg

Bucket: CrNi steel

4 double swivel castors

Available versions:

- 4901.3503
trolley without push bar
- 4902.3503
trolley with push bar,
height: 800 mm



4906.9003
OR step stool

- OR step stool for better access to the surgical field
- 1 step with slip-proof ridged tread
- Step seamlessly welded

Technical information

Dimensions width x length:
603 x 350 mm

Max. load: 125 kg

Material: CrNi steel

Electrostatically conductive
rubber feet

Available versions:

- 4906.9003, height: 130 mm
- 4906.9103, height: 220 mm



4907.9003
OR step stool

- OR step stool for better access to the surgical field
- 2 steps with slip-proof ridged tread
- Steps seamlessly welded

Technical information

Dimensions width x length:
605 x 650 mm

Height: 220 / 440 mm

Max. load: 100 kg

Material: CrNi steel

Electrostatically conductive
rubber feet



**4805.91A3
IV pole MRI**

- Infusion stand suitable for use in MRI areas up to 3 Tesla
- Minimum distance from MRI: 50 cm
- IV pole consisting of antimagnetic parts only
- 4 bottle hooks
- 1 drip glass
- Suitable for gravity infusions

Technical information

- Height adjustment 1,300 – 2,090 mm
- Max. load: 2 kg per bottle hook
- Material: aluminum and plastic
- Castors are antistatic
- 5 double swivel castors (Ø 50 mm) with brakes



**4805.91B3
Infusion stand**

- IV pole suitable for support of infusion bags, bottles and 1 pump
- One-hand height adjustment is possible
- 4 bottle hooks
- 1 drip glass
- Suitable to support gravity infusions

Technical information

- Height adjustment: 1,350 – 2,350 mm
- Max. load: 2 kg per bottle hook
- Material: CrNi steel
- Castors antistatic
- 5 castors (Ø 80 mm) with brakes



**4001.59A0
Bowl**

- 5 liter bowl with interior scale for use with the heated bowl holder
- Can be sterilized

Technical information

- Ø 320 mm
- Height: 120 mm
- Material: CrNi steel
- Mechanical cleaning is possible
- Recommended combination with**
 - heated bowl holder 4904.95/96A0
 - clamp ring 4001.60A0



**4001.60A0
Clamp ring**

- Clamp ring for fixation of sterile covers in combination with a bowl and the heated bowl holder

Technical information

- Ø 305 mm
- Height: 70 mm
- Material: CrNi steel
- Mechanical cleaning is possible
- Recommended combination with**
 - heated bowl holder 4904.95/96A0
 - bowl 4001.59A0



**4904.95A0
Heated bowl holder**

- Warming of physiological saline solutions
- Dampening or warming compresses and towels
- Temperature range 35 – 60 °C set by button
- Temperature monitoring on LED display
- Safety thermostat

Technical information

- Height adjustment: 1,050 – 1,350 mm
- Overall width: 650 mm
- Cable length: 5,700 mm
- Max. load: 6 kg
- Material: CrNi steel
- Teflon-coated heating plate
- 5 double swivel castors (Ø 75 mm), 3 of which electrically conductive, 2 of which with brakes
- Recommended combination with**
 - bowl 4001.59A0
 - clamp ring 4001.60A0

Variant:

- 4904.96A0 with mains plug for Switzerland



**4632.01A0
OR swivel stool pneumatic**

- Mobile swivel stool for flexible use in the operating theatre
- Round padded seat
- Covered castors
- Central brake system
- Foot pedal for convenient height-adjustment

Technical information

- Pneumatic height-adjustment: 570 – 740 mm
- Padded seat: Ø 350 mm, black
- Pad thickness: 100 mm
- Max. load: 200 kg
- Not electrically conductive
- 5 double castors (Ø 50 mm)
- Available versions:**
 - 4632.01A0 without back rest
 - 4632.01E0 with padded back rest



4634.01A0
OR swivel stool hydraulic

- Mobile swivel stool for flexible use in the operating theatre
- Round padded seat
- Covered castors
- Central brake system
- Foot pedal for convenient height-adjustment

Technical information

Hydraulic height-adjustment:
600 – 840 mm

Padded seat: Ø 350 mm, black

Pad thickness: 100 mm

Max. load: 200 kg

5 double castors (Ø 50 mm)

Electrically conductive

Available versions:

- 4634.01A0 without back rest
- 4634.01C0 with padded back rest



4602.3402
OR swivel stool

- Mobile swivel stool for flexible use in the operating theatre
- Round padded seat
- 5 double castors for stability

Technical information

Height adjustment manual:
460 – 670 mm

Padded seat: Ø 330 mm

Pad thickness: 70 mm

Max. load: 200 kg

Electrically conductive

5 double swivel castors (Ø 50 mm)

10

Overview

Accessories matrix and index



The most important information at a glance: Use the accessories matrix to quickly and easily find out which accessory can be used with any given OR table.

The accessories are sorted by article number, and symbols are used to identify the most important information. Firstly, the OR tables that the accessory can be used with and, secondly, what type of mounting point it has. You can also see what page the accessory is on.

If you do not know the article number, the index will help. It lists all the accessories in the catalog alphabetically. At the end of the chapter there is a page for making notes and an order form.

Accessories matrix and index

for easy orientation

Explanation of
the accessory matrix

		Page	Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0	Alphamaquet 1150.30B0/D0	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0
0802.0783	Wire basket	32	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b
1000.24A0	Padded roll for countertraction post	121	e	e	e		e	e	e	e	e	e	e	e		e	e	

Article number

Article name

Combination option

Name of OR table

a	Positioning aid
b	Accessory trolley
c	Medical furniture
d	Can be fitted to OR table mounting point
e	May be only used on the OR table mounting point in conjunction with other accessories
f	Accessories for side rails
g	Side rail accessories, only suitable for use in combination with other accessories

		Page	Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0*	Alphamaquet 1150.30B0/D0*	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0
0802.0783	Wire basket	32	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b
1000.24A0	Padded roll for countertraction post	137	e		e	e			e	e	e	e	e	e	e	e	e	e
1000.5600	Pad	136	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1000.5700	Head ring	128	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1000.68C0	Pad	136	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1000.6900	Plexus cushion	135	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1000.77A0	Tunnel cushion	135	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1000.78A0	Vacuum mattress	137	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1001.40A0	Arm protector ¹	53		a	a		a		a	a	a	a	a	a	a	a	a	a
1001.40B0	Arm protector ² , curved	53		a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1001.42A0	Arm support ³	53	e	d	e	e	g	d	e	e	a	e	e	e	e	e	e	e
1001.44E0	Arm support	47	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f
1001.44F0	Arm support	47	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f
1001.44D0	Arm positioning device PUR	47	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f
1001.45C0	Upper arm posturing plate	104	f		f	f			f	f	f	f	f	f	f	f	f	f
1001.4600	Arm restraint cuff	51	g	g	g	g	g	g	g	g	g	g	g	g	g	g	g	g
1001.4700	Leg restraint cuff	84	a		a	a			a	a	a	a	a	a	a	a	a	a
1001.48A0	Weinberger arm restraint cuff	103	e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	e
1001.54A0	Side rail extension	28	f		f	f			f	f		f	f	f	f	f	f	f
1001.56A0	Body strap	29	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1001.56B0	Body strap	29	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1001.57A0	Leg strap	84	a		a	a			a	a	a	a	a	a	a	a	a	a
1001.58A0	Body strap	30	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a
1001.59A0	Body strap	29	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f
1001.59B0	Body strap	30	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f
1001.6000	Arm support	48	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f
1001.63A0	Multi-purpose plate for arm surgery	51	e	g	e	e	g	g	e	e	e	e	e	e	e	e	e	e
1001.64X0	Adapter for CF table top	52	e	d			d											
1001.64C0	Adapter for CF table top (not for 1160.16A2)	52			e	e	d	d			e	e	e	e	e	e		e
1001.65A0	Leg holder	73	g		g	g			g	g	g	g	g	g	g	g	g	g
1001.66A0	Leg support	83	g		g	g			g	g	g	g	g	g	g	g	g	g

¹ suitable 1160.10/30 with SFC, 7.700.01B0/B2

² suitable 1160.10/30 with SFC, 7.700.01B0/B3

³ suitable 1160.10/30 with SFC, 7.700.01B0/B4

*Previously sold product

[illegible]

*Previously sold product

[illegible]

*Previously sold product

[illegible]

*Previously sold product

a b c d e f g Positioning aid Accessory trolley Medical furniture Can be fitted to OR table mounting point May be only used on the OR table mounting point in conjunction with other accessories Accessories for side rails Side rail accessories, only suitable for use in combination with other accessories		Page	Accessories															
			Magnus				Otesus				Alphamaquet				Meera			
			1180.16A1/A2/A3/A5/A6				1160.10A0/D0				1150.20B0/D0*				7200.01B0/B2			
			1160.16A0/A1/A2				1160.30A0/D0				1150.30B0/D0*				7000.01B0			
Accessories	Description	Page	e		e	e		e	e		e	e		e	e		e	e
1004.42B0	Joint-supporting arm	96	e		e	e		e	e		e	e		e	e		e	e
1004.75A0	Wrist support	123	e		e	e		e	e		e	e		e	e		e	e
1004.8000	Humerus counter traction post	103	f		f	f		f	f		f	f		f	f		f	f
1004.84A0	Knee positioning device	108	f		f	f		f	f		f	f		f	f		f	f
1004.84B0	Knee positioning device	109	f		f	f		f	f		f	f		f	f		f	f
1004.85B0	Countertraction post for femur	93	e		e	e		d	e		e	e		e	e		e	e
1004.86B0	Universal support	98	e		e	e		e	e		e	e		e	e		e	e
1004.8700	Leg plate clamp	94	e		e	e		e	e		e	e		e	e		e	e
1004.8800	Side rail clamp	99	e		e	e		e	e		e	e		e	e		e	e
1004.91A0	Side rail extension	28	f		f	f		f	f		f	f		f	f		f	f
1004.91B0	Side rail extension	28	f		f	f		f	f		f	f		f	f		f	f
1004.93B0	Condyle fixation	95	e		e	e		e	e		e	e		e	e		e	e
1004.9800	Humerus positioning device	103	f		f	f		f	f		f	f		f	f		f	f
1005.01B0	Knee crutches	76	e		e	e		e	e		e	e		e	e			e
1005.22A0	Extension rail	29						g										
1005.2300	Cross-bar attachment	114	g		g	g		g	g		g	g		g	g		g	g
1005.30A0	Fisso Arm support lateral	51	f		f	f		f	f		f	f		f	f		f	f
1005.40A0	Locking device for side rails	28	f		f	f	f	f	f		f	f		f	f		f	f
1005.48B0	Skull clamp holder (not 1160.16A1)	114	e	d	e	e	d		e		e	e		e	e		e	
1005.49B0	DORO® Skull clamp radiolucent	114	e	e	e	e	e	e	e		e	e		e	e		e	
1005.42A0	DORO® LUCENT Skull Clamp	115	e	e	e	e	e	e	e		e	e		e	e		e	
1005.50A0	DORO® Adjustable Base Unit	112	e		e	e		e	e		e	e		e	e		e	
1005.51A0	DORO® Swivel Adaptor	113	e		e	e		e	e		e	e		e	e		e	
1005.52C0	DORO® QR3 Skull Clamp	112	e		e	e		e	e		e	e		e	e		e	
1005.53A0	DORO® Reusable Skull Pins, Pediatric	113	e	e	e	e	e	e	e		e	e		e	e		e	
1005.53B0	DORO® Reusable Skull Pins, Adult	113	e	e	e	e	e	e	e		e	e		e	e		e	
1005.56A0	DORO® Swivel Horseshoe Head Rest, adult	113	e		e	e	e		e		e	e		e	e		e	
1005.63A0	DORO® Transitional Member Radiolucent, long	115	e		e	e	e	e	e		e	e		e	e		e	
1005.63B0	DORO® Transitional Member Radiolucent, short	115	e		e	e	e	e	e		e	e		e	e		e	
1005.64A0	DORO® Swivel Adaptor Radiolucent	116	e		e	e		e	e		e	e		e	e		e	
1005.67A0	DORO® Sterile Disposable Skull Pins, adult	113	e	e	e	e	e	e	e		e	e		e	e		e	
1005.71A0	DORO® Table Adaptor Maquet, radiolucent	115	e		e	e		d	e		e	e		e	e		e	

*Previously sold product

		Page																
			Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0*	Alphamaquet 1150.30B0/D0*	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0
a	Positioning aid																	
b	Accessory trolley																	
c	Medical furniture																	
d	Can be fitted to OR table mounting point																	
e	May be only used on the OR table mounting point in conjunction with other accessories																	
f	Accessories for side rails																	
g	Side rail accessories, only suitable for use in combination with other accessories																	
1005.72A0	DORO® Interface Skull Clamp Radiolucent	116	e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	e
1005.78C0	Pediatric operating table top	67										d	d	d	d	d	d	d
1005.88A0	Safety side rails, pair	20	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1005.89A0	Cover for safety side rails	20	g	g	g	g	g	g	g	g	g	g	g	g	g	g	g	g
1005.95A0	Carbon-fiber arm surgery table	52	f		f	f			f	f	f	f	f	f	f	f	f	f
1007.04BC	Kneeling frame	122			d	d						d	d	d	d	d	d	
1007.04D0	Kneeling frame	122	d															
1007.05A0	Sitting bracket with pelvis support	123	g		g	g				g	g	g	g	g	g	g	g	g
1007.06A0	Accessory trolley	32	b															
1007.07A0	Trolley	123			b	b			b	b	b	b	b	b	b	b	b	b
1007.07D0	Trolley	122	b															
1007.12A0	Guiding roller for head-side traction	116	e		e	e			e	e	e	e	e	e	e	e	e	e
1007.16A0	Head rest adapter	36	e		e	e			e	e	e	e	e	e	e	e	e	e
1007.18B0	Back plate for shoulder operations	106			d	d					d	d	d	d	d	d	d	d
1007.19BC	Back plate for shoulder operations	107								d								
1007.20A0	Trolley	106			b	b			b	b	b	b	b	b	b	b	b	b
1007.21B0	Back plate for shoulder operations	108			d	d					d	d	d	d				d
1007.24A1	Universal frame for patient positioning	118	e		e	e				e	e	e	e	e	e	e	e	e
1007.25A0	Side rail clamp	120	e		e	e				e	e	e	e	e	e	e	e	e
1007.26A0	Short support plate	119	e		e	e				e	e	e	e	e	e	e	e	e
1007.26B0	Long support plate	119	e		e	e				e	e	e	e	e	e	e	e	e
1007.27A0	Adapter for universal frame for Magnus	119	e															
1007.28B0	Adjustable positioning elements	120	e		e	e				e	e	e	e	e	e	e	e	e
1007.30A0	Adapter for radiolucent spine surgery frame	120	e		e	e				e	e	e	e	e	e	e	e	e
1007.40A0	Traction bar with ball joint	91																d
1007.40B0	Traction bar with ball joint	90			f	f	f				f							
1007.41A0	Traction bar	92																d
1007.41B0	Traction bar	90			f	f	f				f							
1007.42AC	Pair of leg plates	79																d
1007.43A0	Screw tension device	97	e		e	e			e	e	e	e	e	e	e	e	e	e
1007.44A0	Sliding clamp with side rail	92																e
1007.45A0	Support plate	99																e
1007.46A0	Trolley for traction bar	92																b
1007.46B0	Trolley for traction bar	90			f	f	f				f							

*Previously sold product

		Page																
			Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0*	Alphamaquet 1150.30B0/D0*	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0
a	Positioning aid																	
b	Accessory trolley																	
c	Medical furniture																	
d	Can be fitted to OR table mounting point																	
e	May be only used on the OR table mounting point in conjunction with other accessories																	
f	Accessories for side rails																	
g	Side rail accessories, only suitable for use in combination with other accessories																	
1007.47A0	Accessory trolley	92																b
1007.48A0	Frame for sterile drape	99	f		f	f			f	f	f	f	f	f	f	f	f	f
1007.49A0	Extension frame	99	e		e	e			e	e	e	e	e	e	e	e	e	
1007.50AC	Adapter clamp	94																e
1007.51A0	Vertical adjustment, long	98																e
1007.51B0	Vertical adjustment, short	98	e		e	e			e	e	e	e	e	e	e	e	e	e
1007.52A0	Bar extension	97																d
1007.53A0	Kneeling frame	122			d	d					d	d	d	d	d		d	d
1007.54A0	Trolley	123			b	b					b	b	b	b	b		b	b
1008.01A0	KTEK Curve spine frame	120	e	e	e	e	e			e	e	e	e	e	e	e	e	e
1008.02A0	KTEK Curve cart	121	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b
1008.02B0	KTEK Stirrup Cart	73	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b
1008.03A0	Curve comfort covers	121																
1008.03B0	Curve standard covers	121																
1008.03C0	Head block	121																
1008.05A0	Secure release clamp	73	f		f	f			f	f	f	f	f	f	f	f	f	f
1008.06B0	KTEK Nissen Straps, standard	83	f		f	f			f	f	f	f	f	f	f	f	f	f
1008.06D0	KTEK Nissen Straps, bariatric	84	f		f	f			f	f	f	f	f	f	f	f	f	f
1008.07E0	Clean catcher non sterile	22	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1008.07F0	Clean catcher sterile	22	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1008.08A0	KTEK 250 Stirrup – Junior	72	g		g	g			g	g	g	g	g	g	g	g	g	g
1008.09A0	KTEK Comfort Stirrup 800	72	g		g	g			g	g	g	g	g	g	g	g	g	g
1008.10A0	Knee crutch lift assist	74	g		g	g			g	g	g	g	g	g	g	g	g	g
1008.11A0	Radiation protection	26	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1008.12A0	Patient safety strap, medium	30	f		f	f			f	f	f	f	f	f	f	f	f	f
1008.14C0	Arthroscopic leg holder	109	f		f	f			f	f	f	f	f	f	f	f	f	f
1008.14B0	Inlay for arthroscopic leg holder	109																
1008.15A0	KTEK Trigger Armboard	49	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1008.17A0	KTEK Armboard Swivel, large	47	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1008.17B0	KTEK Armboard Swivel, small	48	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f	f
1009.01C0	Infusion stand	20	f	g	f	f	g	f	f	f	f	f	f	f	f	f	f	f
1009.02C0	Meniscus positioning device	108	f		f	f			f	f	f	f	f	f	f	f	f	f
1009.69A0	Sensor Drive	27									f	f	f					f
1118.50B0	Pair of leg plates	78			d	d						d	d	d	d	d	d	d

*Previously sold product

a b c d e f g Positioning aid Accessory trolley Medical furniture Can be fitted to OR table mounting point May be only used on the OR table mounting point in conjunction with other accessories Accessories for side rails Side rail accessories, only suitable for use in combination with other accessories		Page	Accessories for OR tables																															
			Maquet Magnus 1180.10A0		Maquet Magnus 1180.16A1/A2/A3/A5/A6		Maquet Otesus 1160.10A0/D0		Maquet Otesus 1160.30A0/D0		Maquet Otesus 1160.16A0/A1/A2		Alphamaquet 1150.16B0		Alphamaquet 1150.20B0/D0*		Alphamaquet 1150.30B0/D0*		Corin 7700.01A0/B0/A2/B2		Alphamaxx 1133.22B4/B5		Maquet Meera 7200.01B0/B2		Maquet Meera ST 7100.01B0/B2		Maquet Meera CL 7000.01B0		Maquet Lyra 1118.10B0/K0		Betaclassic 1118.01B0/K0		Maquet Yuno II 1433.02B0	
1118.60B0	Pair of leg plates	77																																
1130.53B0	Head rest SFC	40				e	e						e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	d	e	e				
1130.53D0	Head rest PUR	40				e	e						e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	d	e	e				
1130.54B0	Connection bracket	36				e	e						e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	d	e	e				
1130.55B0	Seat plate extension	66																											d					
1130.64C0	Head rest	40				e	e						e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	d	e	e				
1130.64G0	Head rest	40				e	e						e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	d	e	e				
1130.65A0	Transfer board	74				d	d								d	d	d	d	d	d	d	d	d	d	d	d	d	d	d	d	d			
1130.67B0	Head rest	41				e	e						e	e	e	e	e	e	e	e	e	e	e	e	e	e	e	d	e	e				
1130.70AC	Maquetmatic leg plate	76				d	d										d	d	d	d	d	d	d	d	d	d	d			d				
1130.71B0	Filler piece	76																d	d	d	d	d	d	d	d	d	d			d				
1130.81A0	Head rest adapter	39				d	d										d	d	d	d	d	d	d	d	d	d	d			d				
1131.31BC	Extension plate	63				d	d										d	d	d	d	d	d	d	d	d	d				d				
1131.55BC	Seat plate extension	66				d	d										d	d	d	d	d	d	d	d	d	d				d				
1132.34A0	Back plate for shoulder operations	107				d	d										d	d	d	d	d	d	d	d	d	d				d				
1132.40A0	Accessory trolley	31	b	b		b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b	b			
1132.45AC	Back plate	59				d	d										d	d	d	d	d	d	d	d	d	d				d				
1132.65A0	Transfer board	75				d	d										d	d	d	d	d	d	d	d	d	d	d	d	d	d				
1132.81A0	Head rest adapter	39																										d	d					
1133.53BC	Pair of leg plates	81				d	d										d	d	d	d	d	d	d	d	d	d	d	d	d	d				
1133.58BC	Leg plate	80				d	d										d	d	d	d	d	d	d	d	d	d				d				
1133.67BC	Pair of carbon-fiber leg plates	78				d	d										d	d	d	d	d	d	d	d	d	d	d	d	d	d				
1133.73AC	Pair of leg plates	83				d	d										d	d	d	d	d	d	d	d	d	d	d	d	d	d				
1133.73BC	Pair of leg plates	82				d	d										d	d	d	d	d	d	d	d	d	d	d	d	d	d				
1150.31B0	Back plate	62														d																		
1150.31D0	Back plate	62												d	d																			
1150.45BC	Carbon-fiber multi-purpose plate	59													d																			
1150.55HC	Seat plate extension	65											d	d																				
1150.59AC	Extension device	88												d																				
1150.64B0	Pair of leg plates	80	e			e	e						d	e	e	e	e	e	e	e	e	e	e	e	e	e	e	e						
1150.64E0	Pair of leg plates	80																												d				
1150.70A0	Accessory adapter	60		e						d							e	e	e	e	e	e	e	e	e	e								

*Previously sold product

[illegible]

*Previously sold product

a b c d e f g Positioning aid Accessory trolley Medical furniture Can be fitted to OR table mounting point May be only used on the OR table mounting point in conjunction with other accessories Accessories for side rails Side rail accessories, only suitable for use in combination with other accessories <th rowspan="7">Page</th> <th colspan="18"></th>		Page																		
			Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0*	Alphamaquet 1150.30B0/D0*	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0		
1180.11A1	Padded plate	64	e																	
1180.11B0	Motorized joint module	63	d																	
1180.11B1	Padded plate	64	e																	
1180.11C1	Urology seat plate	64	d																	
1180.19A0	Extension device	87	e																	
1180.31A0	Back plate	60	e																	
1180.32A0	Extension plate	63	d																	
1180.34A0	Shoulder module	106	e																	
1180.34T0	Trolley	106	b																	
1180.36A0	Connection bracket	36	d																	
1180.37A0	Carbon-fiber back plate	58	e																	
1180.38A0	Universal extension plate	62	d																	
1180.41A0	Pad for rectal positioning	82	e																	
1180.44AC	Carbon-fiber multi-purpose plate for lithotomy position	58	d																	
1180.45A0	Carbon-fiber multi-purpose plate	58	d																	
1180.46A0	Extension plate	43	e		e	e											e			
1180.47A0	Catheter holder	27		f			f													
1180.48B0	Lightweight leg plate	77	e																	
1180.49A0	Side rail adapter pair	27		f			f													
1180.50A0	Head rest	42	d																	
1180.51A0	Adapter clamp with transfer board	75	e																	
1180.53A0	Head rest	43	d																	
1180.54A0	Pair of leg plates	81	e																	
1180.54B0	Pair of leg plates	81	d																	
1180.55A0	Universal module	64	e																	
1180.56A0	Transfer board pad	75	d																	
1180.57A0	Transfer board	75	d																	
1180.72X0	Accessory adapter	60	e	d	e	e	d				e						e			
1419.01B0	Traction device	88										e	e	e						
1419.01C0	Traction device	88													e	e				
1419.01X0	Traction device	88																		
1419.01HC	Traction device	88									e									

*Previously sold product

		Page																		
a	Positioning aid		Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0*	Alphamaquet 1150.30B0/D0*	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0		
b	Accessory trolley																			
c	Medical furniture																			
d	Can be fitted to OR table mounting point																			
e	May be only used on the OR table mounting point in conjunction with other accessories																			
f	Accessories for side rails																			
g	Side rail accessories, only suitable for use in combination with other accessories																			
1419.01JC	Traction device	88											e	e	e					
1433.33AC	Carbon-fiber back plate	59			d	d					d	d	d	d	d			d		
1433.34AC	Carbon-fiber back plate for shoulder operations	107			d	d					d	d	d	d	d			d		
1433.41AC	Adapter plate for traction	89																d		
1433.42A0	Femur hook	93																f		
1433.45AC	Carbon-fiber plate for insertion	57																d		
1433.45BC	Carbon-fiber insertable plate	57																d		
1433.66AC	Carbon-fiber extension plate with 3 holes	91																d		
1433.66BC	Carbon-fiber extension plate with 1 hole	91																d		
1433.67A0	Tibia adapter	94																e		
4000.05A0	Pad	142	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4001.59A0	Bowl	144	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4001.60A0	Clamp ring	144	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4006.01A0	Oasis Elite Closed Head Ring Adolescent	126	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.02A0	Oasis Elite Closed Head Ring Adult	126	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.03A0	Oasis Elite Open Head Ring Adolescent	126	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.04A0	Oasis Elite Open Head Ring Adult	126	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.05A0	Universal operating table pad	129	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.06A0	Oasis ¾ Length Table Pad	129	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.07A0	Oasis Full Length Table Pad	129	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.08A0	Oasis Torso Section	129	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.09A0	Oasis Perineal Table Pad	130	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.10A0	Oasis Universal Patient Positioner	134	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.11A0	Oasis Universal Patient Positioner	134	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.12A0	Oasis Armboard Pad	131	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.13A0	Oasis Armboard Pad	131	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.14A0	Oasis Crutch Stirrup Pads	132	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.15A0	Oasis Crutch Stirrup Pads	132	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.16A0	Oasis Elite Heel Pad	133	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.17A0	Oasis Thigh Pad	132	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.18A0	Oasis Elite Lateral Positioner	135	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.19A0	Oasis Elite Prone Head Rest, large	127	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		

*Previously sold product

a b c d e f g		Page																		
			Maquet Magnus 1180.10A0	Maquet Magnus 1180.16A1/A2/A3/A5/A6	Maquet Otesus 1160.10A0/D0	Maquet Otesus 1160.30A0/D0	Maquet Otesus 1160.16A0/A1/A2	Alphamaquet 1150.16B0	Alphamaquet 1150.20B0/D0*	Alphamaquet 1150.30B0/D0*	Corin 7700.01A0/B0/A2/B2	Alphamaxx 1133.22B4/B5	Maquet Meera 7200.01B0/B2	Maquet Meera ST 7100.01B0/B2	Maquet Meera CL 7000.01B0	Maquet Lyra 1118.10B0/K0	Betaclassic 1118.01B0/K0	Maquet Yuno II 1433.02B0		
4006.21A0	Oasis Elite Supine Head Rest Adult	127	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.22A0	Oasis Elite Supine Head Rest Adolescent	127	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.23A0	Azure Head Protector	128	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.24A0	Azure Sacral Protector, small	130	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.25A0	Oasis Elite Flat Bottom Chest Roll	134	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.26A0	Azure Arm Protector, short	131	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.27A0	Azure Arm Protector, long	131	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.28A0	Azure Leg Protector	132	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.29A0	Azure Heel Protectors	133	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.3100	Prone positioning pad	136	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4006.3200	Lateral positioning pad	135	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a	a		
4009.3000	Infusion stand	21	f	g	f	f	f	g	f	f	f	f	f	f	f	f	f	f		
4402.9002	Hand and arm positioning table	142	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4403.9302	Instrument stand	140	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4404.32A3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4404.32B3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4404.32C3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4404.32X3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4405.32A3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4405.32B3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4405.32C3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4405.32X3	Instrument table	141	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4406.32A3	Instrument table	142	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4406.32B3	Instrument table	142	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4406.32C3	Instrument table	142	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4406.32X3	Instrument table	142	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4464.9500	Instrument stand	140	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4464.9600	Instrument stand	140	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4464.98A0	Mayo instruments table	140	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4464.9X00	Instrument stand	140	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4602.3402	OR swivel stool	146	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		
4632.01A0	OR swivel stool pneumatic	145	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c		

*Previously sold product

a b c d e f g		Page																	
4632.01E0	OR swivel stool pneumatic	145	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4634.01A0	OR swivel stool hydraulic	146	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4634.01C0	OR swivel stool hydraulic	146	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4805.91A3	IV pole MRI	144	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4805.91B3	Infusion stand	144	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4809.9003	Bowl stand	140	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4901.3503	Disposal trolley	143	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4902.3503	Disposal trolley	143	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4904.95A0	Heated bowl holder	145	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4904.96A0	Heated bowl holder, swiss version	145	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4906.9003	OR step stool	143	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4906.9103	OR step stool	143	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
4907.9003	OR step stool	143	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	c	
6005.25A0	Accessory adapter	119	e		e	e												e	
SU01.04990001	Swivel arm board	49	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f	
SU01.04990002	Swivel arm board	49	f	g	f	f	g	g	f	f	f	f	f	f	f	f	f	f	

*Previously sold product

Index

A

Accessory adapter
1150.70A0 60
1150.72A0 60
1150.72B0 59
1180.72X0 60
6005.25A0 119

Accessory stand
for carbon-fiber back plate
1003.87A0 59

Accessory trolley
1003.48A0 31
1007.06A0 32
1007.47A0 92
1132.40A0 31
1160.57A0 91

Accessory trolley for
carbon-fiber insertable plates
1003.86A0 58

Adapter clamp
1007.50AC 94

Adapter clamp with transfer
board
1180.51A0 75

Adapter for CF table top
1001.64X0 52

Adapter for CF table top
(not for 1160.16A2)
1001.64C0 52

Adapter for radiolucent
spine surgery frame
1007.30A0 120

Adapter for universal frame for
Magnus
1007.27A0 119

Adapter plate for extension
1433.41AC 89

Adapter plate for traction
1433.41XC 76

Adapter Trimano Fortis
1002.31A0 104

DORO® Adjustable Base Unit
1005.50A0 112

Adjustable head rest adapter
1002.10A0 38

Adjustable positioning elements
1007.28B0 120

Anesthesia frame
1002.54A0 18

Anesthesia screen
1002.57A0 19

Anesthesia screen extension
1002.59A0 19

Arm positioning device PUR
1001.44D0 47

Arm protector
1001.40A0 53
1002.25A0 54

Arm protector, curved
1001.40B0 53

Arm rest
1002.49A0 50

Arm rest for beach chair position
1002.43B0 105

Arm rest for lateral position
1002.43A0 105

Arm restraint cuff
1001.4600 51
1002.24C0 54

Arm support
1001.42A0 53
1001.44E0 47
1001.44F0 47
1001.6000 48
1002.32A0 50

Armpit support
1001.92C0 102

Arthroscopic leg holder
1008.14C0 109

Azure Arm Protector, long
4006.27A0 131

Azure Arm Protector, short
4006.26A0 131

Azure Head Protector
4006.23A0 128

Azure Heel Protectors
4006.29A0 133

Azure Leg Protector
4006.28A0 132

Azure Sacral Protector, small
4006.24A0 130

B

Back plate
1132.45AC 59
1150.31B0 62
1150.31D0 62
1160.35AC 62
1160.35BC 62
1180.31A0 60

Back plate for shoulder
operations
1007.18B0 106
1007.19BC 107
1007.21B0 108
1132.34A0 107

Back/ buttocks support
1002.11A0 23

Bar extension
1003.54A0 96
1007.52A0 97

Bariatric set, EU
1001.74B0 66

Body strap
1001.56A0 29
1001.56B0 29
1001.58A0 30
1001.59A0 29
1001.59B0 30

Body support
1002.97A0 24

Bowl
4001.59A0 144

Bowl stand
4809.9003 140

Buttocks support
1002.11E0 24

C

Carbon-fiber arm surgery table
1005.95A0 52

Carbon-fiber back plate
1433.33AC 59
1180.37A0 58

Carbon-fiber pair of leg plates
1133.67BC 78

Carbon-fiber back plate
for shoulder operations
1433.34AC 107

Carbon-fiber extension plate with
3 holes
1433.66AC 91

Carbon-fiber extension plate with
1 hole
1433.66BC 91

Carbon-fiber insertable plate
1160.45AC 57
1160.45BC 57
1160.45XC 57
1433.45BC 57

Carbon-fiber multi-purpose plate
1180.45A0 58
1150.45BC 59

Carbon-fiber multi-purpose plate
for lithotomy position
1180.44AC 58

Carbon-fiber plate for insertion
1433.45AC 57

Catheter holder
1180.47A0 27

Clamp ring
4001.60A0 144

Clamp with blade interface
1003.24A0 71

Clean catcher non sterile
1008.07E0 22

Clean catcher sterile
1008.07F0 22

Clip-on device for
minimally invasive surgery
1003.60D0 87

Column protection device
1004.07A0 31

Condyle fixation
1004.93B0 95

Connection bracket
1130.54B0 36
1180.36A0 36
1004.05A0 37

Countertraction post for tibia and
fibula
1003.50C0 95

Countertraction post for femur
1004.85B0 93

Cover for safety side rails
1005.89A0 20

Cross-bar attachment
1005.2300 114

Curve comfort covers
1008.03A0 121

Curve standard covers
1008.03B0 121

D

Disposable arm rest
1002.52A0 105

Disposable inlay and
mounting brackets, big size
1003.76C0 100

Disposable inlay and
mounting brackets, small size
1003.76A0 100

Disposable mounting brackets
1003.78A0 101

Disposal trolley
4901.3503 143
4902.3503 143

DORO® Interface Skull Clamp,
radiolucent
1005.64A0 116

DORO® Reusable Skull Pins,
Adult
1005.53B0 113

DORO® Reusable Skull Pins,
Pediatric
1005.53A0 113

DORO® Sterile Disposable Skull
Pins, adult
1005.67A0 113

DORO® Table Adaptor Maquet,
radiolucent
1005.71A0 115

DORO® Transitional Member
Radiolucent
1005.63A0 115
1005.63B0 115

Double lateral support
1002.17A0 25

E

Trimano elbow countertraction
post
1002.45A0 105

Elbow support
1004.23C0 103

Extension device
1150.59AC 88
1180.19A0 87

Extension frame
1007.49A0 99

Extension plate
1131.31BC 63
1160.32AC 61
1160.32BC 61
1160.32DC 61
1160.33BC 61
1160.33DC 61
1160.33XC 61
1180.46A0 43
1180.32A0 63

Extension rail
1005.22A0 29

F

Femur hook
1433.42A0 93

Filler piece
1130.71B0 76

Fisso Arm support lateral
1005.30A0 51

Trimano finger fixation
1002.48A0 104

Fixture
1002.65A0 38

Fixture for body supports
1002.19C0 23

Flat head rest
1002.73A0 38

Foot plate
1001.86C0 83
1001.97A0 100

Foot plate, adult
1001.8700 101

Foot plate, child
1001.9000 102

Foot plate mount
1003.49B0 96

Frame for sterile drape
1007.48A0 99

G

Guiding bar
1003.5100 95

Guiding roller for head-side
traction
1007.12A0 116

H

Hand and arm positioning table
4402.9002 142

Head block
1008.03C0 121

Head rest
1002.81D0 39
1130.64G0 40
1130.53B0 40
1130.67B0 41
1150.71B0 44
1150.71C0 44
1160.53AC 42
1160.53BC 42
1180.50A0 42
1180.53A0 43
1160.64A0 41
1160.64BC 40

Head rest that can be angled
1002.83A0 117

Head rest adapter
1002.07A0 36
1007.16A0 36
1130.81A0 39
1132.81A0 39

Head ring
1000.5700 128

Head support for shoulder
operations
1002.15A0 43

Heated bowl holder
4904.95A0 145
4904.96A0 145

Heavy duty leg strap
1001.53A0 83

Horseshoe head rest
1002.03A0 116
1002.04B0 118

Horseshoe head rest – two
section
1002.72A0 38

Horseshoe head rest
with one-handed adjustment
1002.66A0 39

Humerus counter traction post
1004.8000 103

Humerus positioning device
1004.9800 103

I

Infusion stand
1009.01C0 20
4009.3000 21
4464.9X00 140
4805.91B3 144

Inlay for arthroscopic leg holder
1008.14B0 109

Instrument stand
4403.9302 140
4464.9500 140
4464.9600 140

Instrument table
1002.98A0 21
4404.32A3 141
4404.32B3 141
4404.32C3 141
4405.32A3 141
4405.32B3 141
4405.32C3 141
4405.32X3 141
4406.32A3 142
4406.32B3 142
4406.32C3 142
4406.32X3 142

Intermediate adaptor DORO®
1005.72A0 116

IV pole MRI
4805.91A3 144

J

Joint-supporting arm
1004.42B0 96

K

Knee crutches, standard
1003.81A0 71

Knee crutches, adipositas
1003.82A0 71

Knee crutches
1005.01B0 76

Knee crutch lift assist
1008.10A0 74

Knee positioning device
1004.84A0 108
1004.84B0 109

Kneeling frame
1007.04BC 122
1007.04D0 122
1007.53A0 122

KTEK 250 Stirrup – Junior
1008.08A0 72

KTEK Comfort Stirrup 800
1008.09A0 72

KTEK Curve Cart
1008.02A0 121

KTEK Curve Spine Frame
1008.01A0 120

KTEK Nissen Straps, standard
1008.06B0 83

KTEK Nissen Straps, bariatric
1008.06D0 84

KTEK Stirrup Cart
1008.02B0 73

KTEK Armboard Swivel, large
1008.17A0 47

KTEK Armboard Swivel, large
1008.17B0 48

KTEK Trigger Armboard
1008.15A0 49

L

Lateral positioning pad
4006.3200 135

Lateral protection
1004.15A0 19

Lateral support
1002.11C0 23
1002.21C0 24
1002.36A0 25

Leg holder
1001.65A0 73
1001.73A0 74
1002.2300 73
1003.80A0 71

Leg plate
1133.58BC 80

Leg plate clamp
1004.8700 94

Leg rest
1001.94A0 136

Leg restraint cuff
1001.4700 84

Leg strap
1001.57A0 84

Leg support
1001.66A0 83

Lightweight leg plate 1180.48B0.....	77	Oasis Elite Supine Head Rest Adult 4006.21A0.....	127
Locking device for side rails 1005.40A0.....	28	Oasis Elite Closed Head Ring Adolescent 4006.01A0.....	126
Long support plate 1007.26B0.....	119	Oasis Elite Closed Head Ring Adult 4006.02A0.....	126
Maquetmatic leg plate 1130.70AC.....	76	Oasis Full Length Table Pad 4006.07A0.....	129
Mayo instrument stand 1003.32C0.....	21	Oasis Perineal Table Pad 4006.09A0.....	130
Mayo instruments table 4464.98A0.....	140	Oasis Thigh Pad 4006.17A0.....	132
Meniscus positioning device 1009.02C0.....	108	Oasis Torso Section 4006.08A0.....	129
Motorized joint module 1180.11A0 1180.11B0	63 63	Oasis Universal Patient Positioner 4006.10A0..... 4006.11A0.....	134 134
Multi-purpose plate for arm surgery 1001.63A0.....	51	OR step stool 4906.9003 4906.9103 4907.9003	143 143 143
Oasis ¾ Length Table Pad 4006.06A0.....	129	OR swivel stool 4602.3402.....	146
Oasis Armboard Pad 4006.12A0..... 4006.13A0.....	131 131	OR swivel stool pneumatic 4632.01E0.....	143
Oasis Crutch Stirrup Pads 4006.14A0..... 4006.15A0.....	132 132	OR swivel stool hydraulic 4634.01A0..... 4634.01C0.....	146 146
Oasis Elite Flat Bottom Chest Roll 4006.25A0.....	134	P Pad 1000.5600 1000.68C0 1003.7400 4000.05A0.....	136 136 134 142
Oasis Elite Open Head Ring Adolescent 4006.03A0.....	126	Pad for rectal positioning 1180.41A0.....	82
Oasis Elite Lateral Positioner 4006.18A0.....	135	Padded roll for countertraction post 1000.24A0.....	137
Oasis Elite Open Head Ring Adult 4006.04A0.....	126	Padded plate 1180.11A1 1180.11B1	64 64
Oasis Elite Heel Pad 4006.16A0.....	133	Patient safety strap, medium 1008.12A0.....	30
Oasis Elite Prone Head Rest, large 4006.19A0.....	127	Pediatric operating table top 1005.78C0.....	67
Oasis Elite Prone Head Rest, medium 4006.20A0.....	127	Pair of leg plates 1007.42AC.....	79
Oasis Elite Supine Head Rest Adolescent 4006.22A0.....	127		

M

O

1118.60B0 1118.50B0 1133.53BC..... 1133.73AC..... 1133.73BC..... 1150.64B0..... 1150.64E0 1150.73BC..... 1150.73DC..... 1160.50AC..... 1160.50BC..... 1160.50DC..... 1180.54A0..... 1180.54B0.....	77 78 81 83 82 80 80 77 77 79 79 79 81 81	Pelvis plaster cast device 1003.3600.....	100
Pin-joint arm for body supports 1002.40A0.....	22	Plexus cushion 1000.6900.....	135
Prone positioning pad 4006.3100.....	136	Pubis/ sternum support 1002.11B0.....	23
R Radial setting clamp 1003.23C0..... 1003.61A0.....	18 18	Radiation protection 1008.11A0	26
Retrofit kit for standard holder for DORO® Skull Clamps 1160.36A0.....	61	Retrofit kit for standard holder for Sugita Skull Clamps 1160.36B0.....	62
Reusable inlay 1003.77A0	101	Rigid head rest 1002.82A0.....	117
Rotating and tilting clamp 1003.34A0.....	96	S Safety side rails, pair 1005.88A0.....	20
Screw tension device 1003.3700 1007.43A0	97 97	Seat plate extension 1130.55B0 1131.55BC 1150.55HC 1160.55AC	66 66 65 65
1160.55BC 1160.55DC.....	65 65	Seat plate extension, carbon fiber 1003.57B0.....	94
Secure release clamp EU 1008.05A0.....	73	Sensor Drive 1009.69A0.....	27
Short support plate 1007.26A0.....	119	Shoulder module 1180.34A0.....	106
Shoulder support 1002.09C0.....	25	Shoulder adapter pair 1180.49A0.....	27
Side rail clamp 1004.8800..... 1007.25A0.....	99 120	Side rail extension 1001.54A0 1004.91A0 1004.91B0	28 28 28
Sitting bracket with pelvis support 1007.05A0.....	123	Skull clamp holder 1005.48B0.....	114
DORO® QR3 Skull Clamp 1005.52C0.....	112	DORO® Skull clamp radiolucent 1005.49B0.....	114
Skull clamp holder 1005.48B0.....	114	DORO® LUCENT Skull clamp 1005.42A0.....	115
Sliding clamp 1003.58B0.....	93	Sliding clamp with side rail 1007.44A0.....	92
Soiled instrument bowl 1003.44A0.....	21	Standard clamp 1003.22C0.....	18
Support bar 1002.18A0.....	30	Support plate 1007.45A0.....	99
Support plate 1160.59BC.....	89	Supporting roll 1002.11D0.....	24

R

S

DORO® Swivel Adaptor 1005.51A0.....	113	Trolley 1003.46A0..... 1007.54A0 1007.07A0..... 1007.07D0..... 1007.20A0..... 1180.34T0.....	22 123 123 122 106 106
Swivel arm board SU01.04990001 SU01.04990002.....	49 49	Trolley for traction bar 1007.46A0 1007.46B0	92 90
DORO® Swivel Horseshoe Head Rest, adult 1005.56A0.....	113	Trimano Adapter 1002.42A0.....	104
OR swivel stool pneumatic 4632.01A0.....	145	Trimano Fortis 1002.30A0.....	50, 104
T Table with extension, long 1001.75A0.....	25	Tube holder 1002.55A0.....	19
Table with extension, short 1001.76A0.....	26	Tunnel cushion 1000.77A0.....	135
Telescopic bar 1003.59V0.....	98	TUR rinsing set 1003.45D0.....	22
Tibia adapter 1433.67A0	94	U Universal extension plate 1180.38A0.....	62
Traction adapter 1003.72AC	90	Universal frame for patient positioning 1007.24A1	118
Traction bar with ball joint 1007.40B0 1007.40A0	90 91	Universal module 1180.55A0.....	64
Traction bar 1007.41B0 1007.41A0	90 92	Universal operating table pad 4006.05A0.....	129
Traction bar, carbon-fiber version 1003.55A0.....	93	Universal support 1004.86B0.....	98
Traction boot 1003.75A0.....	100	Upper arm posturing plate 1001.45C0.....	104
Traction boot, adult 1001.8800.....	102	Urology seat plate 1180.11C1.....	64
Traction boot, child 1001.9100.....	102	V Vacuum mattress 1000.78A0	137
Traction boot, small 1003.75C0.....	101	Vertical adjustment, long 1007.51A0.....	98
Traction device 1419.01B0..... 1419.01C0 1419.01X0..... 1419.01HC..... 1419.01JC.....	88 88 88 88 88	Vertical adjustment, short 1007.51B0	98
Traction stirrup clamp with rotation 1003.35A0.....	96	W Weinberger arm restraint cuff 1001.48A0	103
Transfer board 1130.65A0 1132.65A0 1180.57A0	74 75 75	Wire basket 0802.0783	32
Transfer board pad 1180.56A0.....	75		

T

V

W

Wrist support 1004.75A0	123
X X-ray cassette holder 1002.94A0.....	26
X-ray cassette slide 1002.93A0	26

U

V

W

This information is intended for an international audience outside the US.

This information is aimed exclusively at healthcare professionals or other professional audiences and are for informational purposes only, is not exhaustive and therefore should not be relied upon as a replacement of the Instructions for Use, service manual or medical advice. Getinge shall bear no responsibility or liability for any action or omission of any party based upon this material, and reliance is solely at the user's risk.

Any therapy, solution or product mentioned might not be available or allowed in your country. Information may not be copied or used, in whole or in part, without written permission by Getinge.

Manufacturer · Maquet GmbH · Kehler Str. 31 · 76437 Rastatt · Germany · info@getinge.com

© 2026 Getinge | Getinge and **GETINGE**  are trademarks or registered trademarks of Getinge AB, its subsidiaries or affiliates.
DMS-0002131 v2 | All rights reserved.

GETINGE 