


Module Spotlight

Loan Management





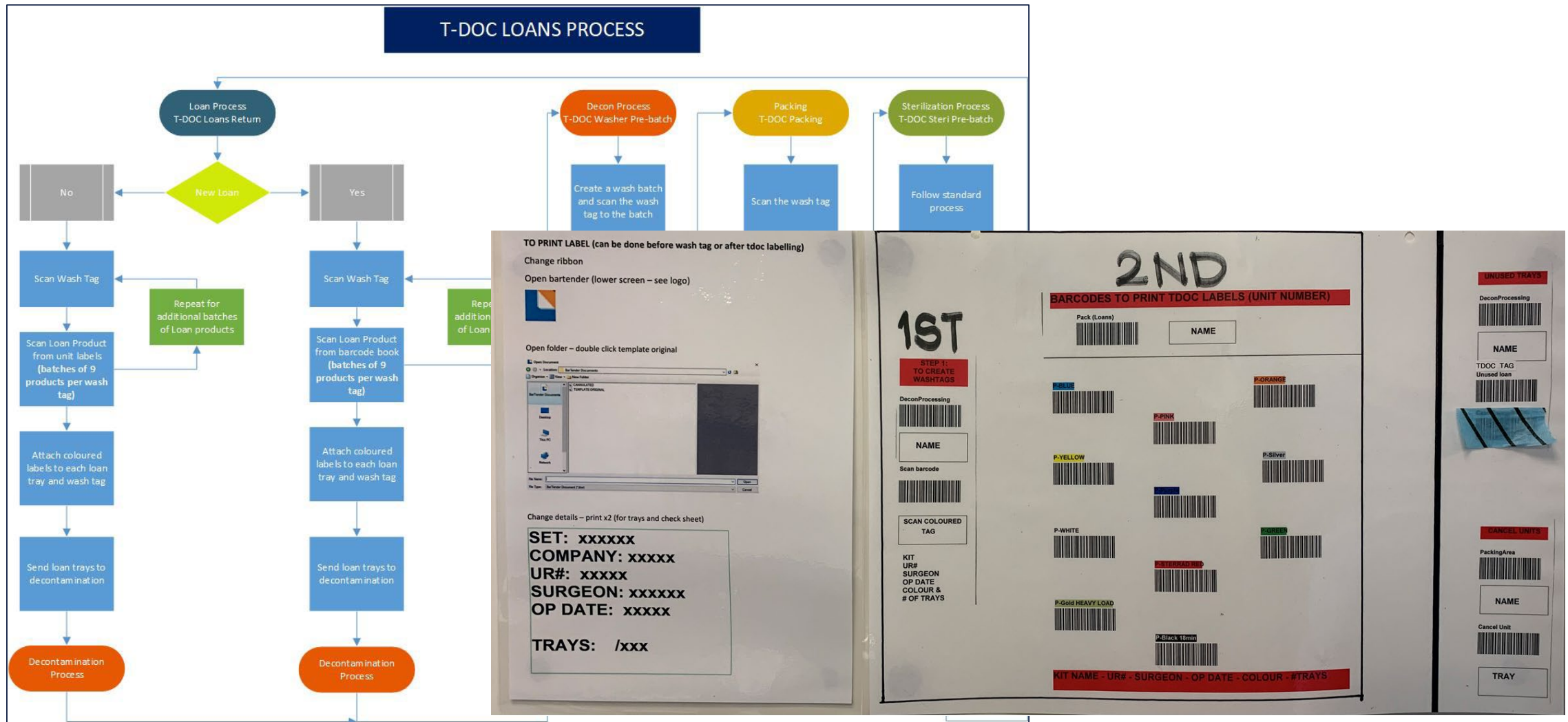









Various ways to currently manage loans



Current Loan Management

Level 1 - Generic Loan Products

PR00025 - External Hospital Loan



Product

PR00025

Type

Tray

Item

EXTLOAN

External Hospital Loan

Traceability

Full traceability

Trace type

Product number

Pack Units > Return/Wash > Pack Units > Sterilize > Store > Patient > Return/Wash > Send to loan company or re-use








Level 2 - Generic Product Serials

Return to company by recording an event or dispatch to company customer



Level 3 - Specific Product Serials

PR00031-01

Serial no.

PR00031-01

⏪

⏩

⏴

⏵

🗨️


Product

PR00031

STRYKER HIP REVISION SET

Unit

Barcodes



Info

Owner

Triggers

User fields

Weight

Fast tracking

Events

Documents

Remark

Status

Normal

Location

Stock

Placement

Return to company by recording an event or dispatch to company customer

Return/Wash > Pack Units > Sterilize > Store > Patient > Return/Wash > Sterilize > Send to loan company or re-use

New Loans Functionality

External loaner interface

Origin


- Growing global demand to handle loaners better
- Dedicated “loaner systems” in the US drive competition

Value


- Improved handling of loaners with integration to the largest loaner systems in the US
- Complete tracking of loaner products in T-DOC with much less work
- Preparation for: A simple stand-alone loaner solution in T-DOC



Workflow

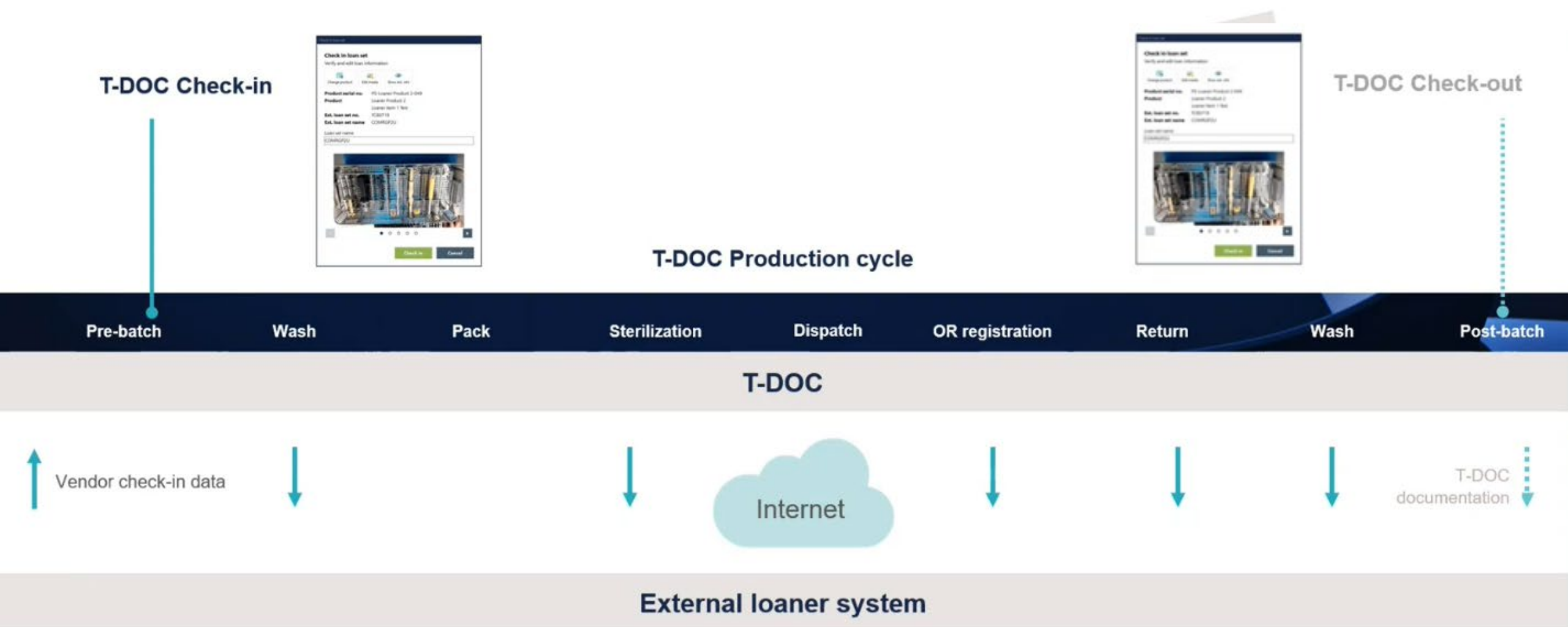

Loan Workflow

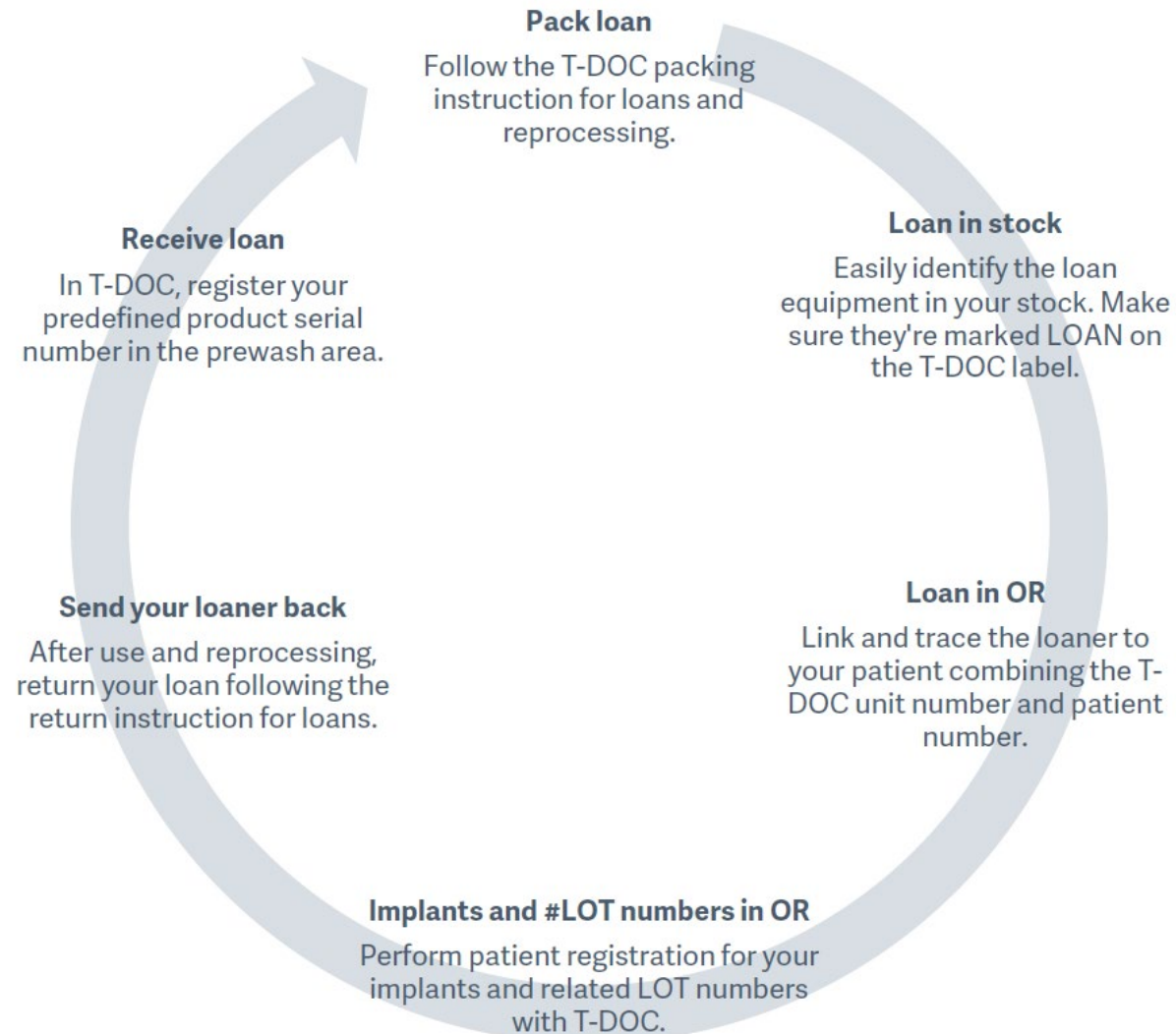



Loaner check-in

Can be performed at:

- Return
- Washer Pre-batch
- Used with Product Serials

Check in loan set


Verify and edit loan information

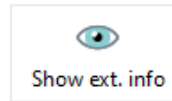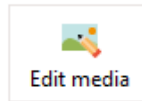

Product serial no.	PR00016-01
Product	PR00016 STRYKER - LOAN TRAYS
Ext. loan set no.	PR00016-01
Ext. loan set name	PR00016-01

Loan set name

STRYKER HIP REVISION SET 01 - DR DONNELLY, PATIENT 569

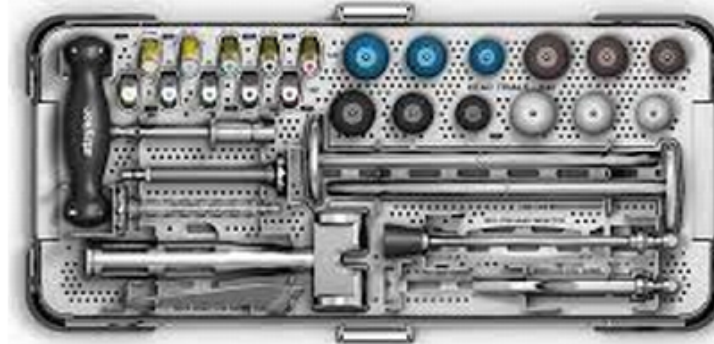

Check in

Cancel


Loaner check-out

Can be performed at:

- Return
- Wash Post Batch
- Washer Pre-batch
- Pack
- Prepare
- Can set to automatically check out after return or wash

Check out loan set


Verify loan information

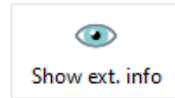

Unit	993
Product serial no.	PR00016-01
Product	PR00016
	STRYKER - LOAN TRAYS
Ext. loan set no.	PR00016-01
Ext. loan set name	PR00016-01

Loan set name

STRYKER HIP REVISION SET 01 - DR DONNELLY, PATIENT 569

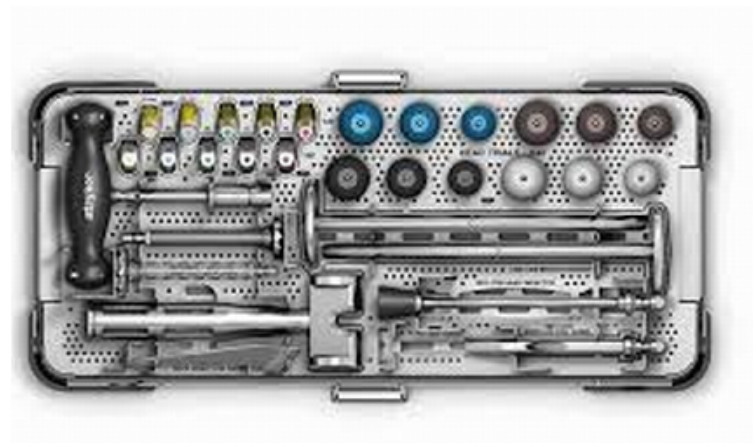

Check out

Cancel

Loaner Management Profiles

Loaner management profiles

Quick search

Full search >>

Recent

Factory

LMPROID

LM001

LMP...

Gen...

LM001 - Generic Loan System

LMPROID

LM001

LMPRONAME

Generic Loan System

General

Remarks

LMPROEXTERNALSYSTYPE

Generic/undefined

LMPROEXTERNALSYSLINKID

LMPROTEMPLATEPRODKEYID

LMPRODEFAULTPRODKEYID

PR00016

LMPROCHECKINLOCAKEYID

Loans Receipt

LMPROAUTOCHECKOUTMODE

Off

LMPROFORCERETURNUNITERROR

LMPROSERINAMECHANGEMODE

Always

LMPROSTATUS

Normal

LMPROEXPARAMS

