

LS3030 DPTE® SS Beta Containers - 201902 - We reserve the right to make technical and construction changes.



DPTE® Beta Containers

Safe bi-directional transfer with re-usable stainless steel

Getinge is a global provider of innovative solutions for operating rooms, intensive care units, sterilization departments and for life science companies and institutions. Based on our firsthand experience and close partnerships with clinical experts, healthcare professionals and medtech specialists, we are improving everyday life for people – today and tomorrow.

Contact · Getinge Infection Control AB · P. O. Box 69 · SE-305 05 Getinge · Sweden · Phone: +46 (0) 10 335 00 00
E-mail: info@getinge.com

www.getinge.com

GETINGE 

DPTE® Stainless Steel Containers

A reliable leaktight sterile transfer solution with steam sterilizable option

Getinge's DPTE® Beta containers are compatible with any DPTE® Alpha port of the same diameter. The patented safety interlock prevents manipulation errors.



Inserts are available for Level 2 containers upwards



Integrity testing using Getinge's wireless, pipeless Transfer Leak Tester

Semi-standard containers and inserts

Level	Ø	Length 300 mm	Length 400 mm	Length 500 mm	Length 600 mm	Length 700 mm
1 Standard autoclavable *	190		•			
	270		•			
	350			•		
2 Semi-standard autoclavable and H ₂ O ₂ sterilizable	190	•	•	•	•	•
	270	•	•	•	•	•
	350	•	•	•	•	•
3 Customizable using pre-designed parts	190	•	•	•	•	•
	270	•	•	•	•	•
	350	•	•	•	•	•
4	Customized containers and inserts, dimensions, parts and finishing to customers' requirements - contact our Sales Team					

* Ø 105 and Ø 460 are part of our range, please contact our Sales Team

Level 1 Standard containers


- Standard lengths
- Container: beta door, beta flange
- Container body 1mm thick
- Handles: attached to flange
- Autoclavable (perforated bottom, filter connection, protective housing)



Autoclavable only

Level 2 Semi-standard containers

- Various lengths
- Container: beta door, beta flange
- Container body 1mm thick
- Handles: selection available, attached to flange
- Autoclavable or H₂O₂ sterilizable (non-perforated bottom, protective housing)
- Option: inserts

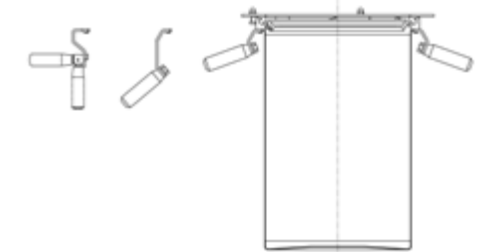


Autoclavable

H₂O₂ sterilizable

Level 3 Customizable containers using pre-designed parts

- Various lengths
- Container: Beta door, beta flange
- Container body 2 mm thick
- Handles: selection, may be attached anywhere on the container
- Autoclavable or H₂O₂ sterilizable
- Option: dismantlable bottom
- Option: additional clamp, may be attached to side and/or bottom
- Option: inserts

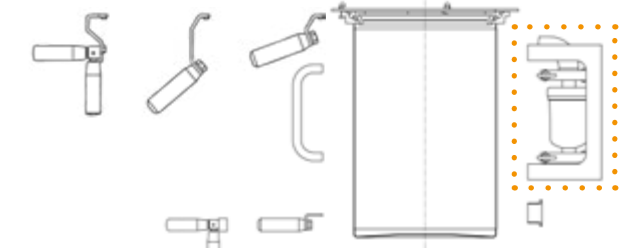

Autoclavable 2 options

H₂O₂ sterilizable

Dismantlable bottom

Level 4 Customized containers

Inserts, dimensions, parts and finishing to customers' requirements



Customized autoclavable stainless steel and plastic insert

Customized container