

Maquet Yuno II for orthopedic and traumatology procedures

Pave the way for high-class techniques


Proper positioning can make all the difference in orthopedic procedures, enabling less invasive techniques that can improve patient outcomes and accelerate recovery. With Maquet Yuno II, Getinge has set a new standard to support specialized orthopedic and traumatology applications.



Maquet Yuno II Mobile Operating Table is designed to facilitate minimally invasive direct anterior approaches for total hip replacements. A newly developed carbon fiber bar with ball joint can be moved both horizontally and vertically for optimal surgical access. When paired with the screw tension device, which enables fine tractions of 12 mm per turn of the handle, the surgeon has the precision required for optimal outcomes.

Key benefits

- **Improved patient safety:**
Automatic lock function prevents injury from accidental release of the carbon fiber bar.
- **Comfortable handling:**
Bars and devices are easy to slide and control, even after the table is fully prepared.
- **Ergonomic comfort:**
Best-in-class height range adjustment for precise and comfortable positioning (26.4" - 50.0").
- **More flexibility:**
Elongation adapter for hip scopes increases the range of possible procedures.
- **Color coded:**
Clearly marked functions reduce risk of error from miscommunication between surgeon and staff.
- **Universal usage:**
The three Maquet Yuno II Table Tops are compatible with our existing accessories for use in all surgical procedures, maximizing OR versatility.

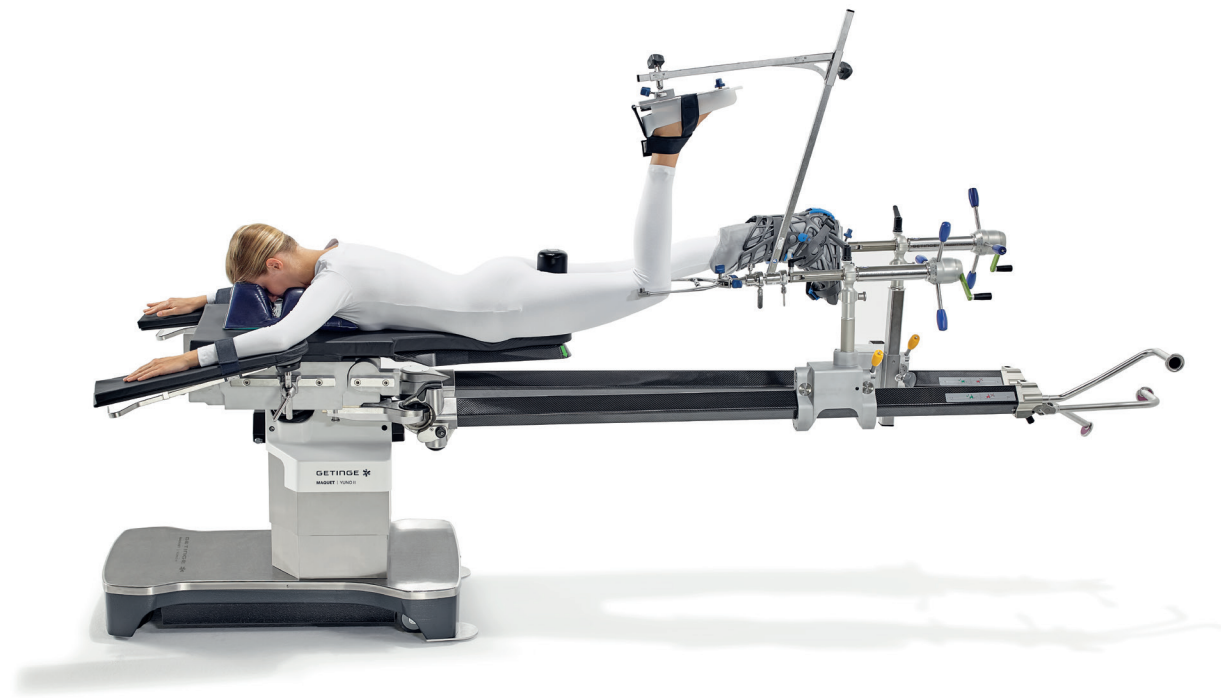

Application scope examples of Maquet Yuno II

- Minimally invasive direct anterior approach for hip surgery
- Hip arthroscopy
- Femur fracture
- Tibia fracture
- Pelvic fracture
- Knee surgery e.g. arthroplasty/arthroscopy
- Spine surgery e.g. lumbar spine
- Shoulder surgery
- Neurosurgery
- General surgery

Getinge is a registered trademark of Getinge AB, its subsidiaries, or affiliates in the United States or other countries. · Copyright 2019 Getinge AB or its subsidiaries or affiliates · All rights reserved

Getinge · 1777 East Henrietta Road · Rochester, NY 14623 · 1.800.475.9040
Sales Office, Canada · 90 Matheson Blvd West, Suite 300 · Mississauga, Ontario L5R 3R3 · Canada

www.getinge.com

GETINGE 

MCV00085406 REVA EXP 02/21